



Auto-Import Images into Metron

Copy photographs from your digital device.



877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com
METRON SUPPORT SERVICES. All rights reserved.

Auto-Import Images into Metron


Copy photographs from your digital device.

(PC on which Metron X-Rays are acquired)


YOUR METRON ADVANCED IMAGING SOFTWARE IS CONTINUALLY BEING ENHANCED. WE LISTEN TO YOU, STRIVING TO ADD AND ENHANCE THE PRODUCT AND SERVICE TO ASSIST YOU WITH THE DIAGNOSIS AND TREATMENT OF YOUR PATIENTS.

Setting up Image Location on your Desktop.


1. Right click anywhere on your desktop...click **"New"**, then scroll to **"Folder"**.

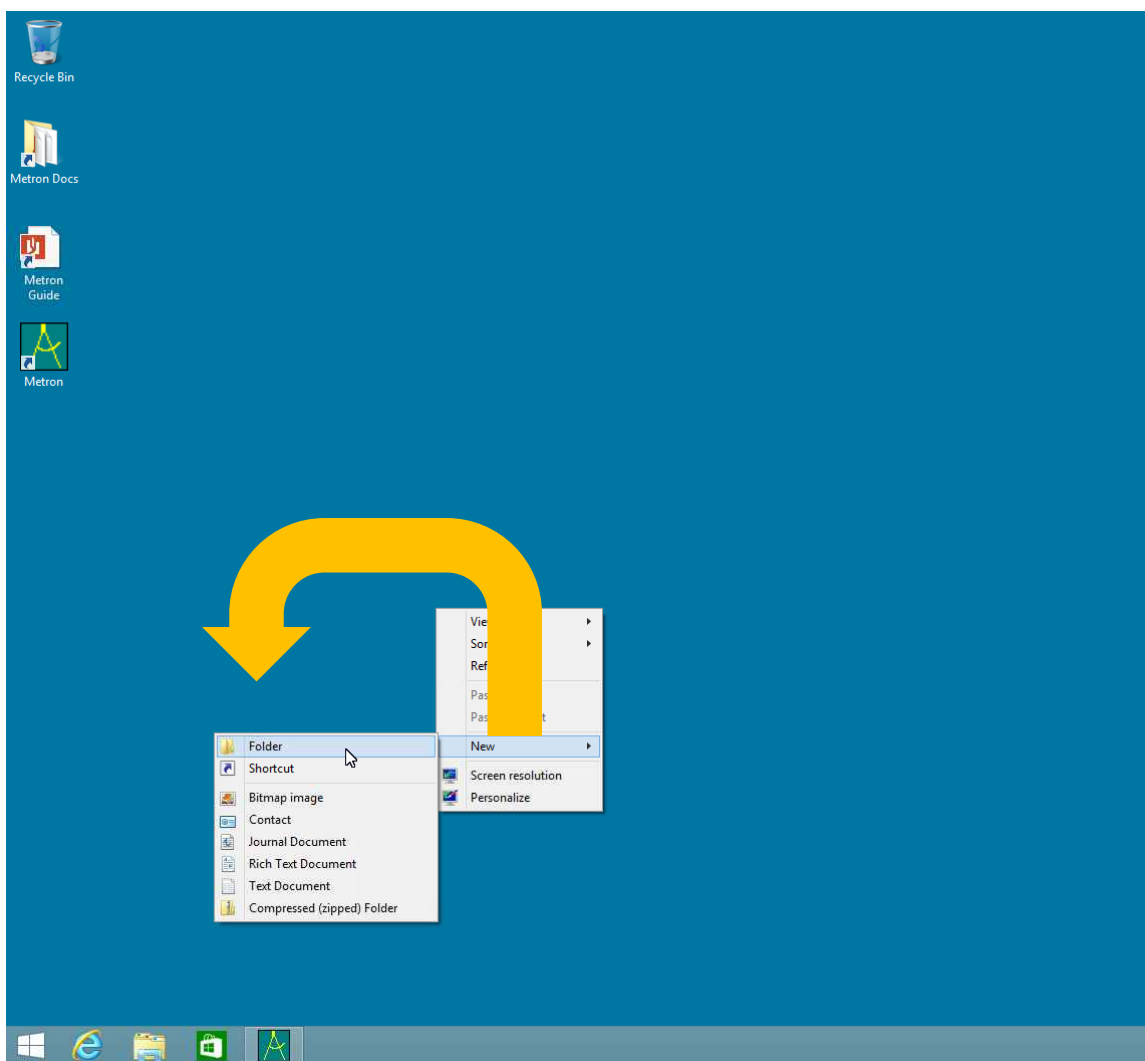

2. Rename the folder. In this example, **Dental Images**, was chosen.

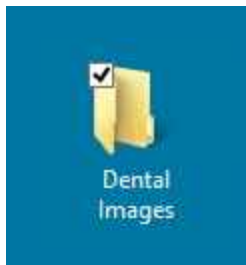


3. Double click the folder to open it. Single click in the path field and copy the contents.

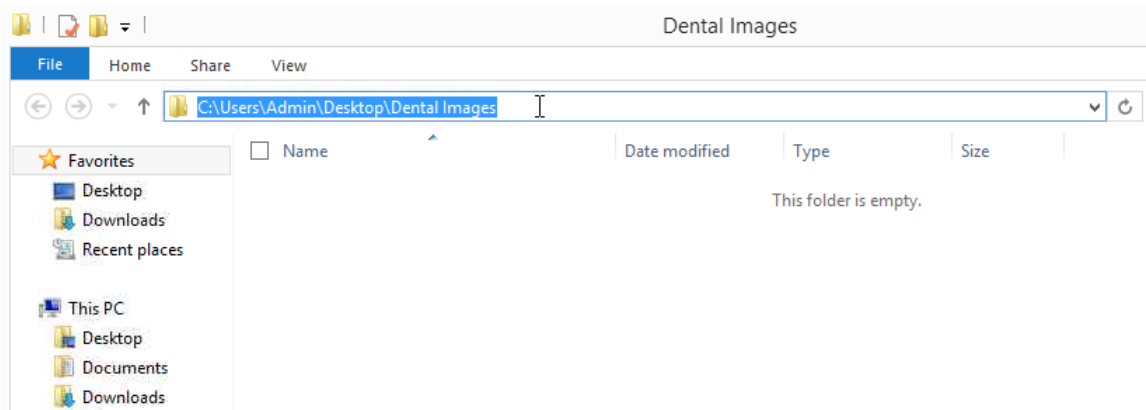


Setting up Metron for Auto-Import.

1. Launch Metron.
2. Click the “**System**” tab, then below “**Preferences**”, select “**Monitored Folders**”.

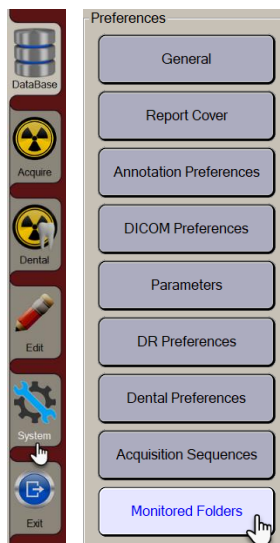

3. In the “**Monitored Folder #1**” area, click the “**Enable**” checkbox and in the “**Nickname**” field enter Dental Images (as an example).

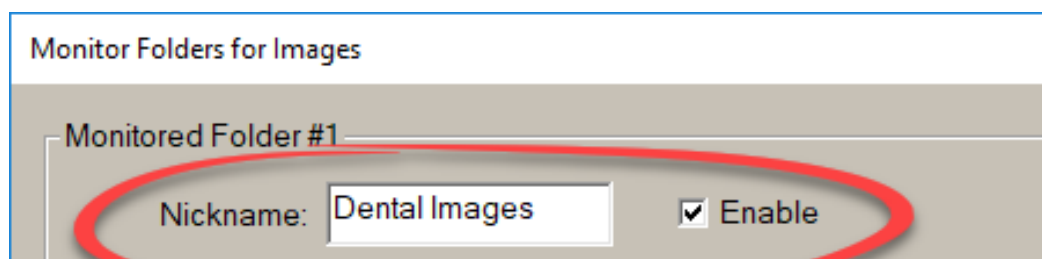

Monitor Folders for Images


Monitored Folder #1

Nickname: ☒ Enable

4. Clear out the contents of the path field below “**Nickname**” and paste path (folder path for this example created on the desktop):

C:\Users\Admin\Desktop\Dental Images

Modality defaults to **Photograph** which is correct.

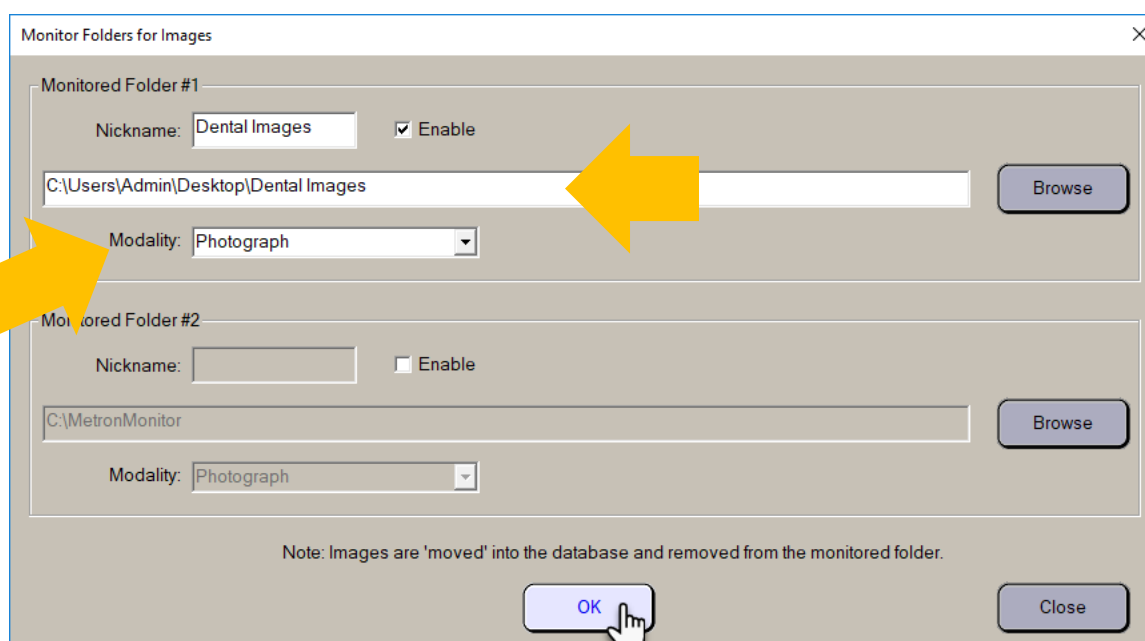

Monitor Folders for Images

Monitored Folder #1

Nickname: ☒ Enable

Modality:

Monitored Folder #2

Nickname: ☐ Enable

Modality:

Note: Images are 'moved' into the database and removed from the monitored folder.

5. Click “**OK**”

Using Auto-Import.

1. Take a photo with your camera device. (Camera, phone, tablet, etc.)
2. Plug in the device and copy the images to the newly created desktop folder “Dental Images”.
3. Create a new study in Metron or locate an existing study for these images to be placed.


Note: When photos are imported into Metron, they will always be imported under the owner, animal, and study date that is currently selected in the Database Browser

4. Once the images are in the desktop folder, the “**Receive Images...**” button in Metron will begin flashing red. Click the “**Receive Images...**” button.




- 5 Click the “**Import all 'Pending'**” button.

Note: Once an image has been imported it is removed from the Desktop folder. Also, if multiple images are available, they will all be imported and removed from this folder.

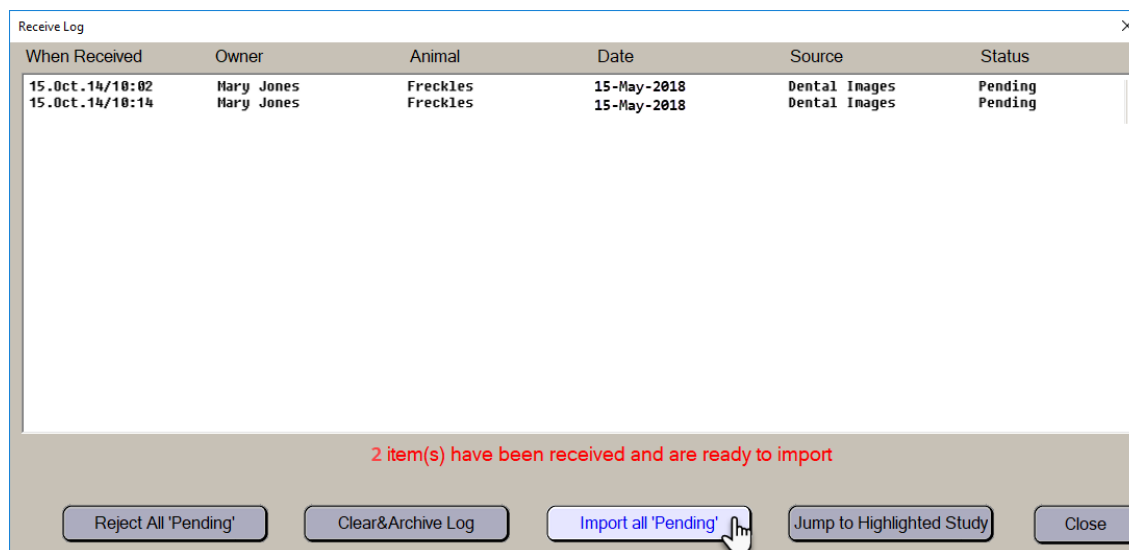

6. The photos you took earlier will now be present as part of the currently selected study.

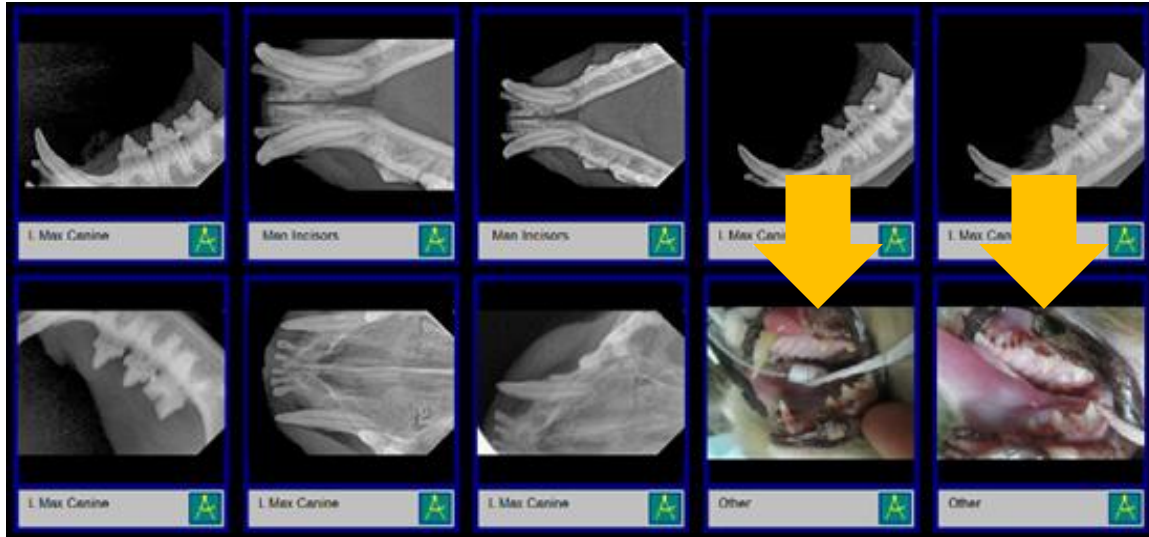

Now you can manipulate the photographs using all the available Metron tools. These photographs are also available to add to your dental report or dental forms.



Contact us if you have additional questions.
We can be reached at metron@metronsupportservices.com or 1-877-638-3868.