



Configure Metron for Fasig-Tipton

Version 8



877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com
METRON SUPPORT SERVICES. All rights reserved.


Configure Metron for Fasig-Tipton


Version 8

YOUR METRON ADVANCED IMAGING SOFTWARE IS CONTINUALLY BEING ENHANCED. WE LISTEN TO YOU, STRIVING TO ADD AND ENHANCE THE PRODUCT AND SERVICE TO ASSIST YOU WITH THE DIAGNOSIS AND TREATMENT OF YOUR PATIENTS.

- 1) Set up some choices by first clicking the “System” tab, in “**Preferences**” select “**DICOM Preferences**”

A) The “Default Tags” Tab

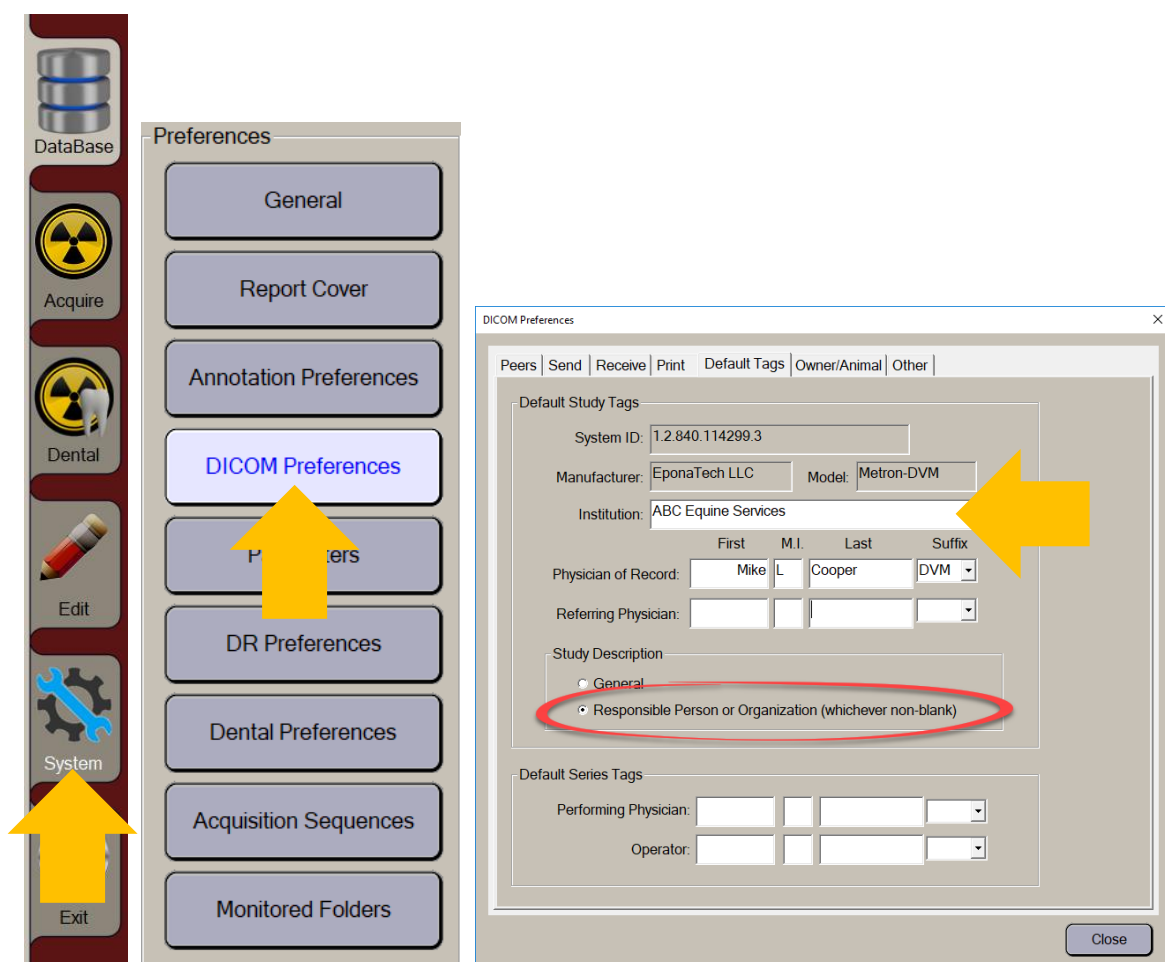

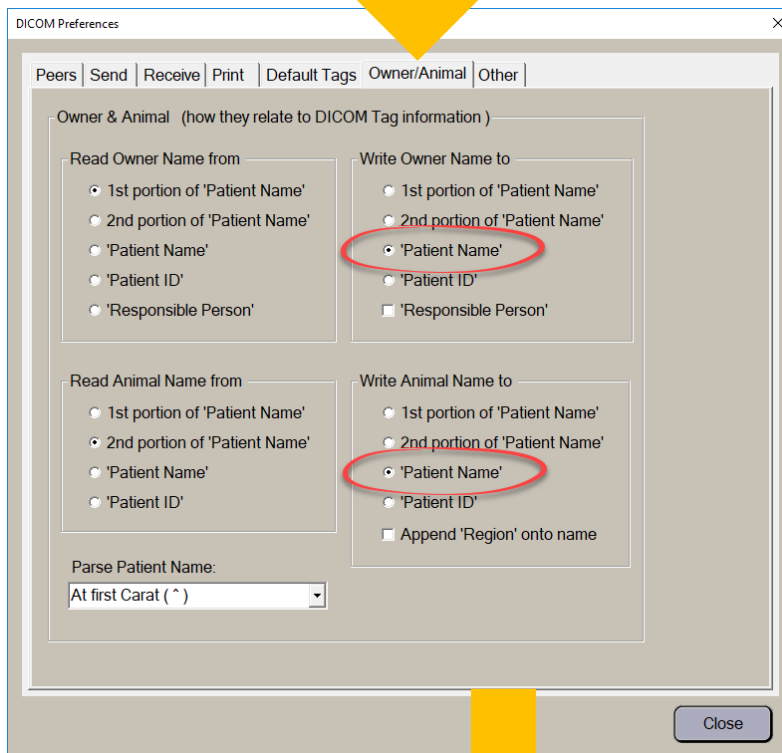

The screenshot shows the 'Owner/Animal' tab of the DICOM Preferences dialog. A yellow arrow points to the top of the dialog. The 'Owner & Animal (how they relate to DICOM Tag information)' section contains four groups of radio buttons. In the 'Write Owner Name to' and 'Write Animal Name to' groups, the 'Patient Name' option is selected and circled in red. The 'Parse Patient Name' dropdown is set to 'At first Carat (^)'.

DICOM Preferences

Peers | Send | Receive | Print | Default Tags | Owner/Animal | Other

Owner & Animal (how they relate to DICOM Tag information)

Read Owner Name from

- ☐ 1st portion of 'Patient Name'
- ☐ 2nd portion of 'Patient Name'
- ☐ 'Patient Name'
- ☐ 'Patient ID'
- ☐ 'Responsible Person'

Write Owner Name to

- ☐ 1st portion of 'Patient Name'
- ☐ 2nd portion of 'Patient Name'
- ☒ 'Patient Name'
- ☐ 'Patient ID'
- ☐ 'Responsible Person'

Read Animal Name from

- ☐ 1st portion of 'Patient Name'
- ☐ 2nd portion of 'Patient Name'
- ☐ 'Patient Name'
- ☐ 'Patient ID'

Write Animal Name to


- ☐ 1st portion of 'Patient Name'
- ☐ 2nd portion of 'Patient Name'
- ☒ 'Patient Name'
- ☐ 'Patient ID'
- ☐ Append 'Region' onto name

Parse Patient Name:

At first Carat (^)

Close

B) The “Owner/Animal” Tab

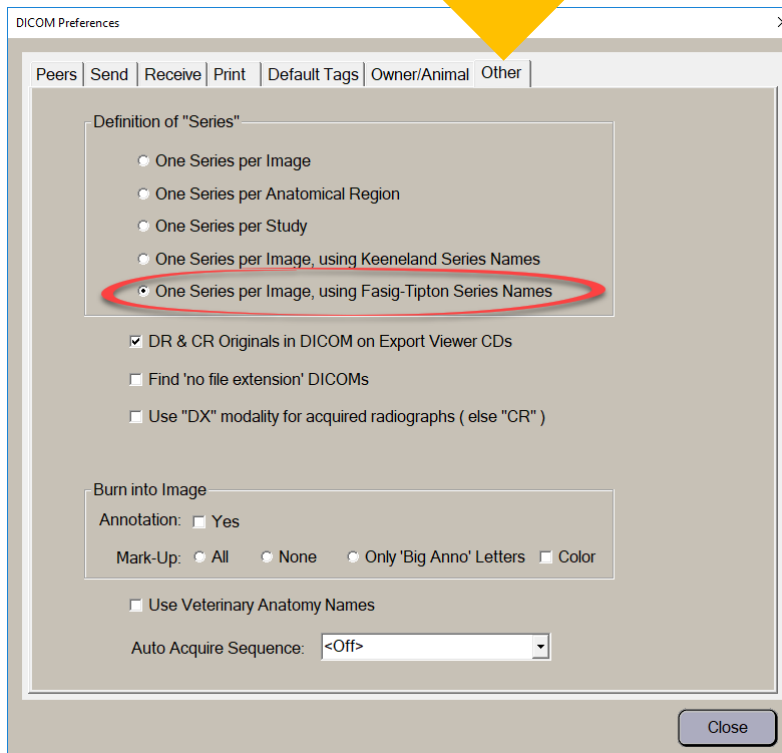

The screenshot shows the 'Other' tab of the DICOM Preferences dialog. A yellow arrow points to the bottom of the previous dialog. The 'Definition of "Series"' section has five radio button options, with the last one, 'One Series per Image, using Fasig-Tipton Series Names', selected and circled in red. Below this are three checkboxes: 'DR & CR Originals in DICOM on Export Viewer CDs' (checked), 'Find 'no file extension' DICOMs' (unchecked), and 'Use "DX" modality for acquired radiographs (else "CR")' (unchecked). The 'Burn into Image' section has an 'Annotation' checkbox (unchecked), a 'Mark-Up' section with three radio buttons ('All' selected), and a 'Color' checkbox (unchecked). There is also a 'Use Veterinary Anatomy Names' checkbox (unchecked) and an 'Auto Acquire Sequence' dropdown set to '<Off>'.

DICOM Preferences

Peers | Send | Receive | Print | Default Tags | Owner/Animal | Other

Definition of "Series"

- ☐ One Series per Image
- ☐ One Series per Anatomical Region
- ☐ One Series per Study
- ☐ One Series per Image, using Keeneland Series Names
- ☒ One Series per Image, using Fasig-Tipton Series Names

☒ DR & CR Originals in DICOM on Export Viewer CDs

☐ Find 'no file extension' DICOMs

☐ Use "DX" modality for acquired radiographs (else "CR")

Burn into Image

Annotation: ☐ Yes

Mark-Up: ☒ All ☐ None ☐ Only 'Big Anno' Letters ☐ Color


☐ Use Veterinary Anatomy Names

Auto Acquire Sequence: <Off>

Close

C) The “Other” Tab


- 2) Set up the “**Patient Name**” as Fasig-Tipton says as the Dam Name.

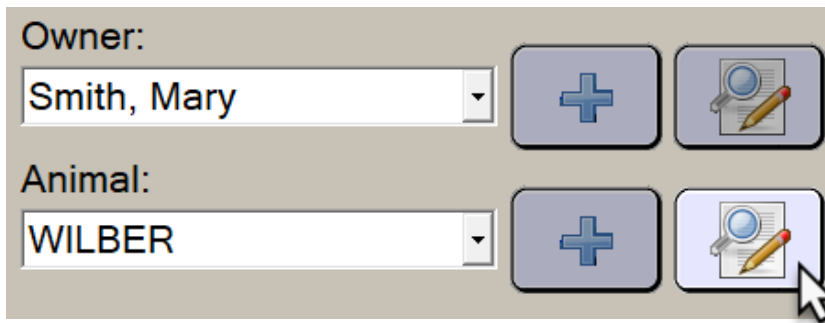

Owner:
Smith, Mary

Animal:
WILBER

- 3) Set up the “**Patient ID**” as Fasig-Tipton wants
(Click the “Edit” button next to the Animal Name)

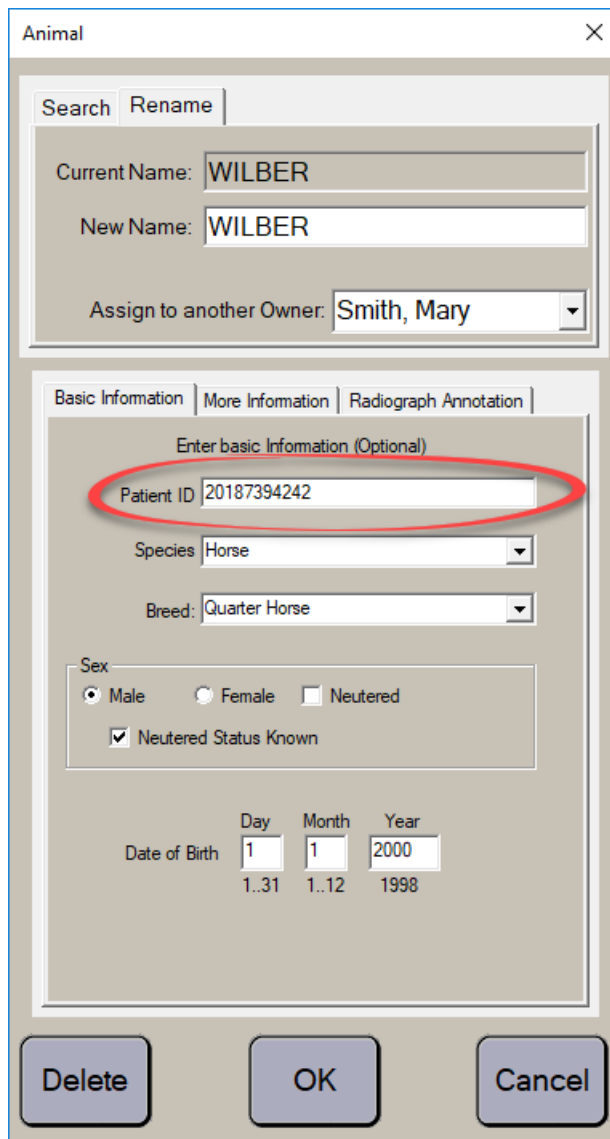

Animal

Search Rename

Current Name: WILBER

New Name: WILBER

Assign to another Owner: Smith, Mary

Basic Information More Information Radiograph Annotation

Enter basic Information (Optional)

Patient ID 20187394242

Species Horse

Breed: Quarter Horse

Sex

☒ Male ☐ Female ☐ Neutered

☒ Neutered Status Known

Date of Birth Day Month Year
1 1 2000
1..31 1..12 1998

Delete OK Cancel

4) Set up the “**Responsible Org**” (the Vet or Practice taking the images) Click the “More Information” tab after you have entered the Patient ID

Animal

Search Rename

Current Name: WILBER

New Name: WILBER

Assign to another user: Smith, Mary

Basic Information More Information Radiograph Annotation

Information used in DICOM file creation (Optional)

Responsible Person

First M.I. Last Suffix

Role: OWNER

Responsible Org: Equine Veterinary Services

Other Patient ID

Patient ID Type: RFID

ID Issuer

Breed Registry Info

Breed Reg

Breed Reg Number

Delete OK Cancel

Note: Items 2, 3, and 4 will need to be set up for each horse. Then burn a DVD in “**DICOM Images Only**” mode for one animal. Put each animal on a separate DVD.



Contact us if you have additional questions.
We can be reached at metron@metronsupportservices.com or 1-877-638-3868.