



Creating a Collection in Metron 8



877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com
METRON SUPPORT SERVICES. All rights reserved.

Creating a Collection in Metron 8


YOUR METRON ADVANCED IMAGING SOFTWARE IS CONTINUALLY BEING ENHANCED. WE LISTEN TO YOU, STRIVING TO ADD AND ENHANCE THE PRODUCT AND SERVICE TO ASSIST YOU WITH THE DIAGNOSIS AND TREATMENT OF YOUR PATIENTS.

The Metron Collection feature makes it easy to display images side by side that are from different moments in time for the same patient or for different patients. For example, if an animal breaks its leg, you could use the Metron Collection feature to compare the animal's injury from onset through the recovery process.


Using collections, you can add an X-Ray from when the leg was initially broken, a postoperative view, and then another X-Ray taken a month later once the bone begins to set.

Using the Collections Feature in Metron 8

1. Navigate to the first Owner/Animal/Study date that has radiograph(s) you would like to compare to radiographs taken on another study date. Click on each image you would like to be part of the "Collection" - now indicated by a red border.

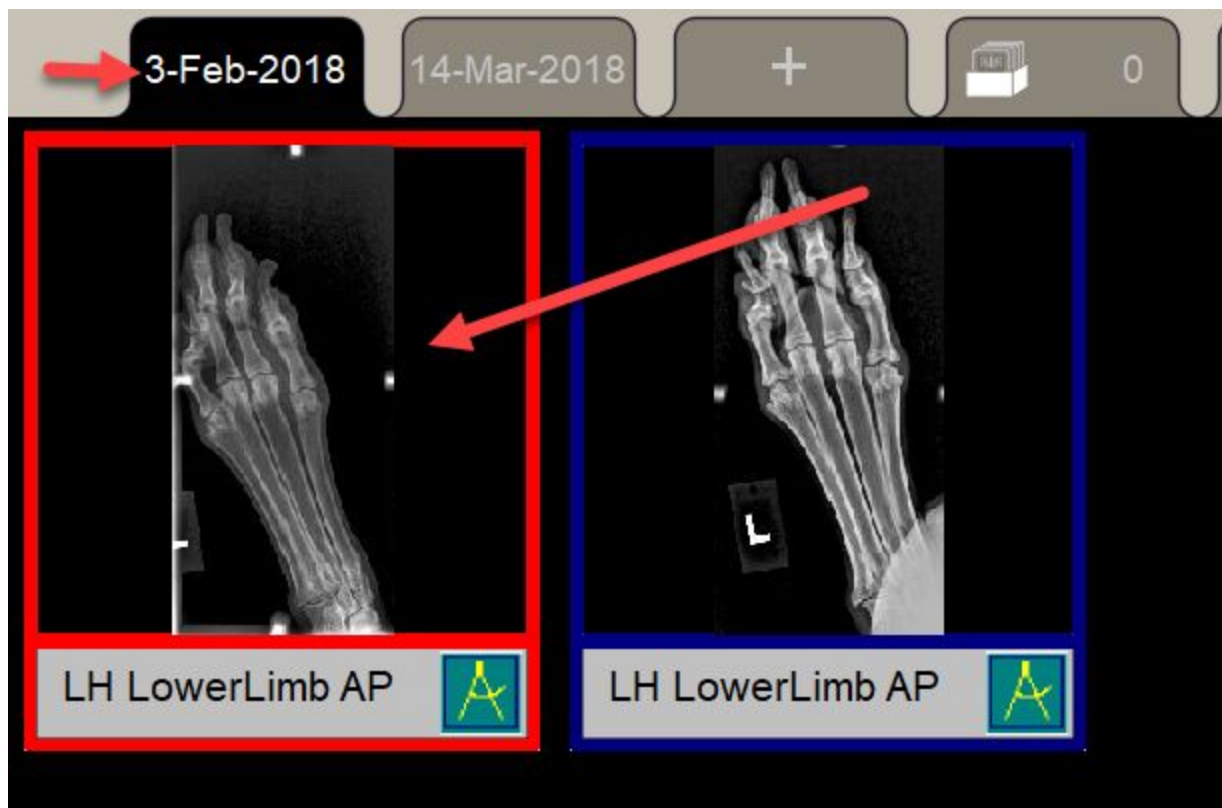


2. Click the **“Add to Collection”** button located on the right-hand side of the Metron button panel to add the image(s) to the collection.

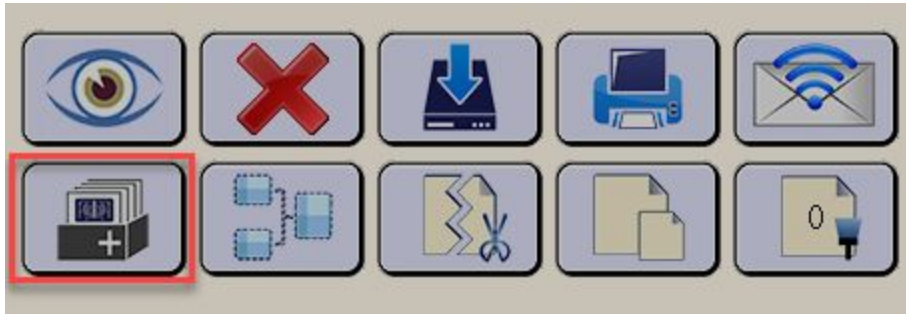

3. Navigate to the second study date with radiograph(s) you would like to compare with the first study date. Click on the image, or images, so they have a red border.



4. Click the **“Add to Collection”** button located on the right-hand side of the Metron button panel to add the image(s) to the collection.

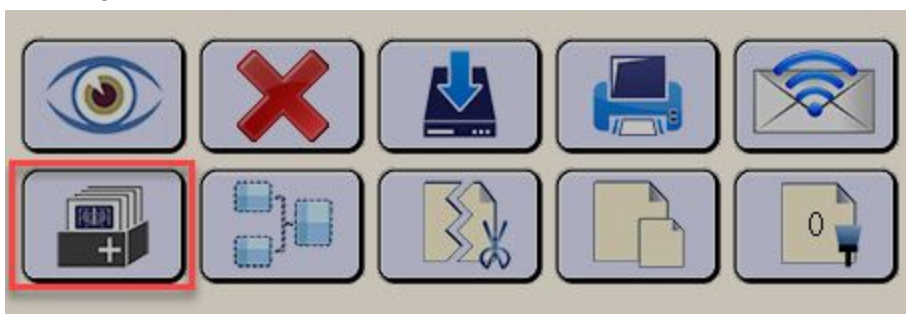


5. Click the **"Collections"** tab at the top of the screen (near the study dates) to enter the Collections view.
(Note: if you have multiple study dates that fill the screen, you may have to click the blue "right" arrow to locate the Collections tab button. Also, note the collections button has a number on it representing the number of images you have added to the collection).



6. You will now see all of the images you added to the collection. In this screen, you may highlight the images so that they have a red border and click the **"View"** button.



7. After you are completed comparing radiographs and are ready to reset the collection, click the “**Reset Collection**” button while in the Collection mode - this will “reset” and remove all images in the current collection.

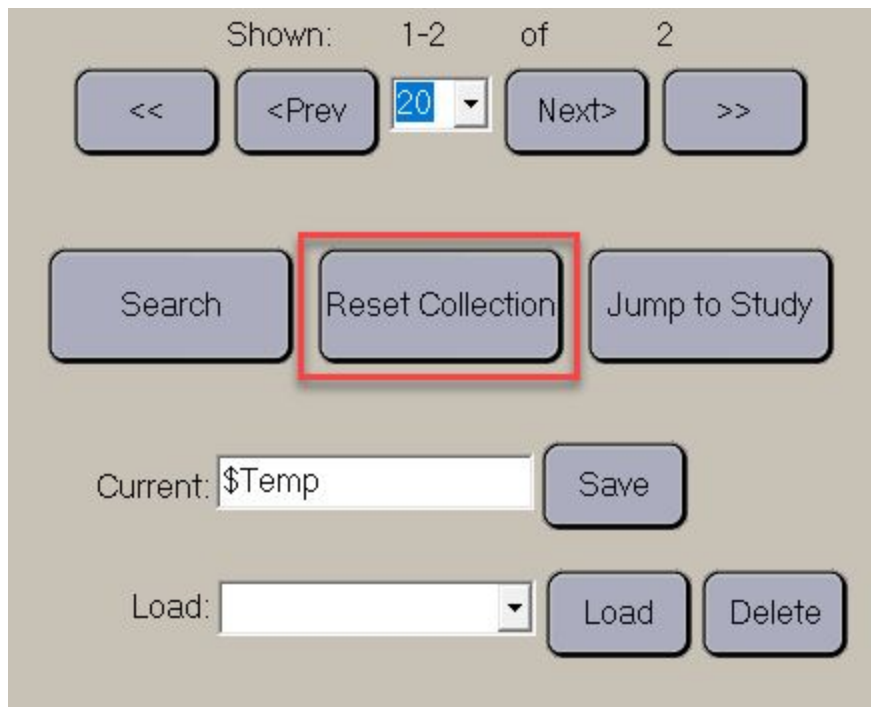

The Metron Collections feature is very easy to use and is a great tool for comparative analysis.

If you have any questions or would like to learn more about some of the advanced Metron 8 Collection features, contact your Support Services Team.



Contact us if you have additional questions.
We can be reached at metron@metronsupportservices.com or 1-877-638-3868.