



Learn More

How to Declare a Metron Database

Provides access to view Metron images on workstations



877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com

METRON SUPPORT SERVICES. All rights reserved.

How to Declare a Metron Database

Provides access to view Metron images on workstations


YOUR METRON ADVANCED IMAGING SOFTWARE IS CONTINUALLY BEING ENHANCED. WE LISTEN TO YOU, STRIVING TO ADD AND ENHANCE THE PRODUCT AND SERVICE TO ASSIST YOU WITH THE DIAGNOSIS AND TREATMENT OF YOUR PATIENTS.


Declaring an existing database in Metron

Do you have multiple copies of the Metron software in use within your practice? Do you need access to digital images located at various workstations? *Reminder:* each Metron software license you own includes up to two viewstation licenses. This is a reason many people would need to be able to declare a database.


Listed below explains how to add an existing database from one computer to another so images can be viewed and utilized at multiple locations.

- Open Metron
- Right Side Tabs, Click “**System**”, Below Database, Click “**New Database**”
- Go to drive letter with new database
 - You can have your IT person map a drive to your acquisition computer's C:\Epona folder
- Highlight database name
- Click “**Declare an Existing Database**”

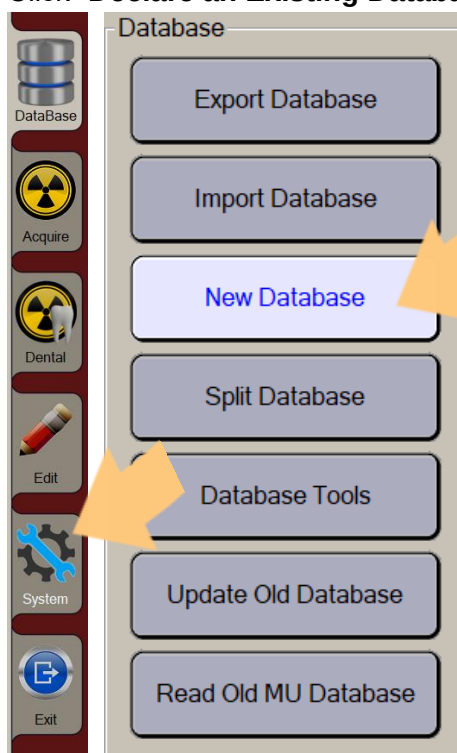


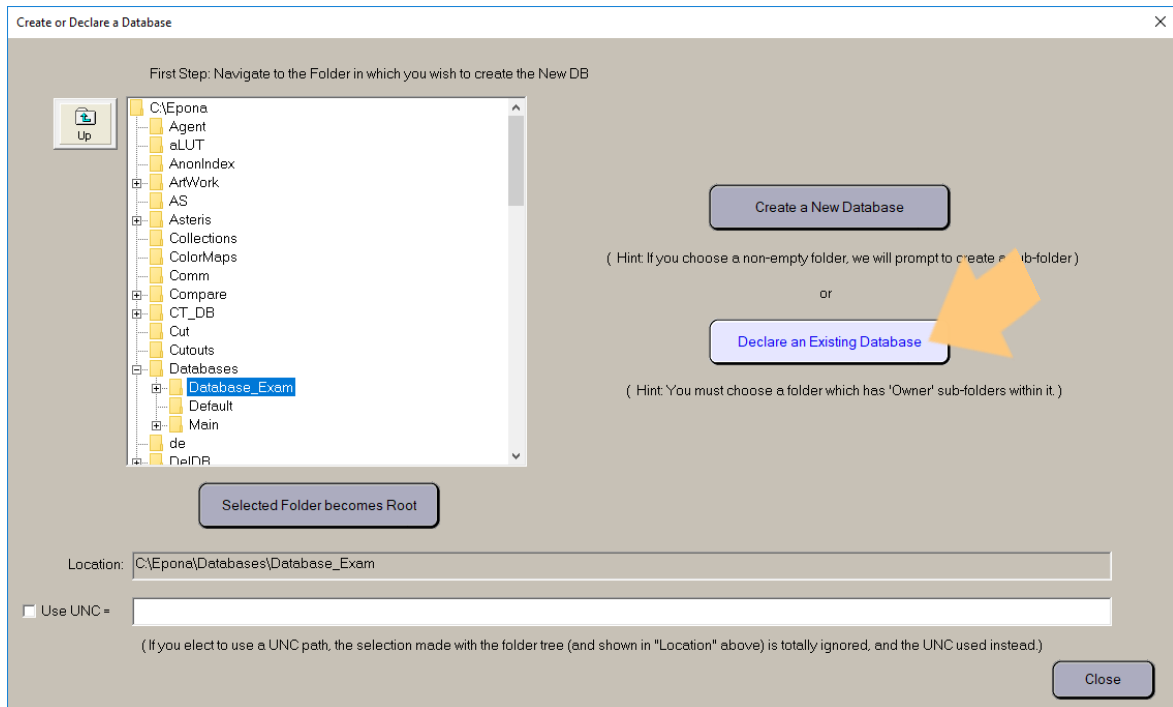


- You should see a list of Owners – Click “Yes”

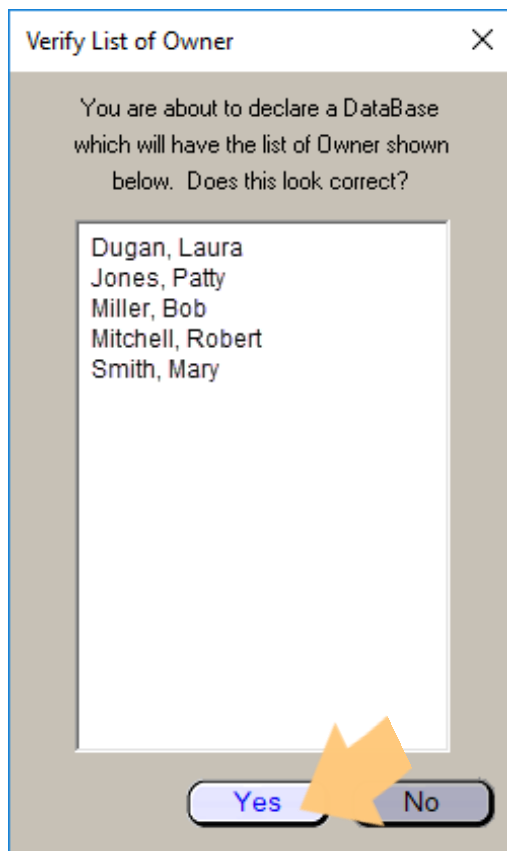

- Enter name of Database – Click “OK”.

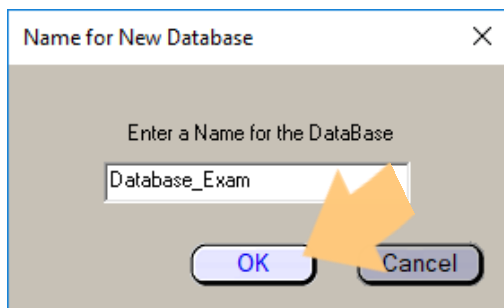

- Click the **Database tab** to return to the Database Panel, then **click the Database tab again** to choose a different database. Dropdown in upper right corner will display choices including your newly named database...Database_Exam. Select and click “Close” to use and view this database.

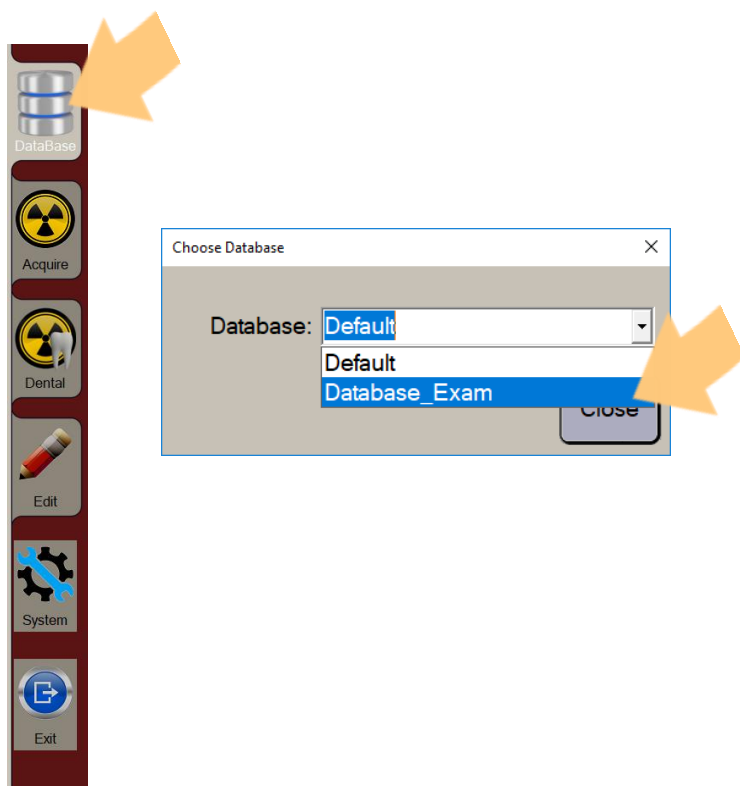

Now that you have your drives mapped, you are able to view and use digital images from multiple workstations in your practice.

Contact us if you have additional questions.
We can be reached at metron@metronsupportservices.com or 1-877-638-3868.