




How to
Acquire in
METRON
Quick Start Guide

To get started using Metron 8, let's begin by acquiring your first image!

Double Click the Metron Icon on the Desktop to Launch Metron 8.

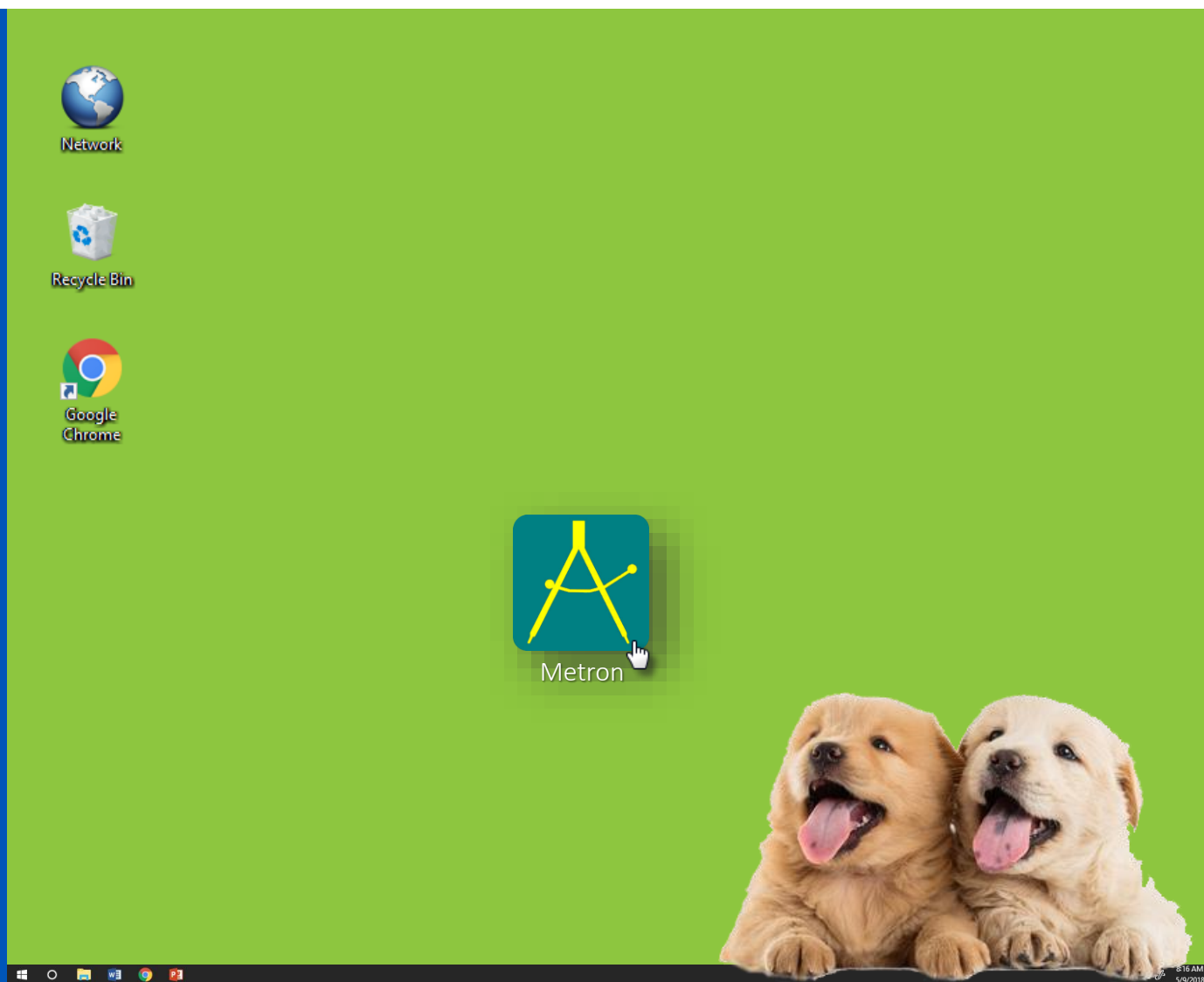

The screenshot displays the Metron 8 software interface. At the top, there's a header bar with a date filter set to '27-Mar-2018' and zoom controls. Below this is a grid of five X-ray images: 'L Abdomen Lateral', 'RF LowerLimb Later', 'RH Pelvis DV', and 'LF LowerLimb Later'. Each image has a small 'A' icon in the bottom right corner. To the right of the grid is a sidebar containing a 'Date Filter' dropdown, 'Owner' and 'Animal' dropdowns (set to 'Smith, Mary' and 'Snickers' respectively), and a grid of icons for various functions like 'Acquire', 'Dental', 'Edit', 'System', and 'Exit'. At the bottom of the sidebar is a red 'EMERGENCY' button. The text 'Metron 8 boasts a user friendly database. Large buttons and tabs for ease of use with touch screen devices.' is overlaid on the left side of the interface. Below this text, the question 'So first, what do those new tabs do?' is displayed.

Metron 8 boasts a user friendly database. Large buttons and tabs for ease of use with touch screen devices.

So first, what do those new tabs do?

The tabs on the right side quickly get you where you need to go with one click or tap.



Database



Acquire



Acquire with Dental Sensor



Edit a Single Image



System Folders and Settings



Exit Metron Program

The top tabs display your **Image Studies**.

The **Collections** tab is also located here.

Create a Collection

to view a group of images from different studies side-by-side.



Collections

Create New Study Date

Image Studies By Date Tab

View Thumbnails



Client Study Creation


To Create a New Client:

1. Select  to the right of Owner.
2. Add the New Owner's Name* exactly how it should appear in the Metron Database.
3. Click **OK**.
4. Select  to the right of Animal.
5. Add New Animal's available information.
6. Click **OK**.

Search or Add to an Existing Client:

- A. Select the Owner using drop down arrow.
- B. Select the Animal. (if multiple)


*We recommend Last Name, First Name.

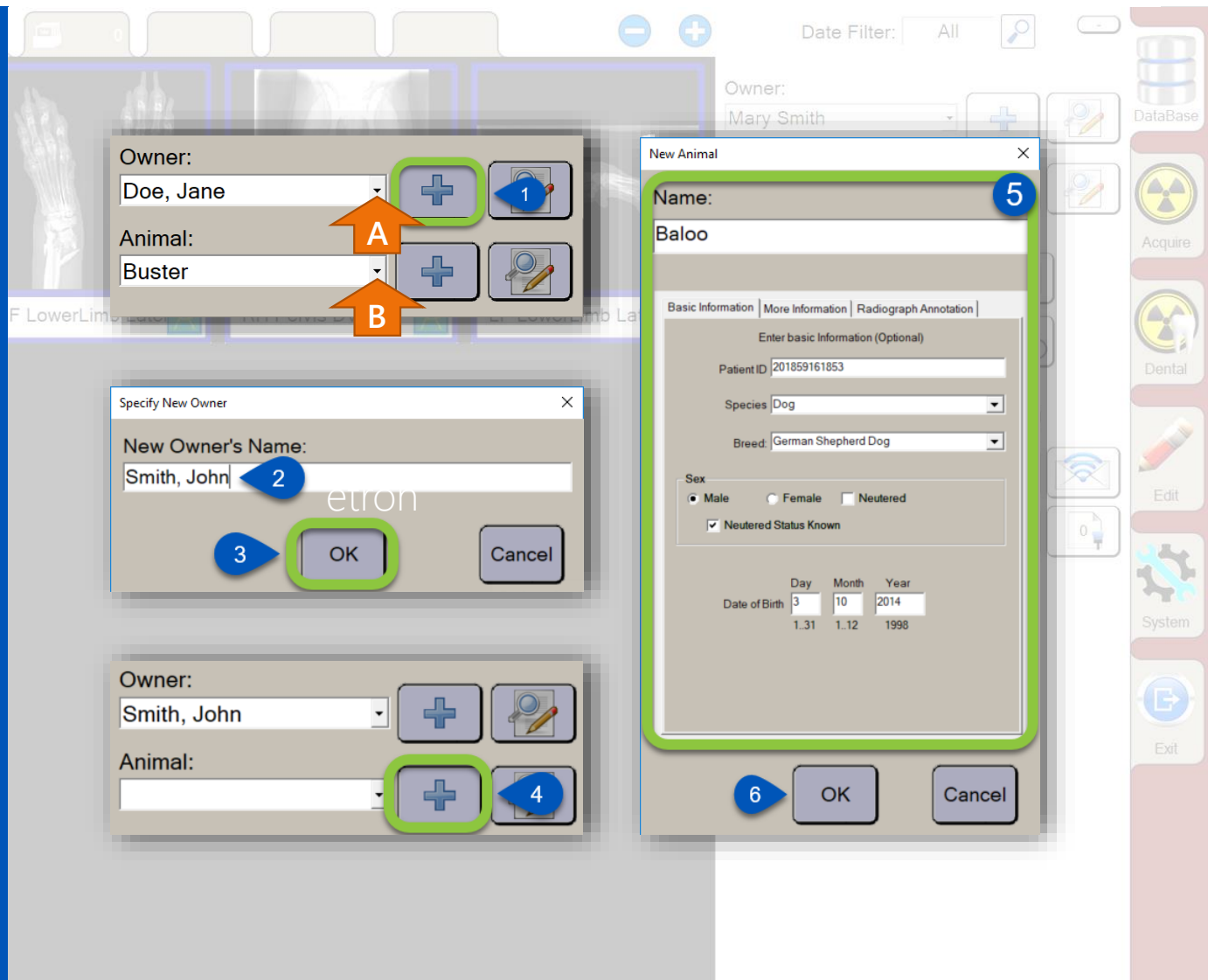

The screenshot illustrates the Metron 8 software interface for creating a new client study. The main window shows the 'Owner' field with 'Doe, Jane' and the 'Animal' field with 'Buster'. To the right of each field is a plus icon. A green circle with a blue '1' highlights the plus icon next to the Owner field. An orange arrow labeled 'A' points to the plus icon next to the Animal field. A green circle with a blue '2' highlights the 'Specify New Owner' dialog box, which contains the 'New Owner's Name' field with 'Smith, John' and 'OK' and 'Cancel' buttons. A green circle with a blue '3' highlights the 'OK' button. An orange arrow labeled 'B' points to the plus icon next to the Animal field. A green circle with a blue '4' highlights the plus icon next to the Animal field. A green circle with a blue '5' highlights the 'New Animal' dialog box, which contains the 'Name' field with 'Baloo', 'Basic Information' tab, 'Patient ID' field with '201859161853', 'Species' dropdown with 'Dog', 'Breed' dropdown with 'German Shepherd Dog', 'Sex' section with 'Male' selected, 'Neutered Status Known' checked, and 'Date of Birth' section with 'Day' 3, 'Month' 10, and 'Year' 2014. A green circle with a blue '6' highlights the 'OK' button at the bottom of the 'New Animal' dialog box.

Acquire an Image

1. Click the **Acquire Tab** to prompt the Acquire Screen.
2. Click **OK** to create the new study under today's date.

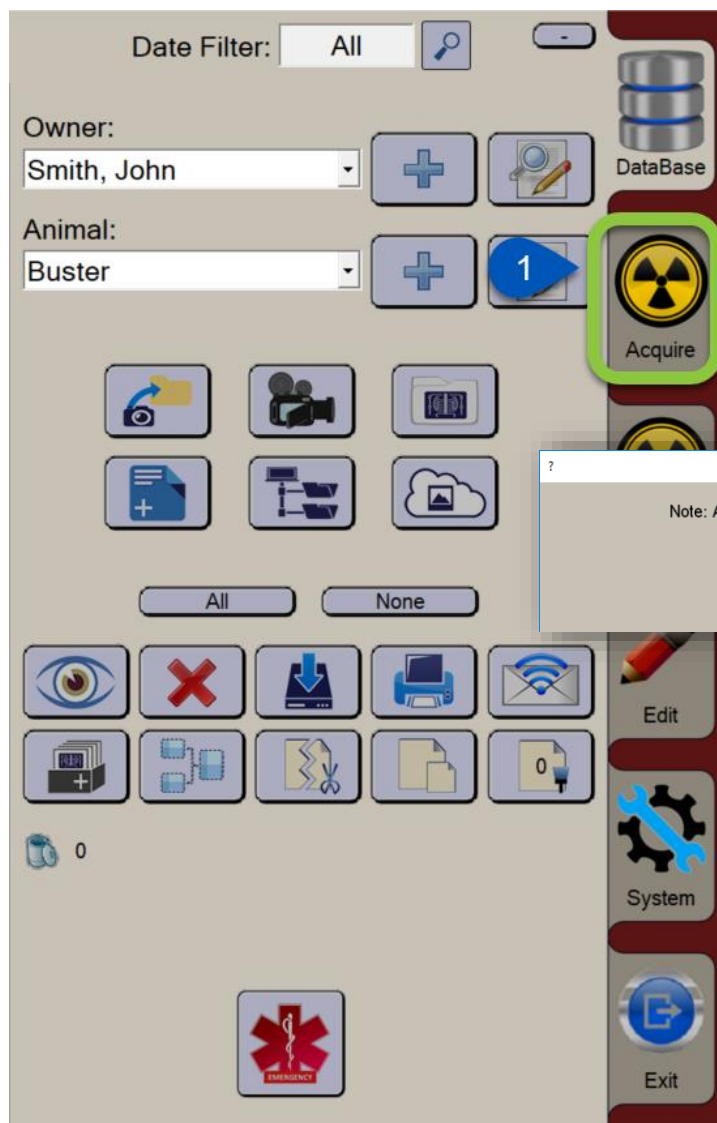


Acquire an Image

Metron will go from  to  when it is able to communicate with a CR/DR system.

1. Best practice, select **Region** to X-Ray*. (click dot-red dot indicates region)
2. Select **View**.
3. When  use your firing switch (hand switch/foot pedal) to fire the generator and acquire the image.

**Follow pop-up prompts if labelling your image after acquisition.*



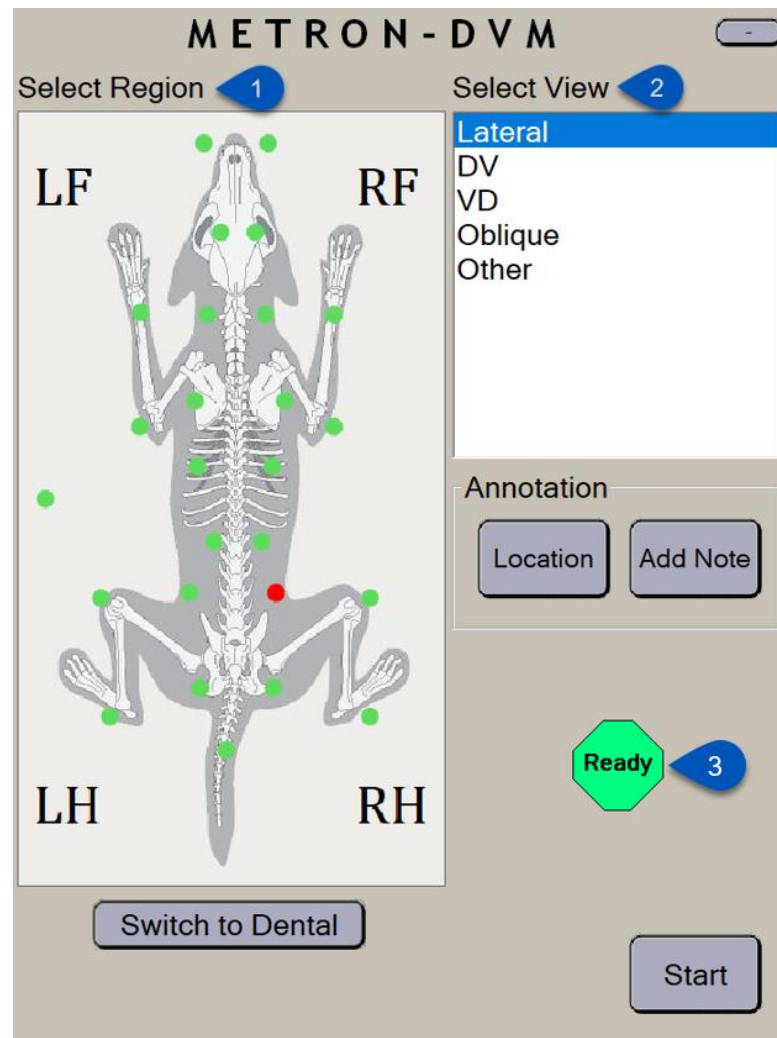


Edit an Image

After the X-Ray has been acquired, the image will be displayed shortly. The image (with any additional user edits) will be saved automatically.

1. Make all adjustments and edits to the image.*
2. Change the location of the annotation to best fit the acquired image (optional). Click **OK** when selected. Auto will choose the preferred option based on the image.
3.  **Accept** the image, then continue to the next shot as needed.
 **Reject** the image. Rejected images go to the trash can.



*Reference your [Metron 8 Button Guide](#) for image editing capabilities.

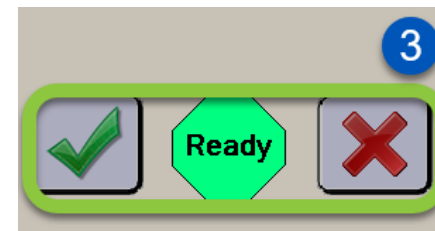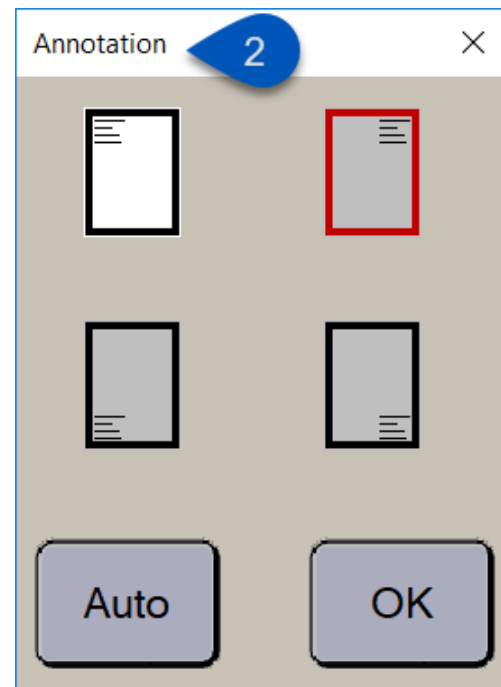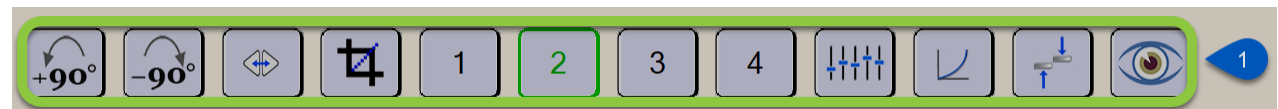

When all images have been acquired, click the **DataBase Tab** to return to all images in that study.

*Reference your [Metron 8 Button Guide](#) for additional features.

METRON8

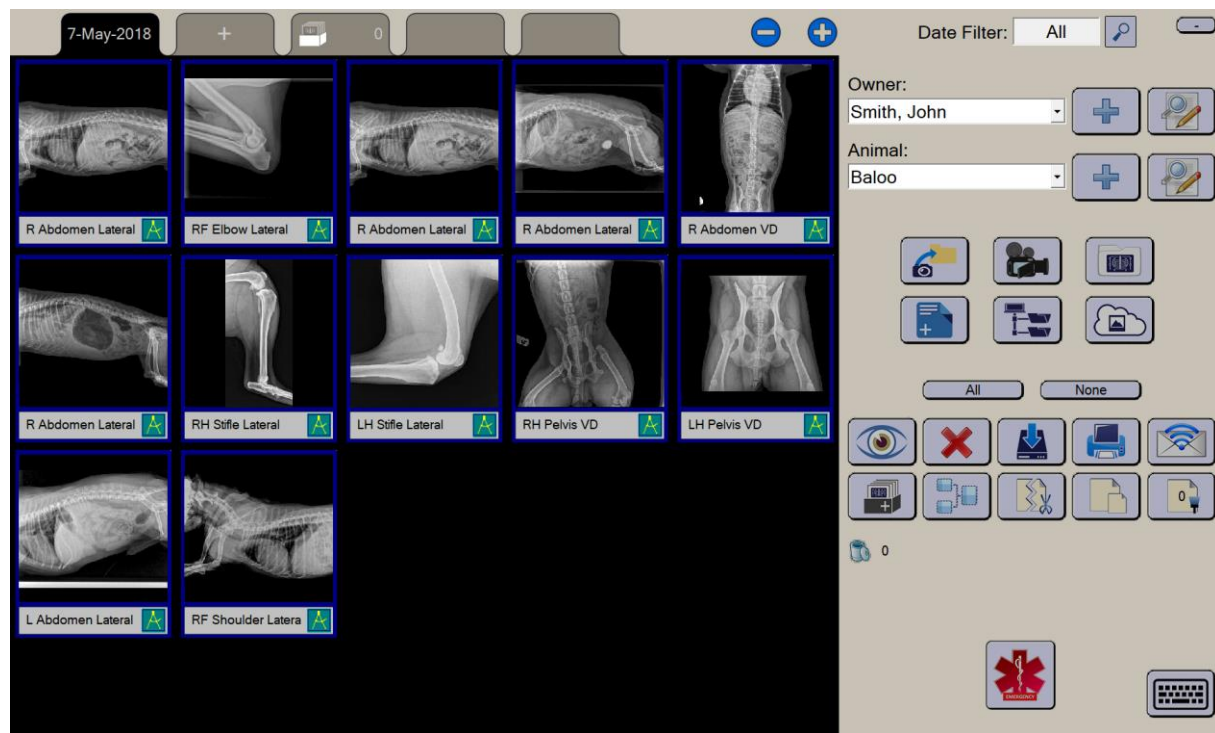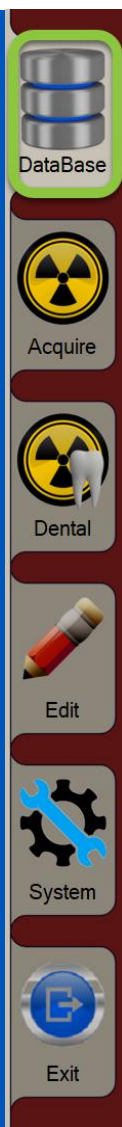

For Emergency Purposes Only



In the event of a situation when time is critical, click the **Emergency Button** to go to the Acquire Panel.

Please Note: Images acquired in Emergency Mode will be created based on the time the images were taken.

Edits to the Owner and Animal information should be completed as time permits to maintain accurate records.



Owner: Owner_15h35m24s
Animal: Animal_15h35m24s
Date: 9-May-2018 12:58:49
RF Foot Lateral
Patient DOB: 1-Jan-2000

This is how the annotation will appear on images in Emergency Mode until edits to the Owner and Animal are made.



Contact Us!



Contact us if you have additional questions.

1-877-638-3868

metron@metronsupportservices.com

Visit us anytime online at:

www.metronsupportservices.com



METR8N