



Import Images from VetRay to Metron

Just a few simple steps to additional features!



877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com
METRON SUPPORT SERVICES. All rights reserved.

Import Images from VetRay to Metron

Just a few simple steps to additional features!

YOUR METRON ADVANCED IMAGING SOFTWARE IS CONTINUALLY BEING ENHANCED. WE LISTEN TO YOU, STRIVING TO ADD AND ENHANCE THE PRODUCT AND SERVICE TO ASSIST YOU WITH THE DIAGNOSIS AND TREATMENT OF YOUR PATIENTS.

1. **Open VetRay** (you will need your VetRay dongle plugged in.)



2. Double Click the study you would like to import into Metron from the Patient List within VetRay.

Patient/Image Administration

Patient List | Owner | Animal | Images | Evaluation


Date: 11/17/2009 Search: Filter: All Sort: Last Name

✓ Status Images Not Proc.


Last Name	First Name	Tel.	Date of Birth	Receiving	Zip/City	Address
Ahl	Valerie			6/24/2008		
Bentz	Kim		12/30/1899	3/8/2006		
Blue Huron Sanctuary			12/30/1899	2/16/2006		
Bluschke	linda		12/30/1899	3/3/2006		
Brenner	Judy		12/30/1899	3/14/2006		
Breunig	Louis			5/16/2008		Mystic Valley Dairy
Brick	Troy		12/30/1899	3/13/2006		
Burcalow			12/30/1899	1/18/2006		
Burcalow	Becky		12/30/1899	2/6/2006		
Butler	Lynn		12/30/1899	3/6/2006		
Chapman	Ruth		12/30/1899	3/6/2006		
Colello	Beth		12/30/1899	3/13/2006		
Dave	Trussel		12/30/1899	2/24/2006		
Diercks	Kim		12/30/1899	2/22/2006		
Diercks	Kim		12/30/1899	2/22/2006		
Doolittle	Pamela			5/15/2008		
Eakins	Linda		12/30/1899	2/17/2006		

Patient: Diercks, Kim, 12/30/1899 711


- Click the **"Select All"** button to highlight all images. However, if you would only like to export a couple of images select the images you would like to export.

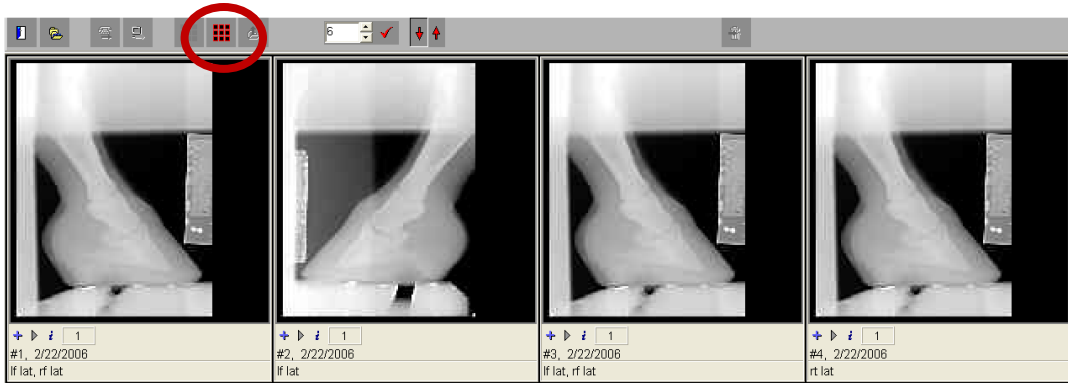

- Click the **"Load Selected Images"** button.

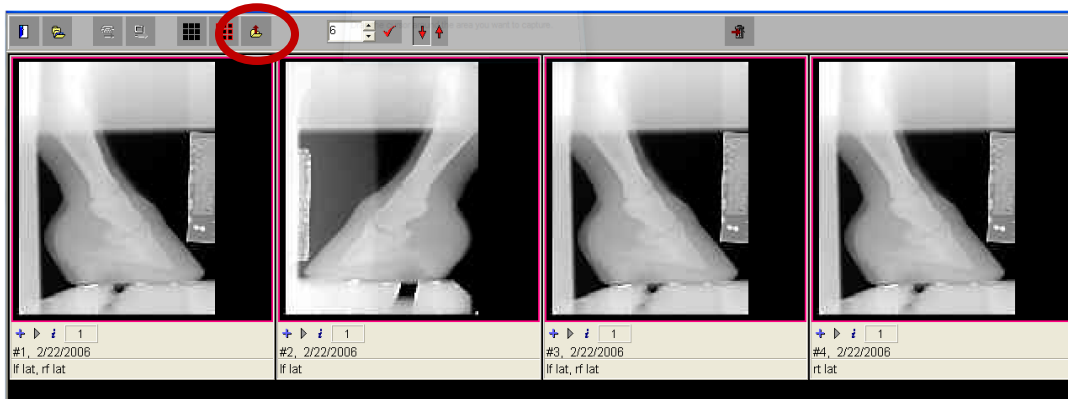


- Click **File**, followed by **Image Export** or press the F2 button on your keyboard.

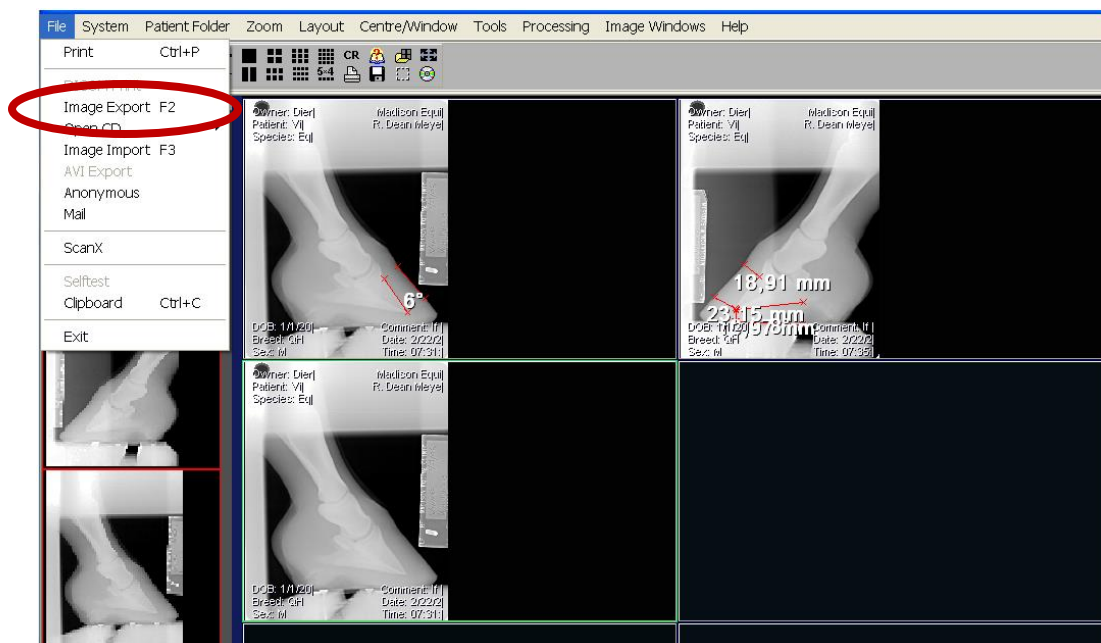


6. Select **"All Images"** and click the **"OK"** Button




7. Select a location you would like to save these images to. We created a folder on the C Drive called VetRay Images. Within that folder we created a folder for the Animal and then another for the Date. By doing this, we will know which images belong to which study. In this example, our location is C:\VetRay Images\Sample Horse#2\ Feb 13 2009

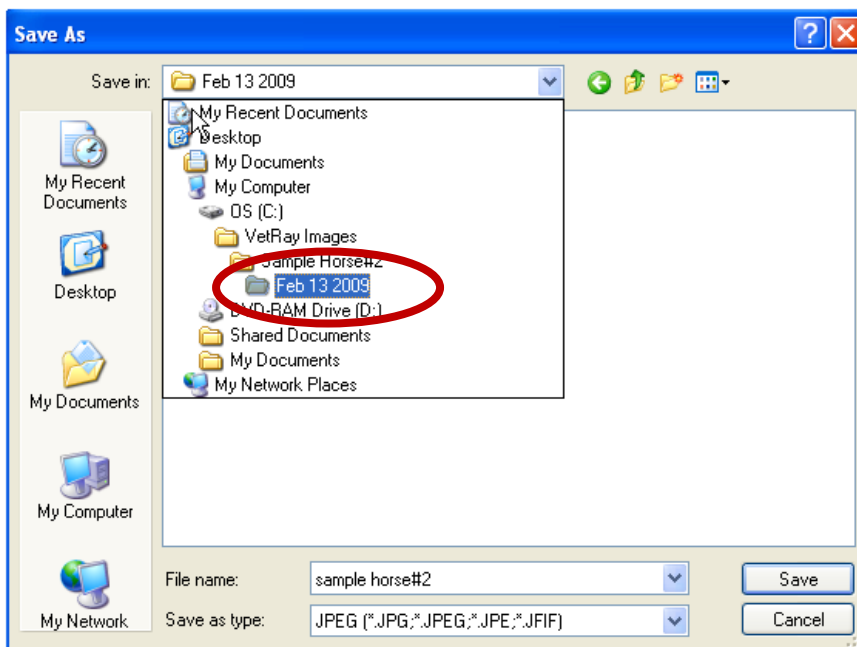


8. You will need to type the File name. We used the name of the animal as our file name. The image name will be the name you typed in plus 3 digit numbers which will change for each image. For example: sample horse#2000 the next image will be sample horse#2001.

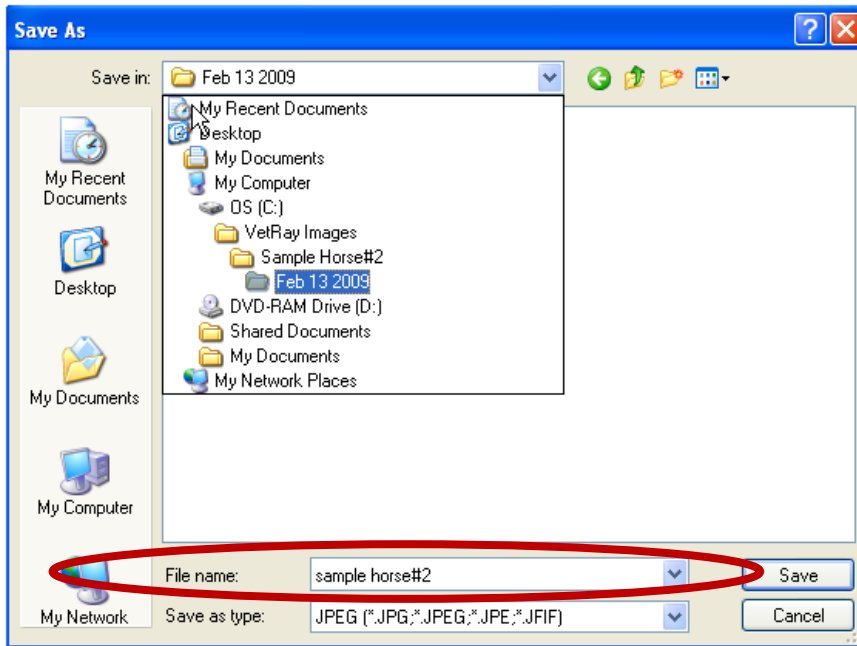

9. Click the “**Save**” button

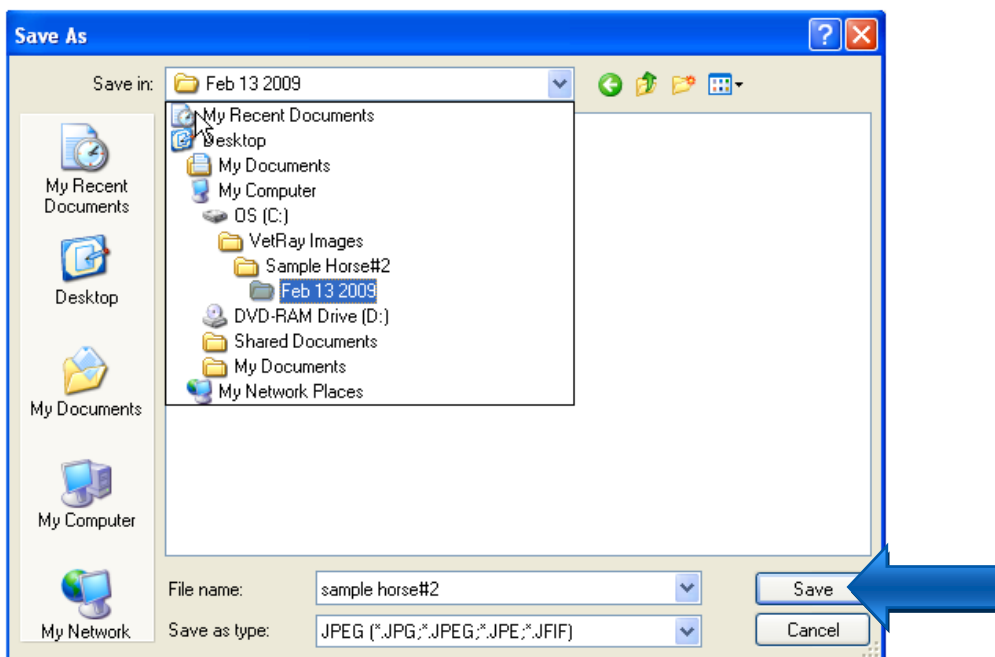

10. Open Metron

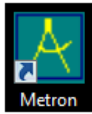

11. Select the Owner, Animal and Date that you will be import the images to.




12. Click the "File(s)..." Button




13. Select the folder that you saved your images in step 7. Our location is C:\VetRay Images\Sample Horse#2\ Feb 13 2009


14. You will then see the images that you exported from VetRay.

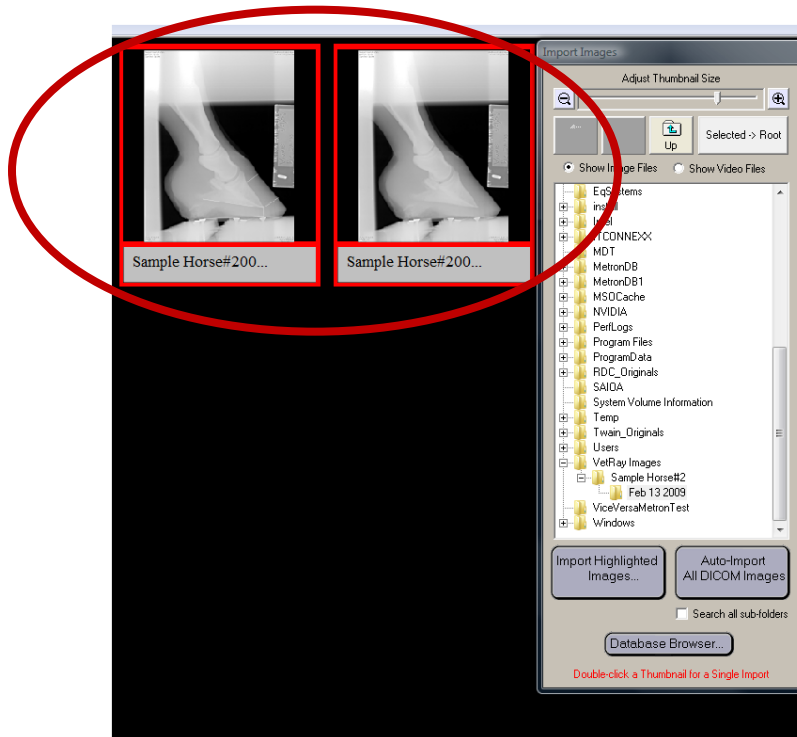

15. Select the images and then click the “Import Highlighted Images...” button.

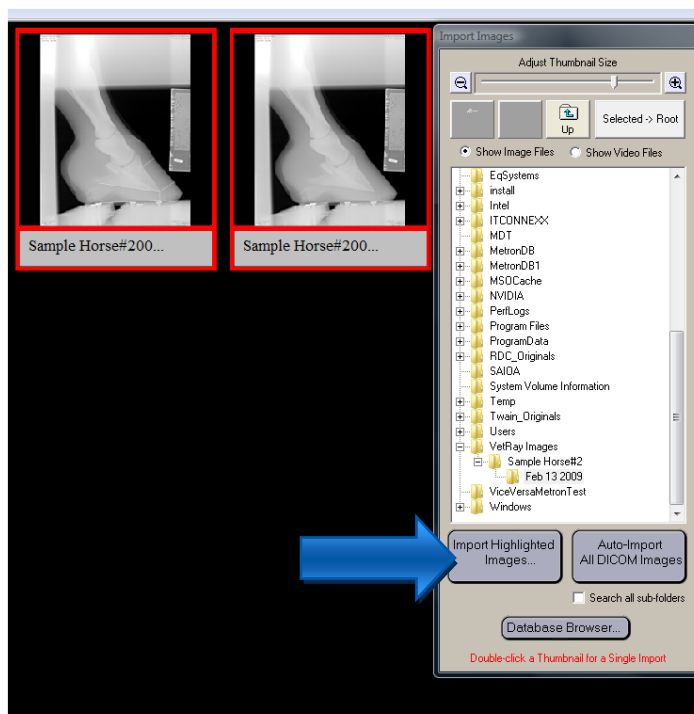


16. Verify that the Owner, Animal and Date fields are correct.

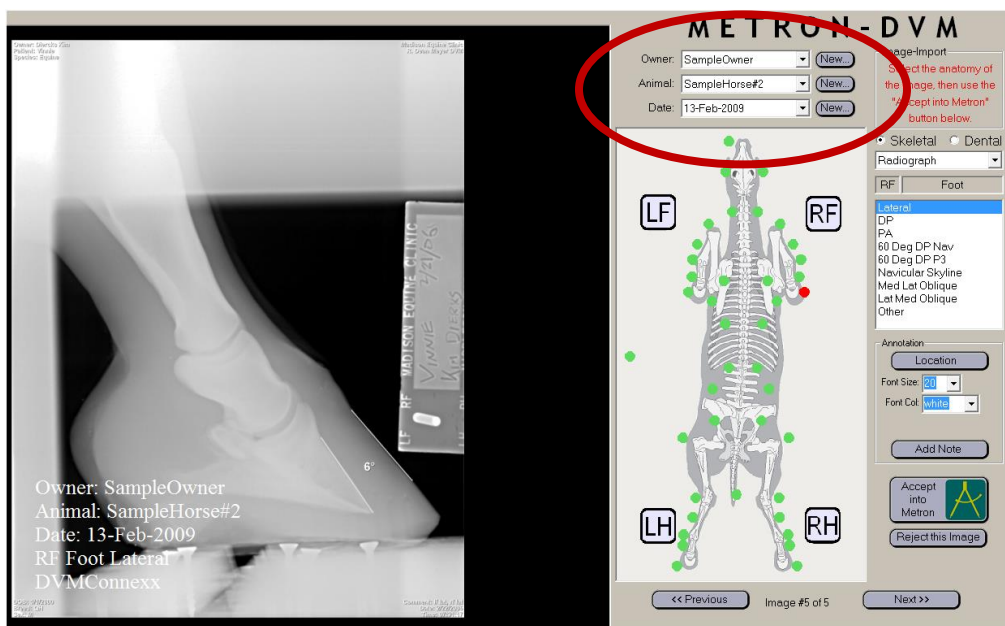


17. Select the **anatomy** of the image



18. Select the “Accept into Metron” button.

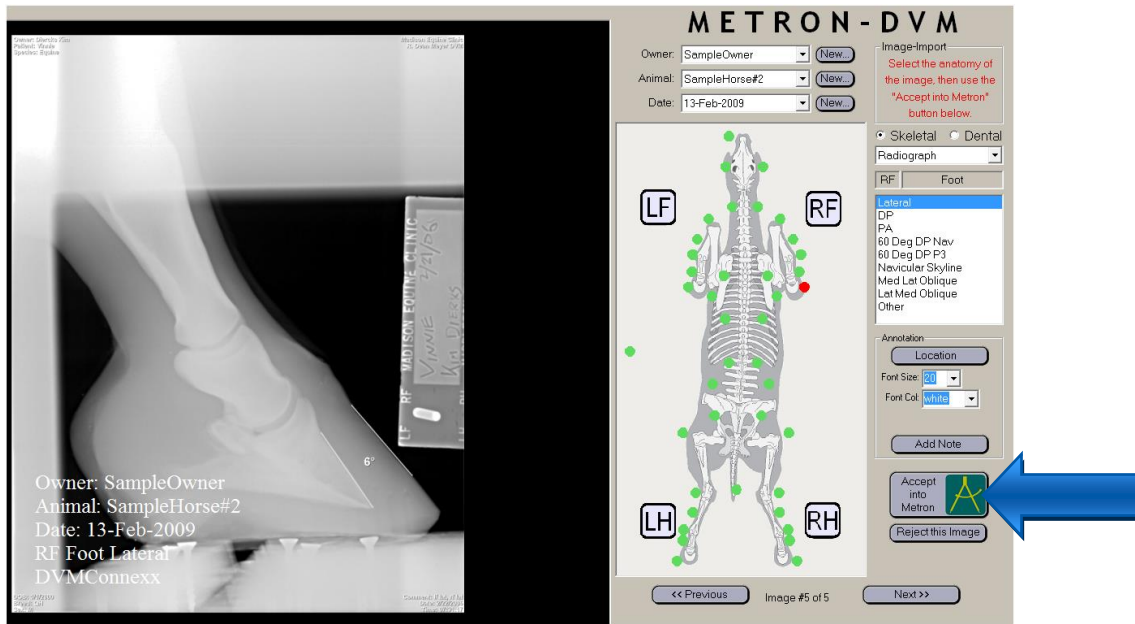

19. Repeat steps 16 through 18 for each image. In the upper right hand corner, you will notice the Multi-Import shows how many more images you will need to import.

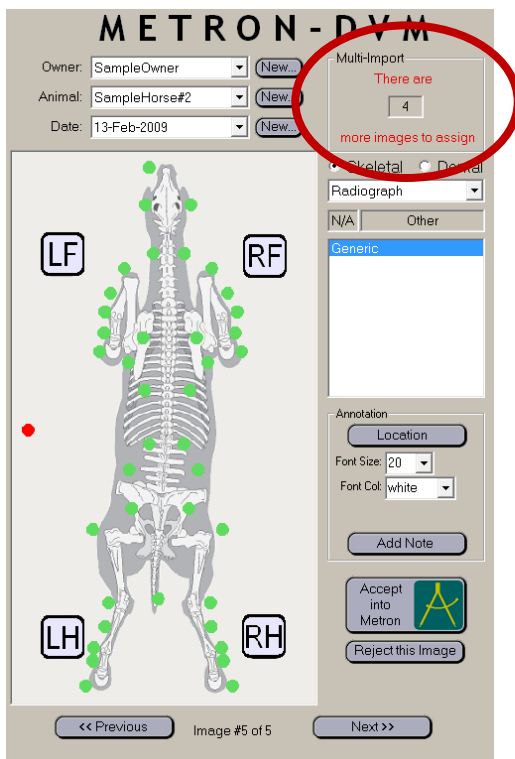



Contact us if you have additional questions.

We can be reached at metron@metronsupportservices.com or 1-877-638-3868.