



Metron Export to Android/Google Drive

Configuring Metron PC for Google Drive/Android Sync




877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com
METRON SUPPORT SERVICES. All rights reserved.

Metron Export to Android/Google Drive


Configuring Metron PC for Google Drive/Android Sync

YOUR METRON ADVANCED IMAGING SOFTWARE IS CONTINUALLY BEING ENHANCED. WE LISTEN TO YOU, STRIVING TO ADD AND ENHANCE THE PRODUCT AND SERVICE TO ASSIST YOU WITH THE DIAGNOSIS AND TREATMENT OF YOUR PATIENTS.

1. In order to setup Google Drive/Sync, you will need an existing Gmail/Google Account. If you don't already have a Gmail account, please visit <http://www.gmail.com> to set one up.
2. On the Metron Acquisition PC (the PC on which the radiographs are taken), open Internet Explorer and browse to <https://drive.google.com>
3. Click the “**Download**” link in the upper right-hand corner.
4. Click the “**Download Drive**” button so that it expands, and then click the “**Mac and PC**” option

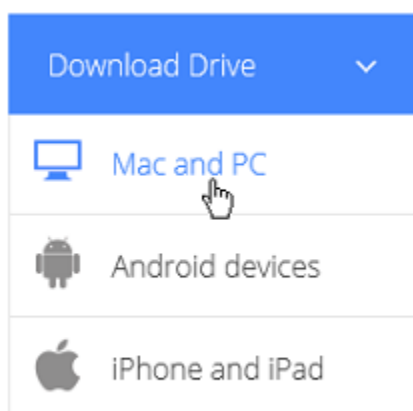

5. Click “Accept and Install”

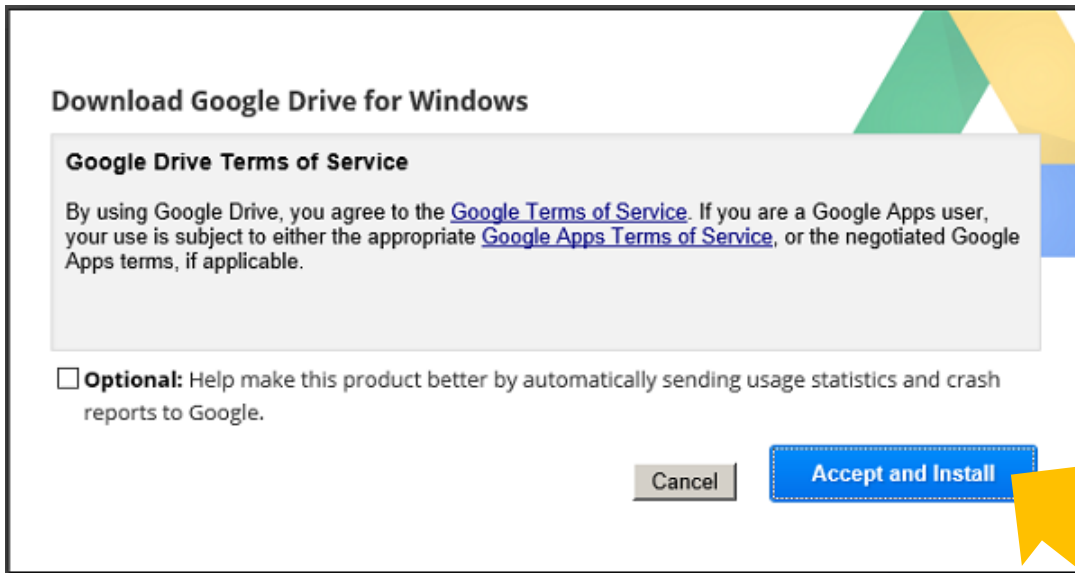


6. Click the “Get Started>” button




7. Enter the Gmail account information and click “Sign In”


8. Click “Next”

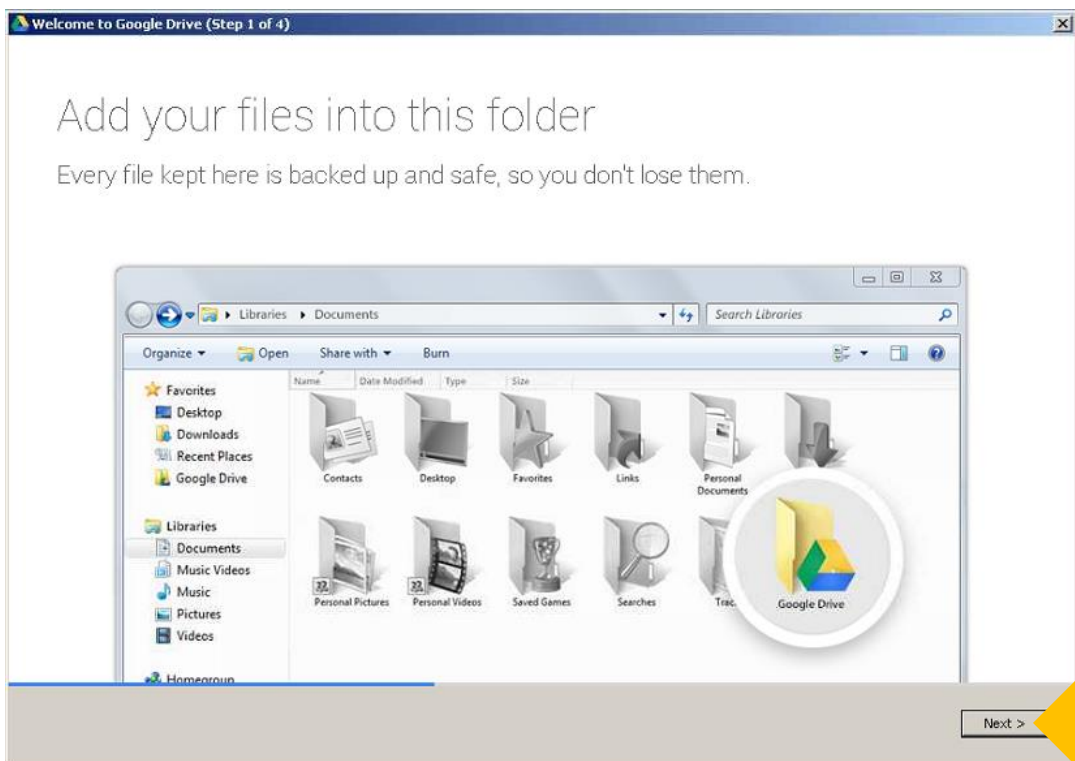


9. Click "Next"

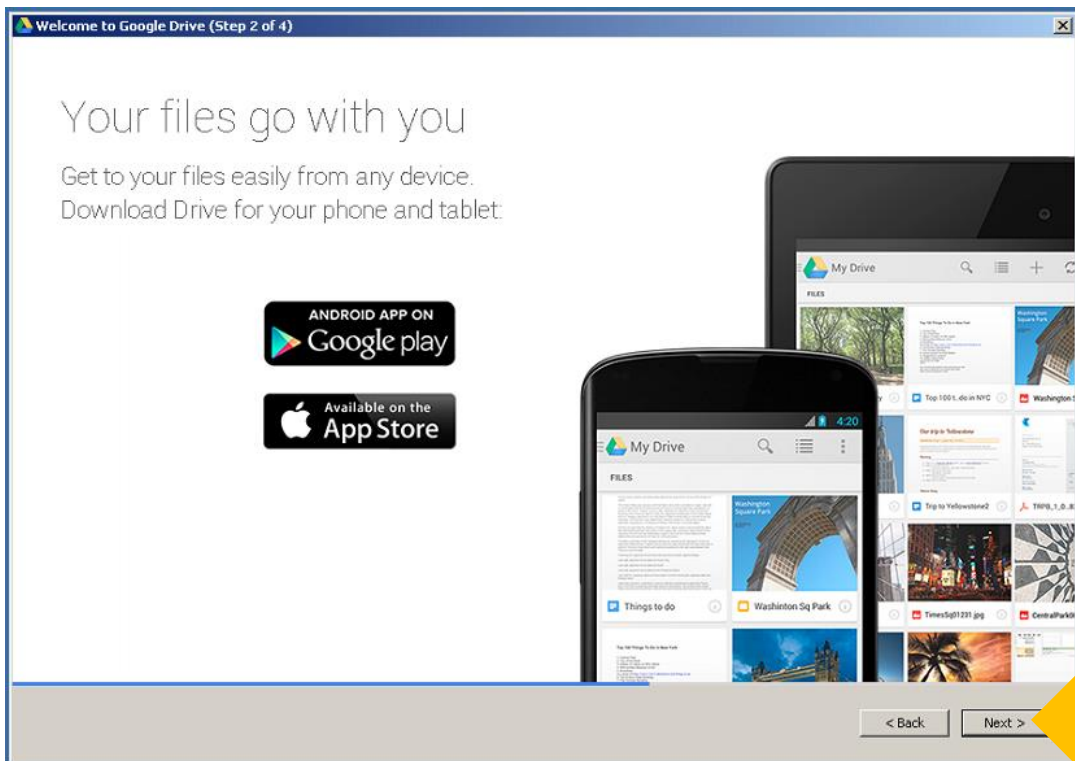


10. Click "Next"

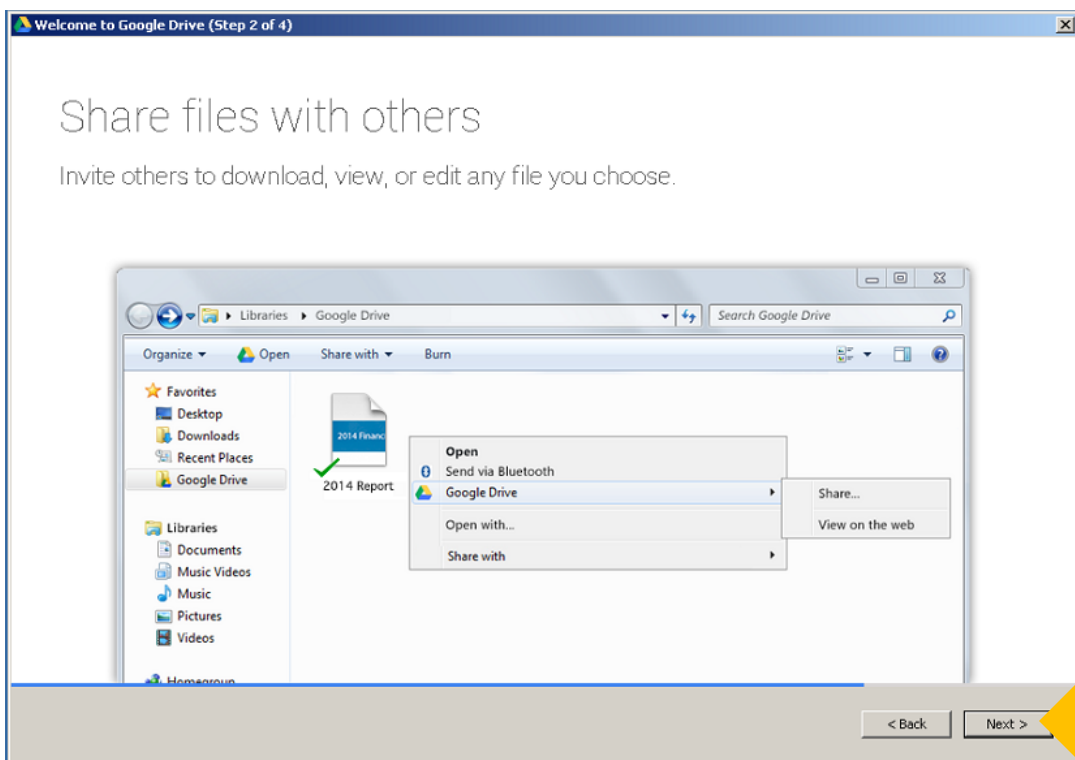


11. Click “**Advanced Setup**”

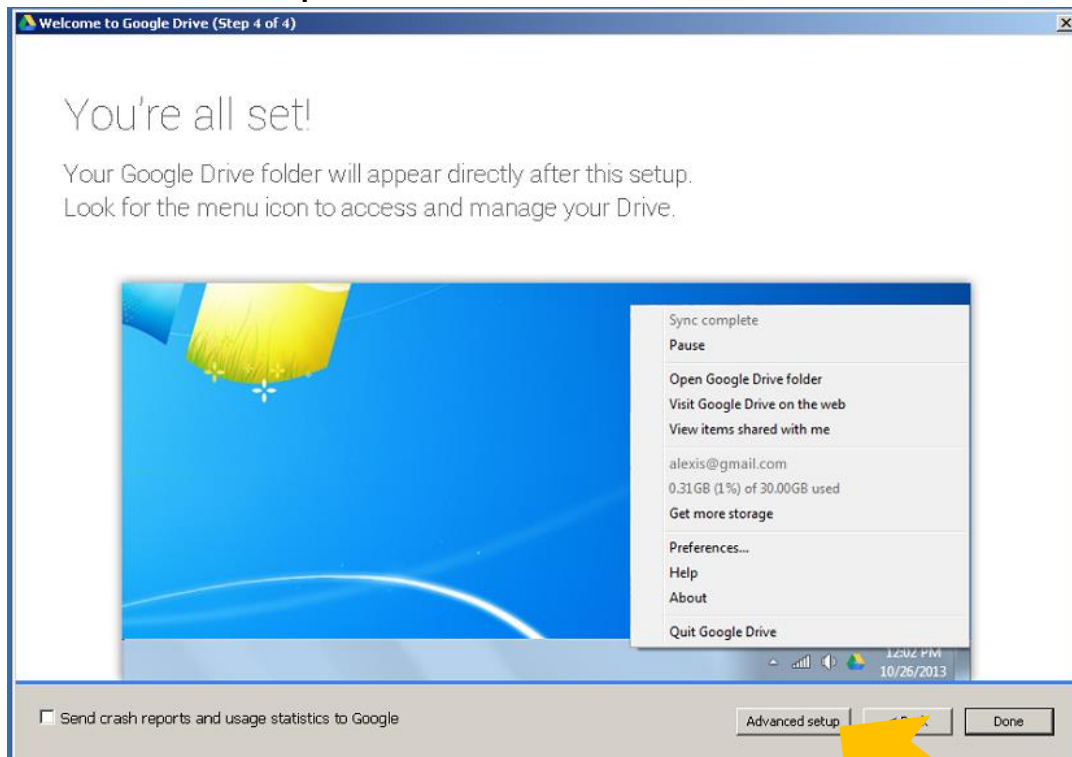


12. Click “**Change**” in the “**Folder location**” Section

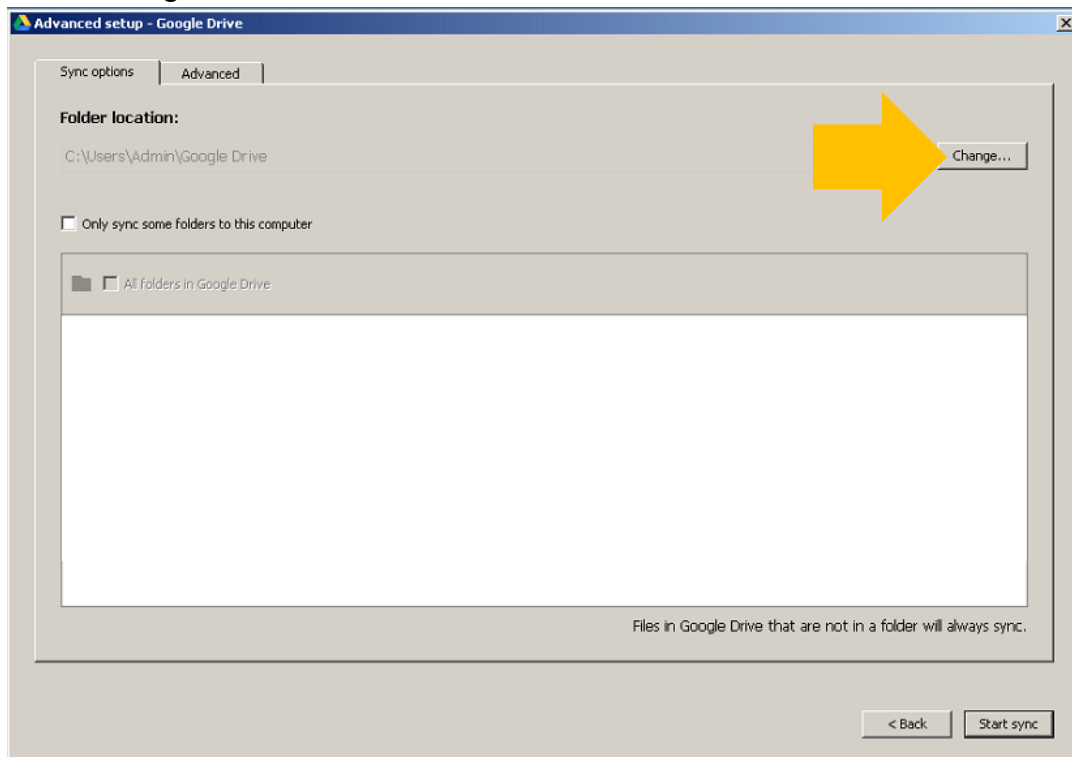

13. Expand “**Computer**” and then click the “**Local Dick (C:)**” Drive

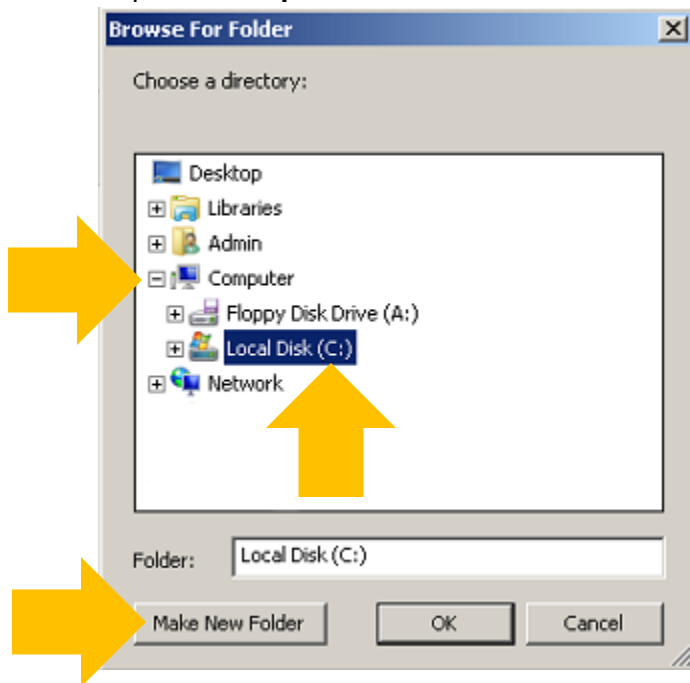


14. Click the “**Make New Folder**” button, name the new folder “**Google Drive**”, and then press the “**enter**” key on your keyboard.

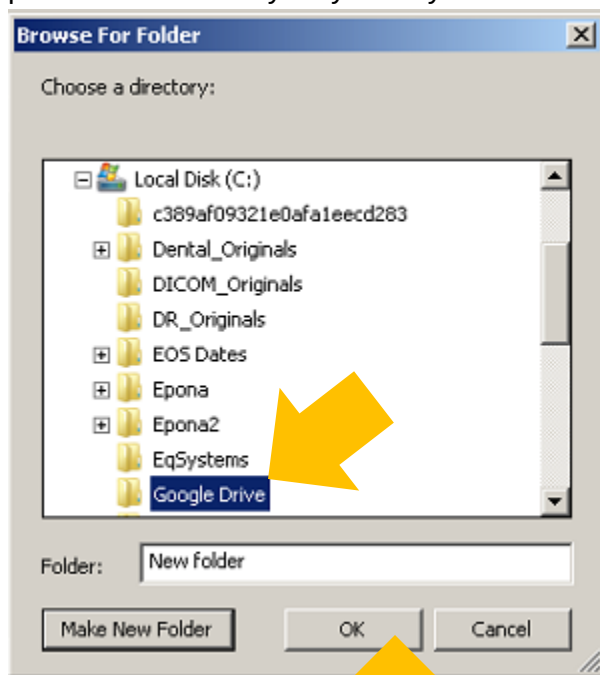

15. Click “**OK**”




16. Click **"Start Sync"**

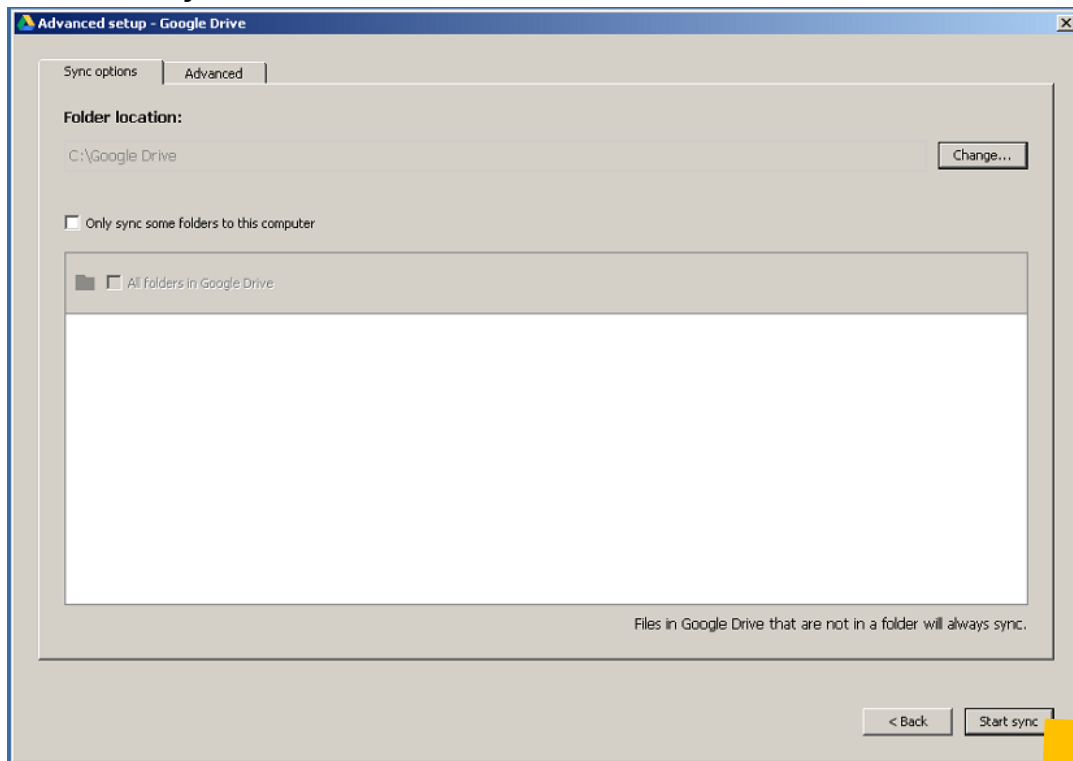


17. In Internet Explorer, browse to <https://drive.google.com> and login with your Gmail credentials.


18. Click the **"Create"** button and select **"Folder"**

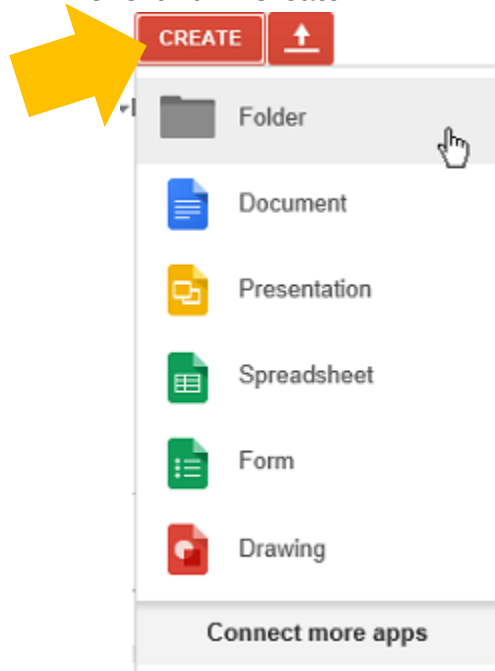


19. Enter the name “**MetronCloud**” and click “**Create**”

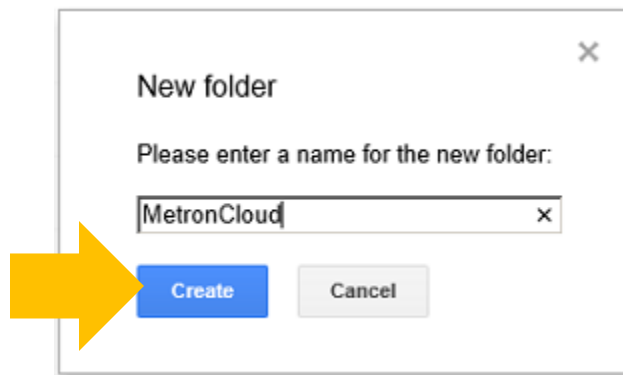

20. Right-click on the Google Drive icon in the system tray and select “**Preferences**”

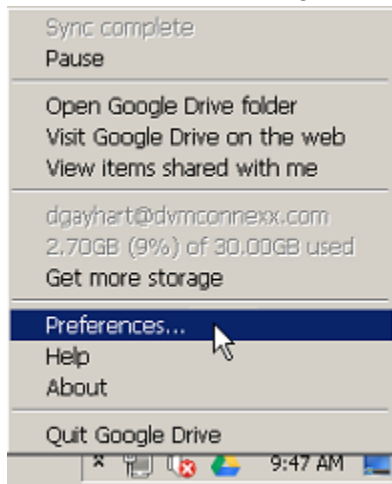

21. On the “**Sync Options**” tab, check the “**Only sync some folders to this computer**” option, and then check the “**MetronCloud**” folder. Click “**Apply Changes**”.

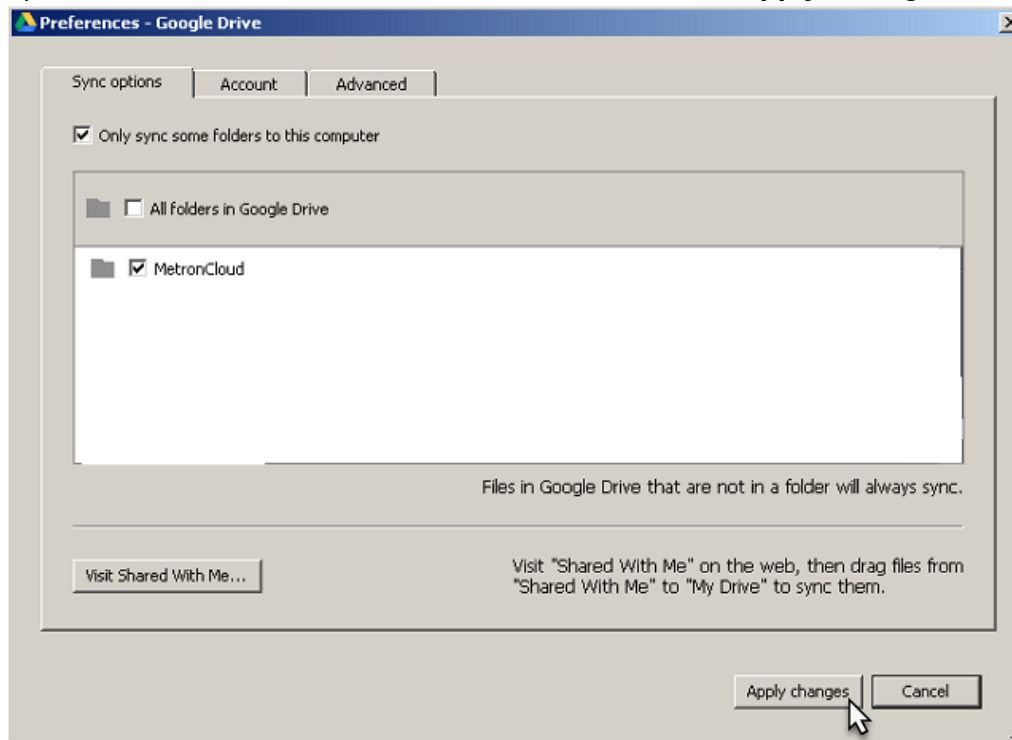


22. Open Metron, highlight any image, and then click the “**Export**” button



23. Click the “**Edit Rules**” button.



24. In the “**Default Folder**” field, enter C:\Google Drive\MetronCloud\<O>\<A>\<D>

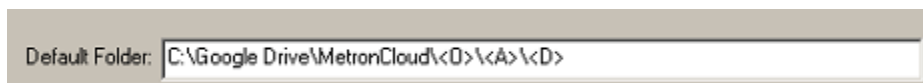

25. Verify the “**Default Format**” drop-down menu is set to “**JPG**”
26. In the “**Auto-Export**” section, verify that the options “**When Image Accepted**” and “**When DICOM Image received**” are checked.
27. Click the “**Close**” button.

Metron is now configured to auto-export .jpg images of all X-Rays Metron acquired to the Google Drive folder!

Installing Google Drive on Android Device

On many android-based tablets, Google Drive will come pre-installed. All you will need to do is open the Google Drive app and login with the user's credentials.

If not installed, below are directions to install it. **Note: Google Drive cannot be installed on Kindle Fire/Amazon-based tablets.**

1. Touch the **Google Play** icon on your phone or tablet.
2. Click on the magnifying glass and search for "**Google Drive**"
3. Locate Google Drive and install.
4. Press the "**Home**" button to return to the home screen on the tablet and, click the  Drive app icon.
5. The **MetronCloud folder** should be present. Perform a test acquisition (either an actual X-Ray or "Demo" X-Ray) in Metron to verify that data is synced successfully to the MetronCloud folder.



Contact us if you have additional questions.

We can be reached at metron@metronsupportservices.com or 1-877-638-3868.