

METR8N


Button Actions Reference Guide

Below is a guide to the Metron 8 buttons...their visual appearance and function.
Format simplified for use with touch screen laptops and tablets.

Database Browser Panel



Emergency: Enter directly to Acquire — give patient info later

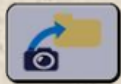

Import an image or video clip from a file



Click this button to add reports, forms or morphs to a study



View highlighted images side by side



Delete highlighted images



Export highlighted images



Print highlighted images



Send to Web-Viewer, iPhone/iPad, to DICOM server, or by email



Add highlighted items to the current collection




Stitch two images together to create one image



Cut images. You can then paste to a different date or animal study



Copy images. You can then paste to a different date or animal study

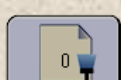

Paste images that you previously cut or copied



Create New Owner or New Animal



Acquire a clip or still frame from a USB video camera



Receive worklists and create studies



Pops On-Screen Keyboard
Touch screen only



Search/Edit Owner or Animal



View log of images received




Use Query Retrieve to pull images from a DICOM server


Single Image View Panel




Rotate image 180 degrees

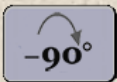

Rotate Clockwise

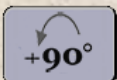

Rotate Counter Clockwise

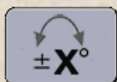


Customize Rotation



Mirror Image



Crop Image

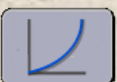

Apply Look Up Table (LUT)



Invert image



Window Level



Image Processing Tuner
Advanced Feature



Pseudo Color Emboss



Toggle in and out of Zoom Mode



Toggle in and out of Magnifier Mode




View DICOM Tags



View Images in Full Screen Mode

Return to Database

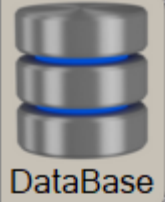

DataBase


Acquire Image



Acquire

(Optional)

Acquire w/Dental Sensor

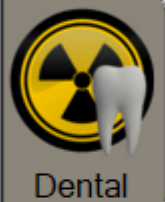

Dental

Edit a Single Image




Edit

System Folders



System

Exit Metron Program

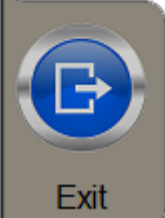

Exit

 METRON
 SUPPORT SERVICES