



Metron Web-Viewer, Android & More

View Metron image files on a Chromebook, Windows tablet, or ANY Android device



877-638-3868 | metron@metronsupportservices.com | www.metronsupportservices.com
METRON SUPPORT SERVICES. All rights reserved.

Metron Web-Viewer, Android & More

View Metron image files on a Chromebook, Windows tablet, or ANY Android device


The document contains detailed steps to send Metron Image files to the Metron Web-Viewer to display on your Chromebook, Windows tablet or any Android device (including your phone). Note: If you are using a phone to display images, it must have an Android OS System.

CHROMEBOOK USERS: If you are using a Chromebook, proceed to Page 3 for instructions.


WINDOWS USERS: If you are using a Windows tablet, proceed to Page 3 for instructions.

ANDROID DEVICE USERS: Proceed with Step 1 and download the Firefox® Browser.

1. Download Firefox® Browser for Android Version 33.0 or higher from Google Play. Click **"Install"**.

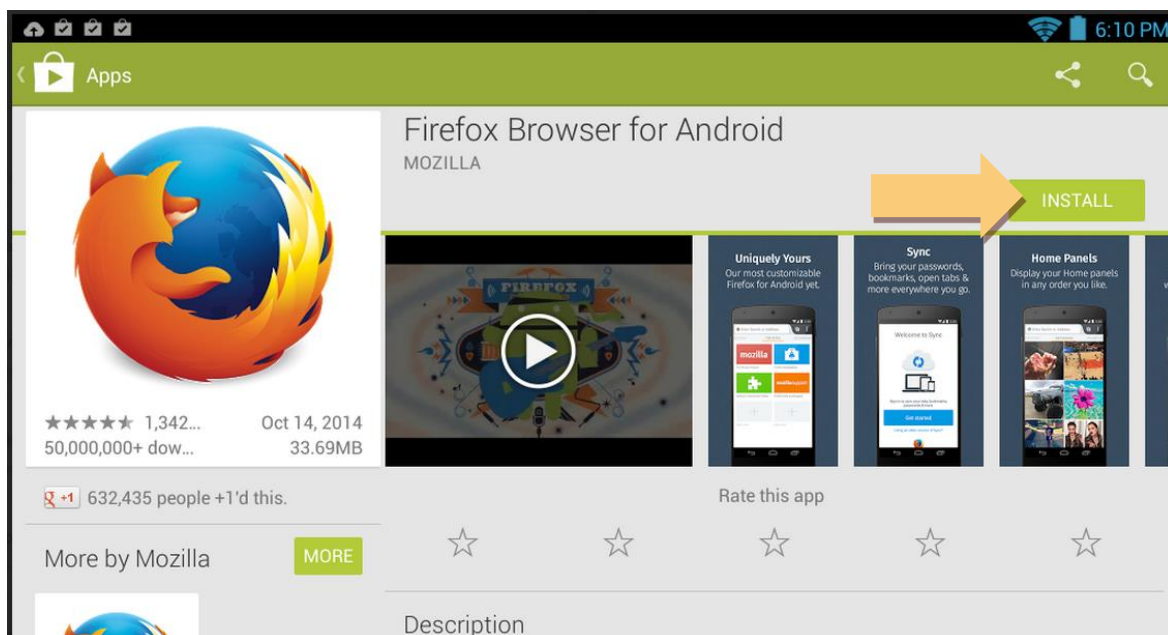

2. Click **"ACCEPT"**.

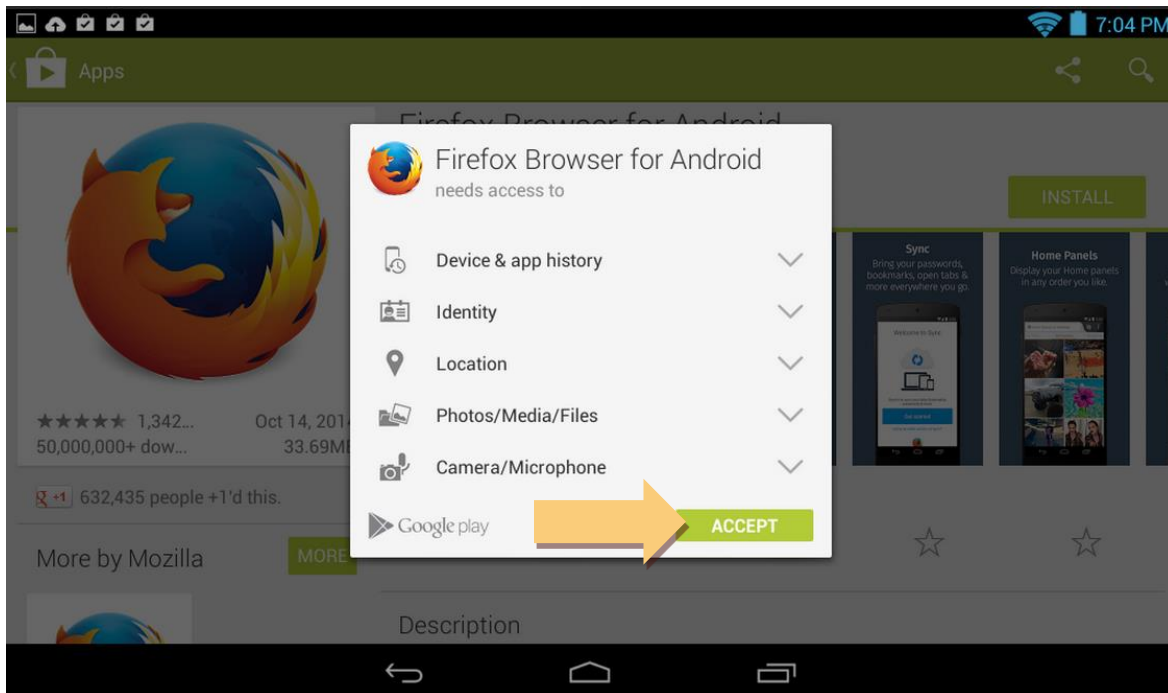


3. The installation will begin. Once completed, click **"OPEN"**.



4. You can also launch Firefox from your Android device.

ANDROID USERS: Open your Firefox Browser.  Go to <http://www.metronwebviewer.com/>


5. Tap anywhere on the screen:

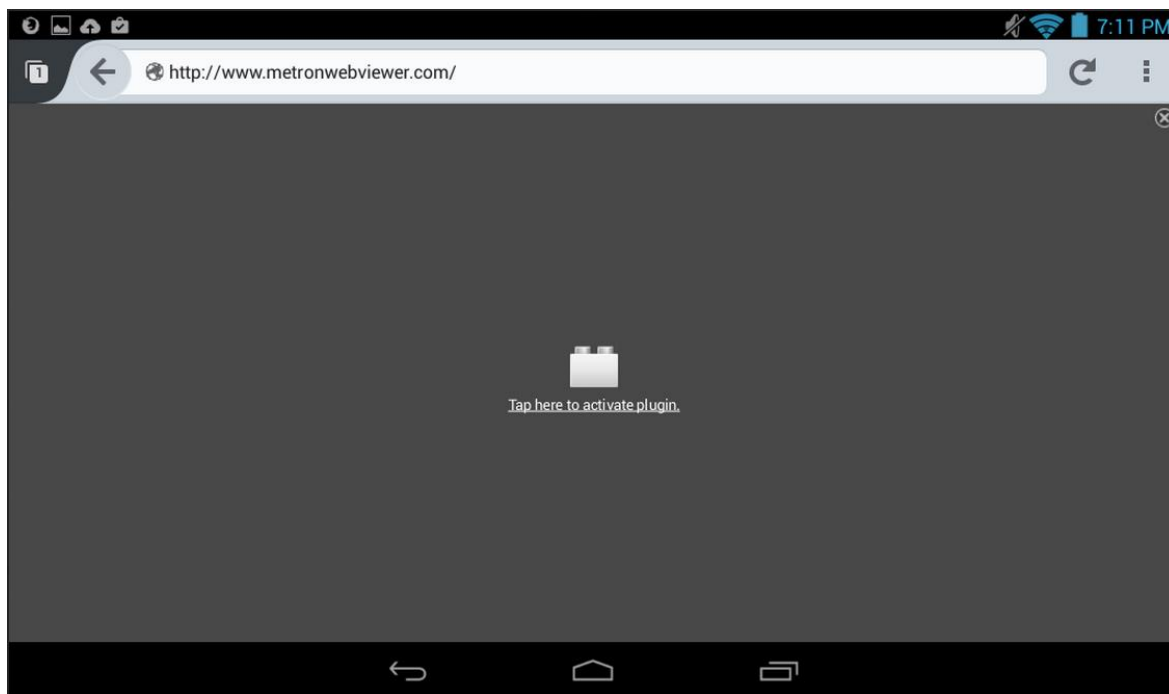


CHROMEBOOK USERS: Open your Chrome Browser.  Go to <http://www.metronwebviewer.com/>

WINDOWS USERS: Open Internet Explorer.  Go to <http://www.metronwebviewer.com/>


6. Enter the same ID for the Metron Web-Viewer.
This ID should match your ID in Metron. Check the box as shown to allow viewers to download images from this study. Click **“Send”**.
7. The User ID should also match Firefox on the Android device. **NO PASSWORD**.
Click **“Log in”**.

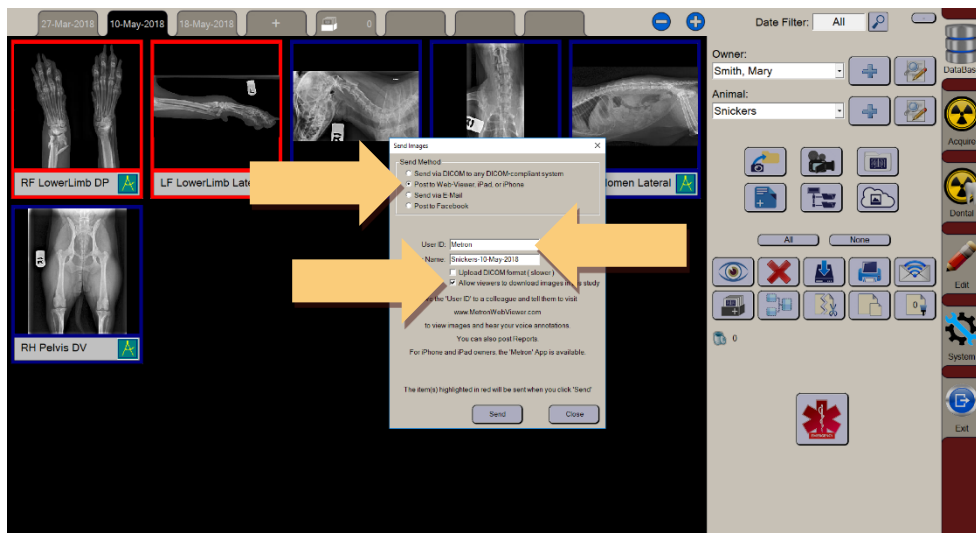

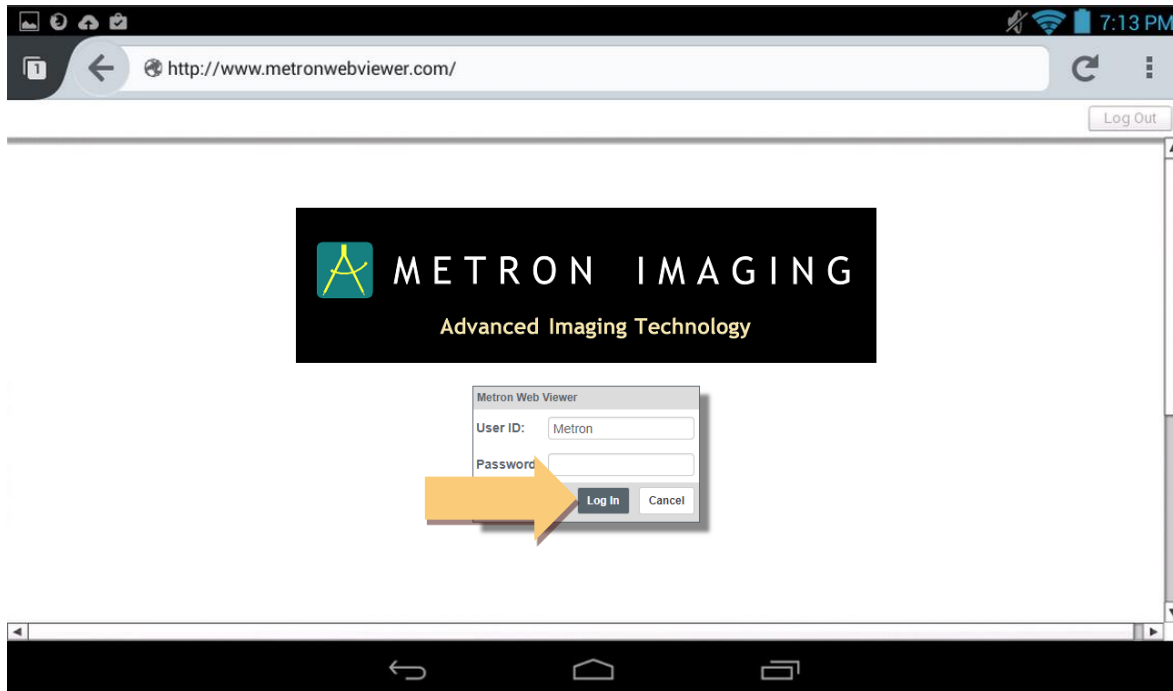

NOTE: If you mistyped the User ID in this screen, click “Cancel” so you can re-enter the User ID. The “Back Space key does not work.

8. Click “View”.

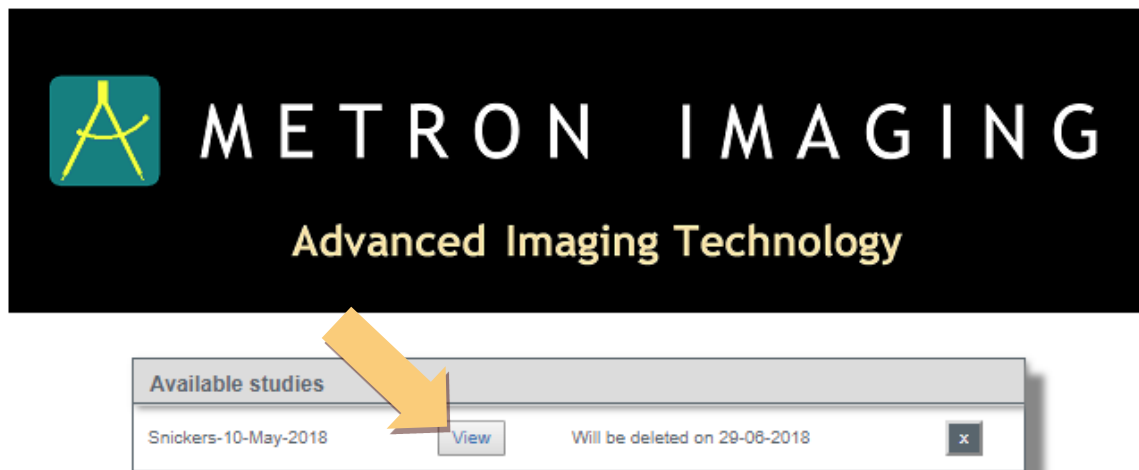

9. Your images are in the Study List:



10. You can also download your studies from the Metron Web-Viewer and save them to your Android device. Click the **“Download”** button.





Contact us if you have additional questions.
We can be reached at metron@metronsupportservices.com or 1-877-638-3868.