



Adding Images from a File into the Metron Database



Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.



Adding Images from a File into the Metron Database

This online training session will provide you with the detailed steps for
"Adding Images from a File into the Metron Database".



Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.



Adding Images from a File into the Metron Database

This online training session will provide you with the detailed steps for "Adding Images from a File into the Metron Database".

The Training Session Objective:

Provide information on how to add an image into the Metron Database from a file on your PC or your local network.

(Video Length: 5 minutes)




Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.

1

Adding Images from a File into the Metron Database

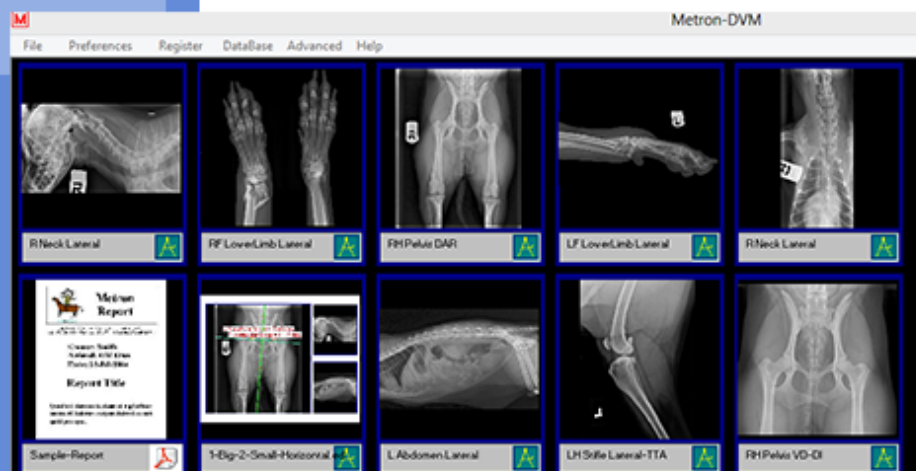


With Metron, you will always start in the 'Database Browser' screen where thumbnails of all the images in the currently selected study are displayed on the screen.



METRON
SUPPORT SERVICES

2

Adding Images from a File into the Metron Database

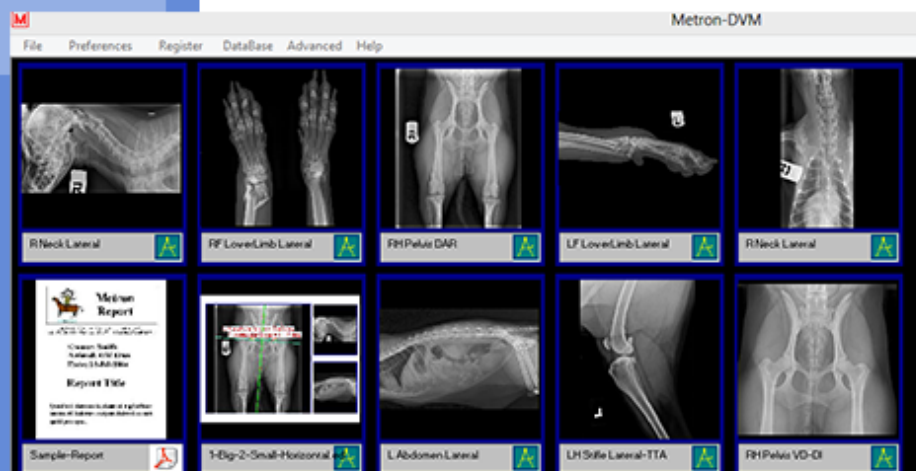

A study is the group of images for a particular owner, animal and date combination. Metron always starts at the study you were viewing when you last used Metron.



METRON
SUPPORT SERVICES

3

Adding Images from a File into the Metron Database



METRON
SUPPORT SERVICES

4

Adding Images from a File into the Metron Database



DATABASE BROWSER

Thumbnail Options... Default

Browse & Acquire | Collection (2 items)

Studies with new images: ☐ Day ☐ Week ☐ Month ☒ All

Owner: Fnh Miller

Animal: Cooper

Date: 9-Apr-2014

Add New Image, Video Clip, Form, or Report



METRON
SUPPORT SERVICES

5

Adding Images from a File into the Metron Database


Here we will present how to add an image into the Metron Database from a saved file on your PC. Any image file on your PC or on your local network can be imported into Metron.



6


Adding Images from a File into the Metron Database

Here we will present how to add an image into the Metron Database from a saved file on your PC. Any image file on your PC or on your local network can be imported into Metron.

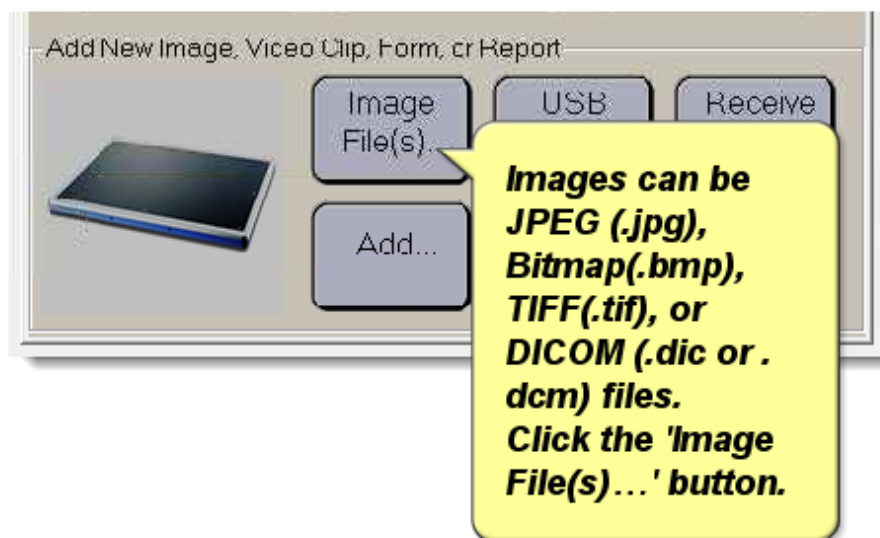

7


Adding Images from a File into the Metron Database

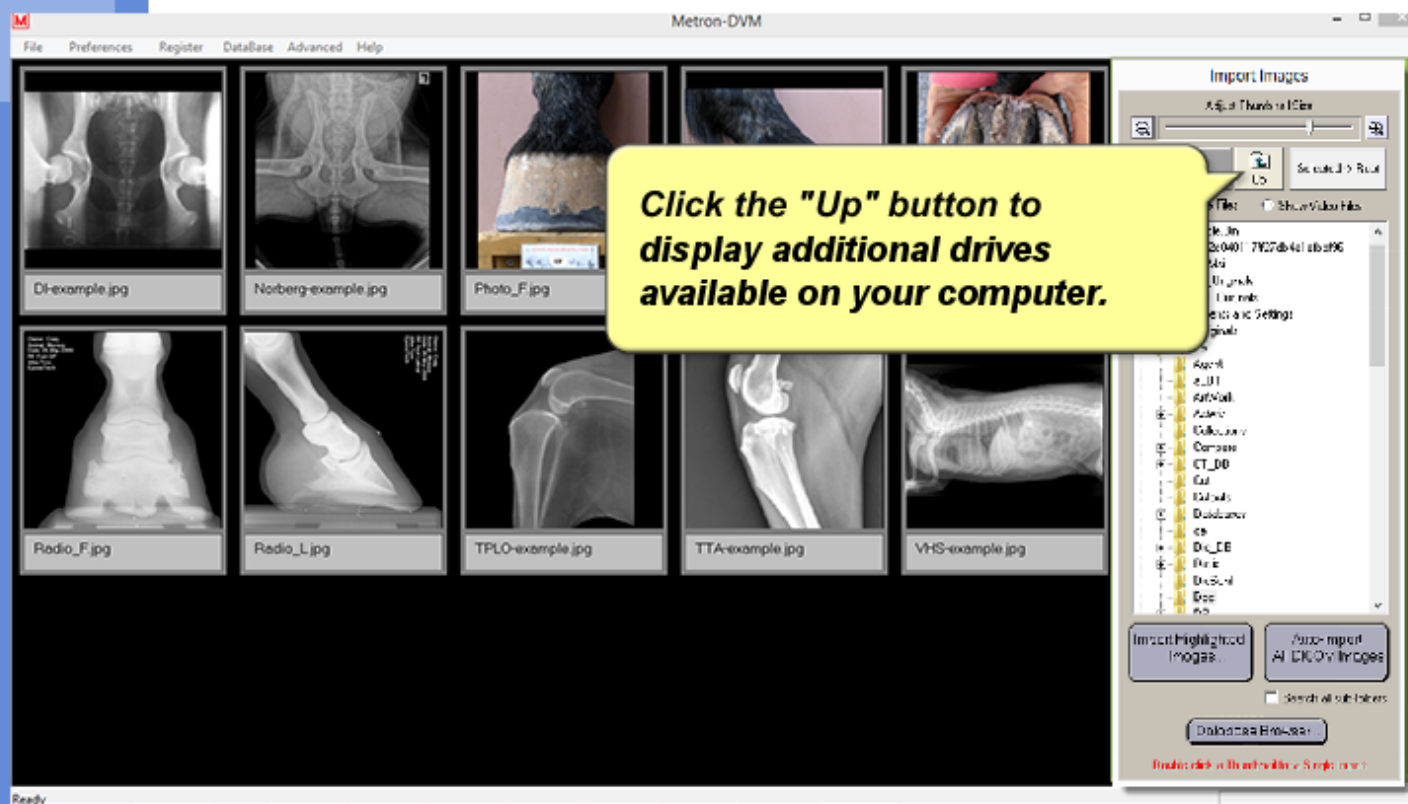

METRON
SUPPORT SERVICES


8

Adding Images from a File into the Metron Database

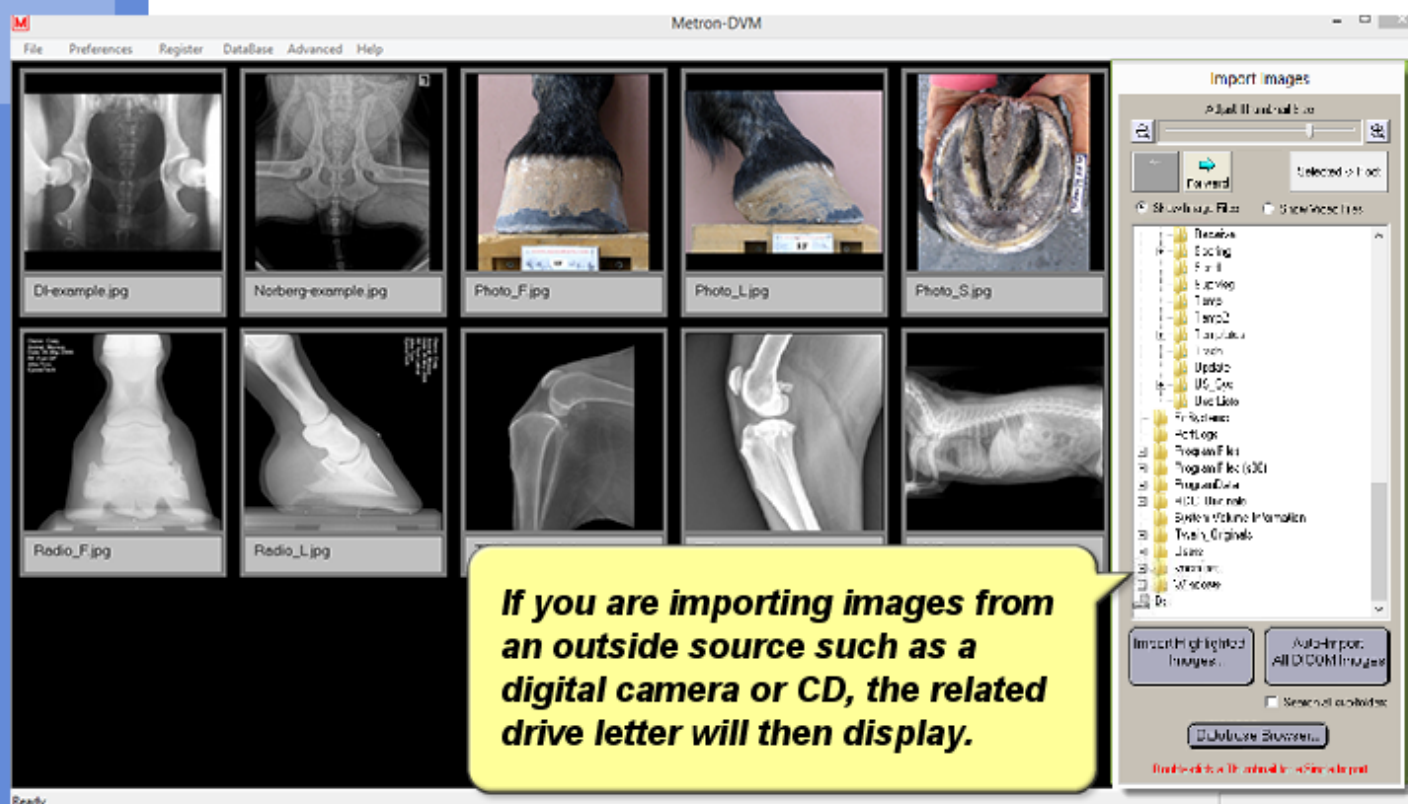

9


Adding Images from a File into the Metron Database

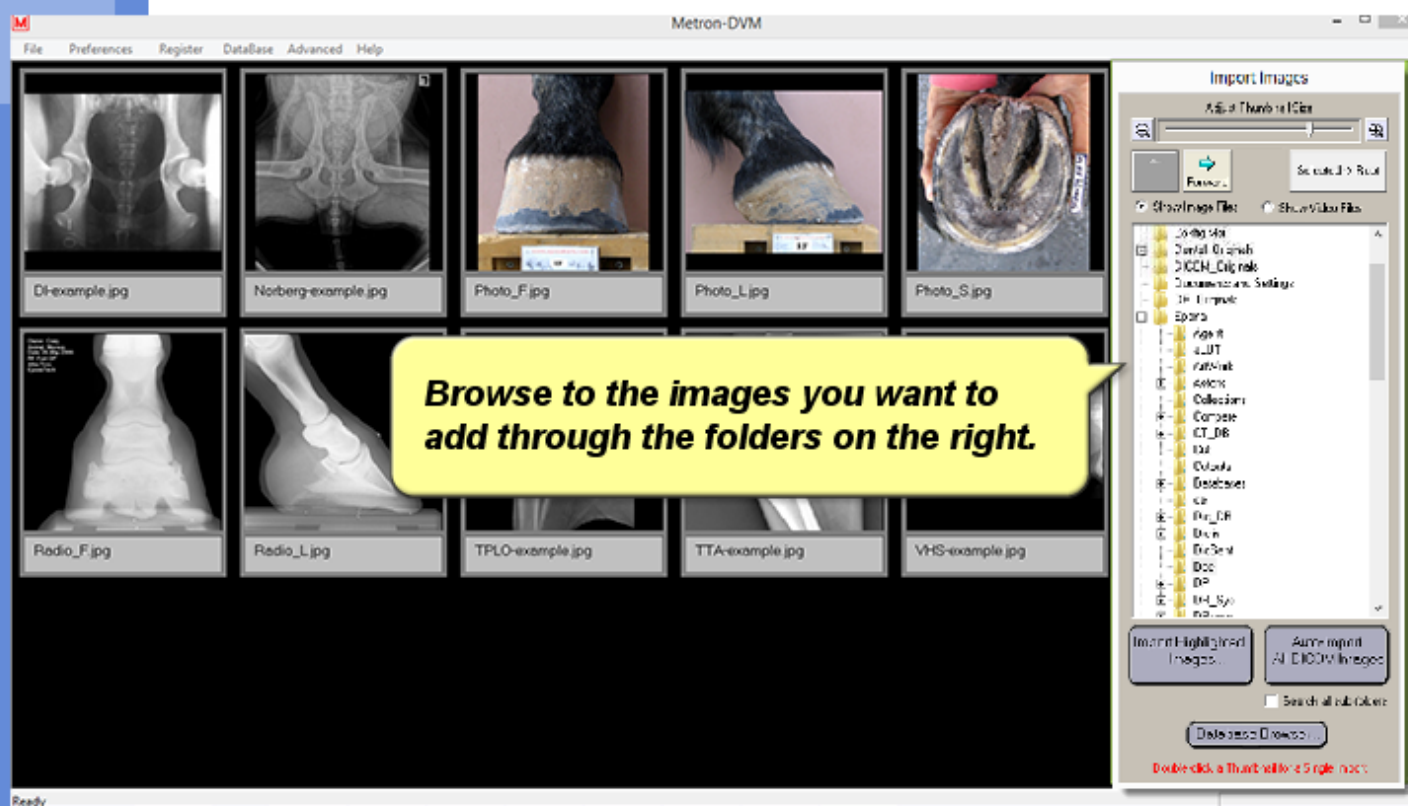

METRON
SUPPORT SERVICES

10

Adding Images from a File into the Metron Database

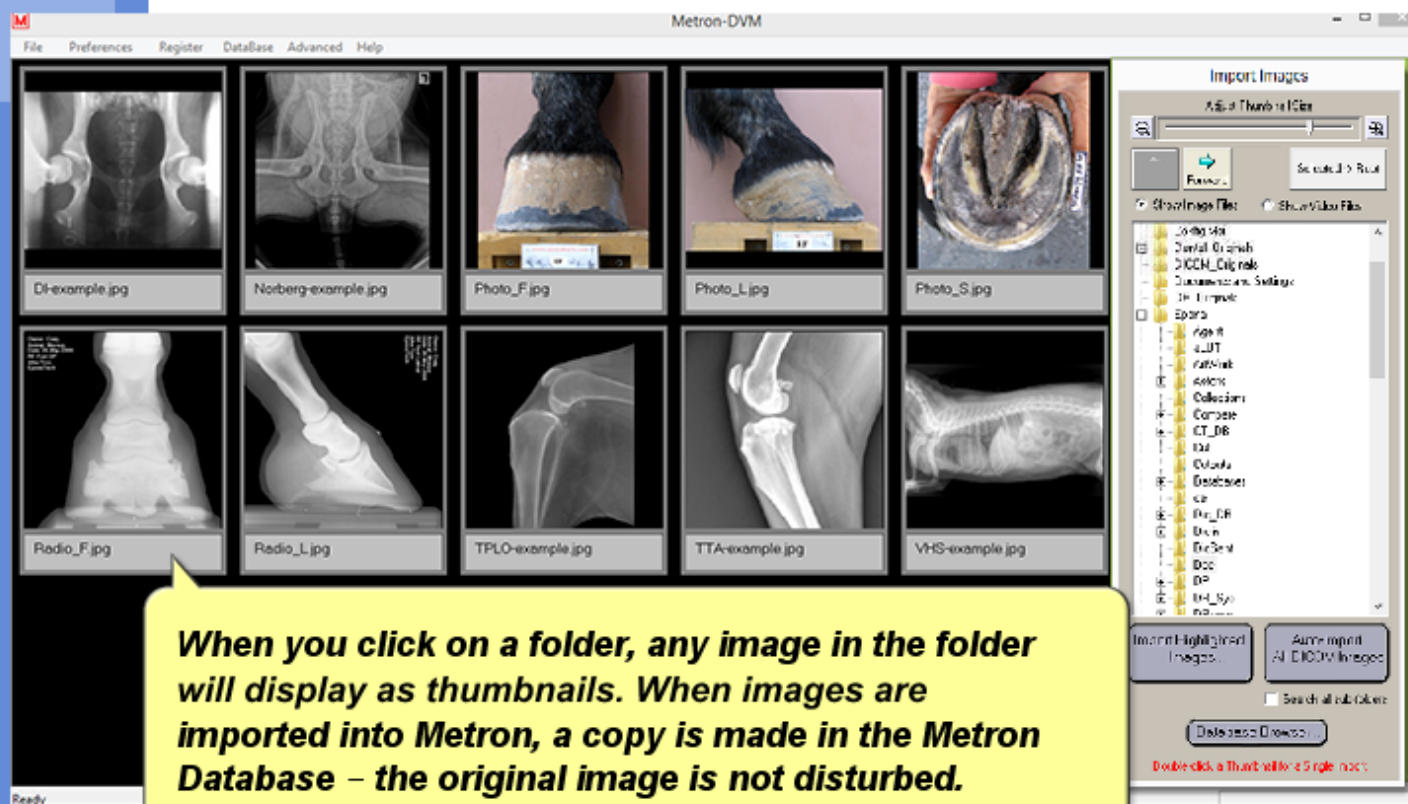

METRON
SUPPORT SERVICES


Adding Images from a File into the Metron Database




METRON
SUPPORT SERVICES

12

Adding Images from a File into the Metron Database

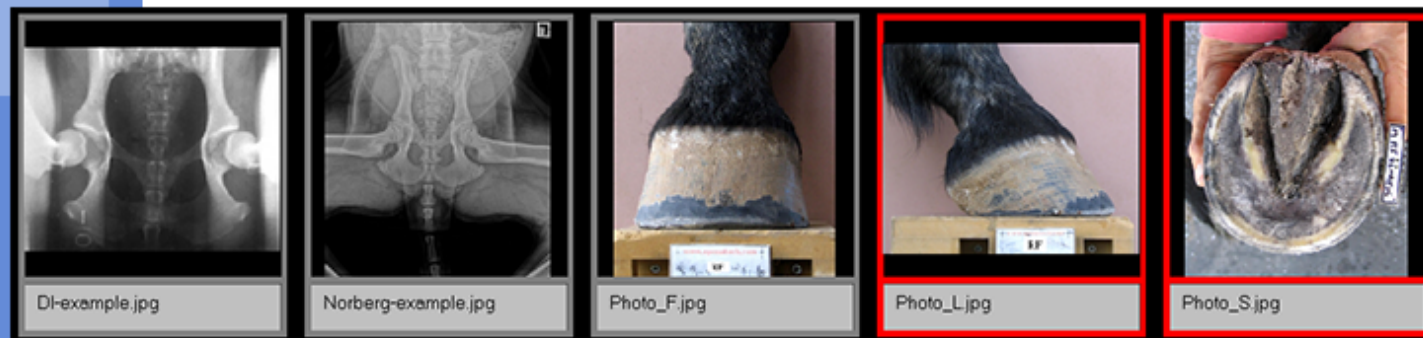


Clicking on multiple images and selecting the 'Import Highlighted Images' button will bring the images to the Single Image Viewer.

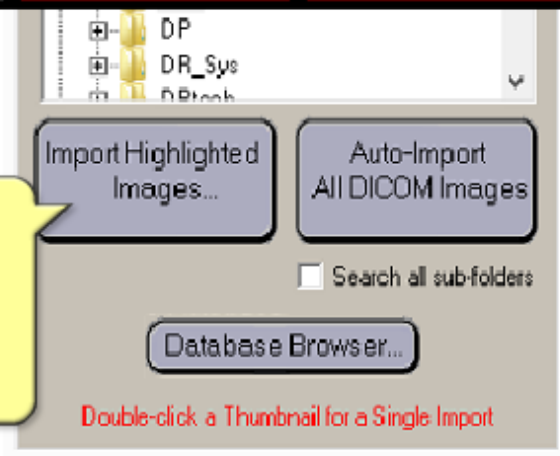


METRON
SUPPORT SERVICES

13

Adding Images from a File into the Metron Database

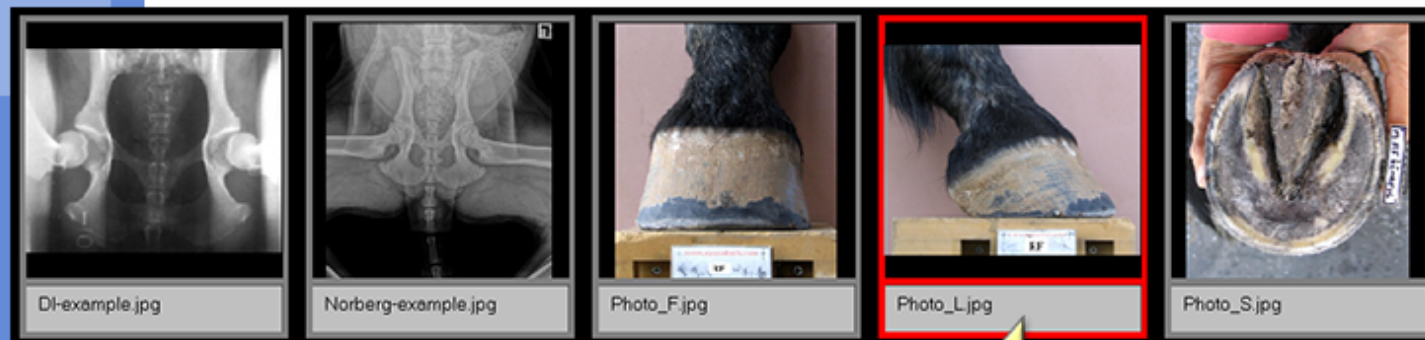

If you would like to import only one image, you can double click the image thumbnail.

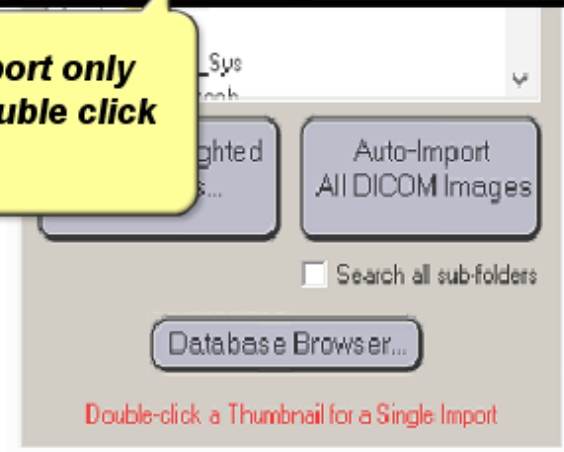

METRON
SUPPORT SERVICES

14

Adding Images from a File into the Metron Database

At the Single Image Viewer, select an owner or add new owner, animal and date.

This will default to the current study you have open.

METRON
SUPPORT SERVICES

Adding Images from a File into the Metron Database

METRON-DVM

Owner:
Animal:
Date:

Image Input
Select the anatomy of the image, then use the "Accept into Database" button below.

☒ Skeletal ☐ Dental


Photograph:

Select View:
Lateral
DP
PA
Sider
Blade
Medial
Other


Annotation
Location:
Sub: White:
☒ B/L ☒ Region ☒ View

Next select the anatomy region, view, and add any annotation changes needed for the specific image.

Adding Images from a File into the Metron Database

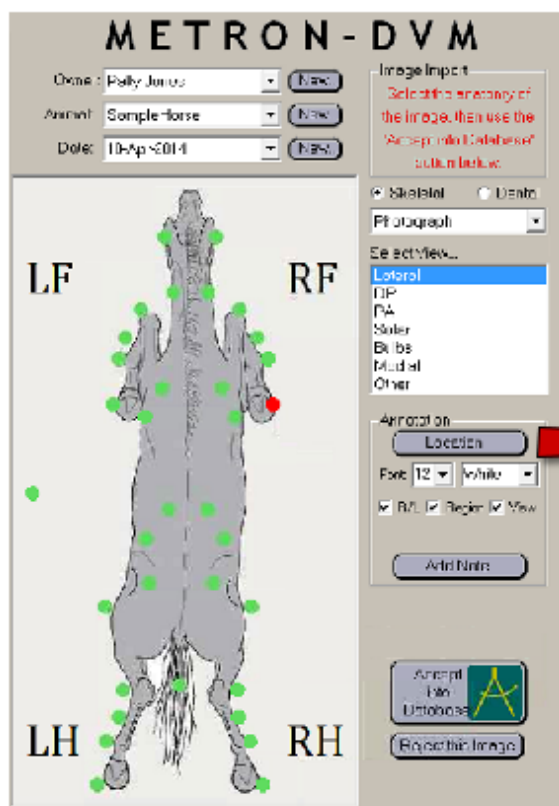

Click on one of four icons in the **Anno Location** panel to place the annotation in the corresponding corner of the image ... or select **'Select Automatically'** to let Metron choose the best corner.

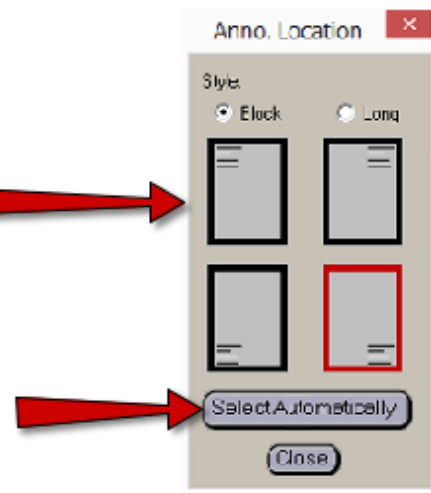

Adding Images from a File into the Metron Database



Adding Images from a File into the Metron Database

METRON-DVM

Owner: Polly Jones

Animal: Sample Horse

Date: 10-Apr-2014

Image Input
 Select the location of the image, then use the 'Accept into Database' button below.

☒ Skeletal ☐ Dental

Photograph:

Select View:
 Lateral
 DP
 PA
 Saker
 Bilbo
 Medial
 Other

Annotation:

 Font: 12 White
☒ B/L ☒ Region ☒ View

When completed, click **'Accept into Database'**. If you have selected more than one item, the next item will automatically display after you click the **'Accept into Database'** button.



Adding Images from a File into the Metron Database

If you decide you do not want to add this image, click 'Reject this Image'.

When completed, click **'Accept into Database'**. If you have selected more than one item, the next item will automatically display after you click the **'Accept into Database'** button.

Adding Images from a File into the Metron Database

METRON-DVM

Owner: Polly Jones

Animal: Sample Horse

Date: 10-Apr-2014

LF RF

LH RH

MultiInput

There are 2 more images to accept

☐ Skeletal ☐ Dental

Photograph

Lateral

DP

PA

Solar

Bite

Mandib

Other

Annotation

Font: 12 White

☒ B/L ☒ Register ☒ View

If you originally selected multiple images, you will be shown each image one at a time for labeling and accepting.

Adding Images from a File into the Metron Database

Metron will remind you how many images you have remaining to label and accept.

The screenshot shows the METRON-DVM software interface. At the top, there are dropdown menus for Owner (Polly Jones), Animal (Sample Horse), and Date (10-Apr-2014), each with a 'New' button. Below these is a diagram of a horse from a dorsal view, with anatomical labels: LF (Left Foreleg), RF (Right Foreleg), LH (Left Hindleg), and RH (Right Hindleg). Green dots are placed on the horse's body to indicate points of interest. To the right of the diagram is a sidebar with various controls: a 'Multi-point' section showing 'There are 2' points and a 'more images to assign' link; a 'Select View...' dropdown menu with options: Lateral, DP, PA, Saker, Blbs, Medial, and Other; an 'Annotation' section with a 'Location' button, a 'Font' dropdown (set to 12), a 'Color' dropdown (set to White), and checkboxes for 'B/L', 'Register', and 'View'; an 'Add Point' button; and a bottom section with 'Accept', 'Use Database', and 'Reject this Image' buttons.

In the following six slides, you will be presented questions to evaluate your knowledge of **"Adding Images from a File into the Metron Database"**.

Let's get started!



METRON
SUPPORT SERVICES

23

Adding Images from a File into the Metron Database

Images may be JPEG, Bitmap, TIFF, or DICOM files. Browse to the images you want to add through the folders on the right of the screen. When you click on a folder, any images in the folder will display as thumbnails:

- ☐ When images are imported into Metron, a copy is made on your PC - the original image is not disturbed.
- ☒ When images are imported into Metron, a copy is made in the Metron database - the original image is not disturbed.
- ☐ When images are imported into Metron, a copy is made in the Metron database - the original image is erased.

Submit

Congratulations, your answer is correct!

Sorry, your answer is incorrect. Please try again.

METRON
SUPPORT SERVICES

Adding Images from a File into the Metron Database

Clicking on multiple images and selecting the “Import Highlighted Images” button will bring the images to the Single Image Viewer. If you would like to import only one image:

- ☒ Double click the image thumbnail.
- ☐ Single click the image thumbnail.
- ☐ Double click “Database Browser”.

Submit

Congratulations, your answer is correct!

Sorry, your answer is incorrect. Please try again.

25

Adding Images from a File into the Metron Database

At the Single Image Viewer, select an owner or add new owner, animal and date. This will default to the current study you have open. Next select the:

- ☐ View and Database browser.
- ☐ Database browser and location.
- ☒ Anatomy region, view, and add any annotation changes for the specific image.

Submit

Congratulations, your answer is correct!

Sorry, your answer is incorrect. Please try again.

METRON
SUPPORT SERVICES

26

Adding Images from a File into the Metron Database

When done, click “Accept into Database”. If you have selected more than one item:

- ☒ The next item will automatically display after you click the “Accept into Metron” button.
- ☐ The next item will not display.
- ☐ The next item will not display until you click on “Reject this Image”.

Submit

Congratulations, your answer is correct!

Sorry, your answer is incorrect. Please try again.

METRON
SUPPORT SERVICES

Adding Images from a File into the Metron Database

If you decide you do not want to add this image:

- ☐ Click "Accept into Database".
- ☐ Click "Add Notes" and type in reject.
- ☒ Click "Reject this Image."

Submit

Congratulations, your answer is correct!

Sorry, your answer is incorrect. Please try again.

Adding Images from a File into the Metron Database

If you originally selected multiple images, you will be shown:

- ☒ Each image one at a time for labeling and accepting.
- ☐ All images at the same time for accepting.
- ☐ All images at the same time for rejecting.

Submit

Congratulations, your answer is correct!

Sorry, your answer is incorrect. Please try again.



Adding Images from a File into the Metron Database

This completes the "Adding Images from a File into the Metron Database" online training. Please feel free to view it again or select another online training session to review.



Click F11 on your keyboard to minimize the browser screen
and click the window closed to exit this training session.

