



Dental Positioning Guides for Canine and Feline



Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.



Dental Positioning Guides for Canine and Feline

This online training session will provide you with the details of how to use the "Dental Positioning Guides for Canine and Feline".



Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.


1

Dental Positioning Guides for Canine and Feline


Dental Positioning Guides provide the ability to understand the ideal location for the X-Ray machine and plate to help produce the best image. The guides show the animal, the film/plate and the X-Ray source.

As you see in the images, the animal position is sometimes shown on its side or sometimes on its back and is able to be rotated.


Maxillary Molars

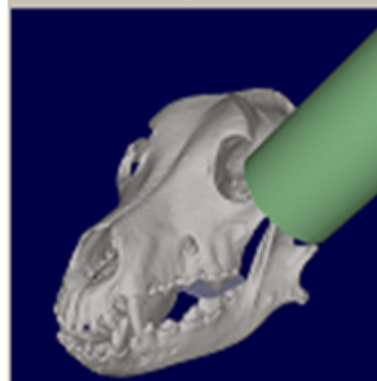

Maxillary Molars

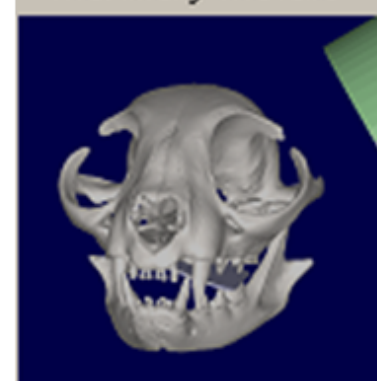

Mandibular Molars

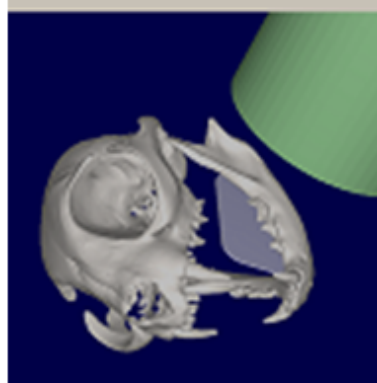

METRON
SUPPORT SERVICES


2

Dental Positioning Guides for Canine and Feline


Dental Positioning Guides provide the ability to understand the ideal location for the X-Ray machine and plate to help produce the best image. The guides show the animal, the film/plate and the X-Ray source.

As you see in the images, the animal position is sometimes shown on its side or sometimes on its back and is able to be rotated.

Mandibular Molars

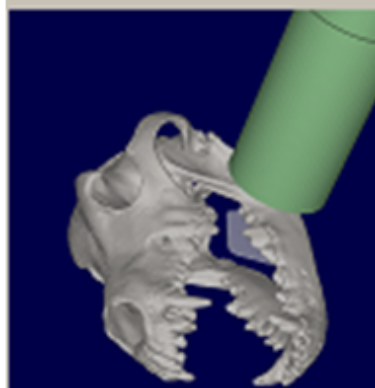


Maxillary Incisors

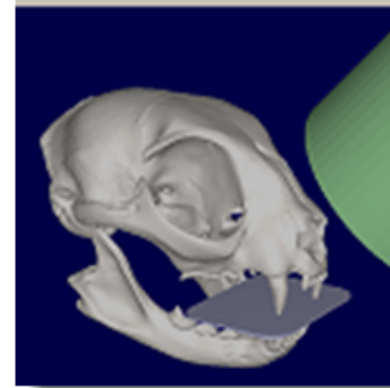


The dental positioning guides are available in Metron for canine and feline. These guides will soon be available for equine.

3

Dental Positioning Guides for Canine and Feline

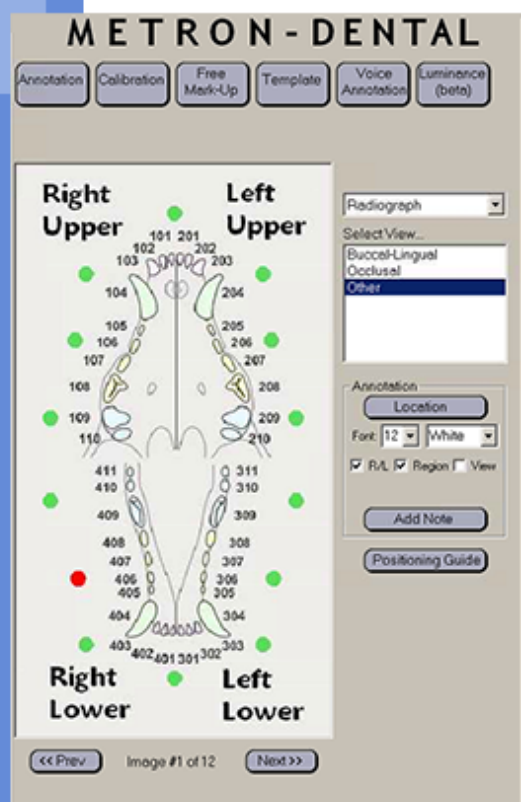


The Positioning Guides are available from the Single Image View mode for canine or feline when acquiring a new image or viewing a current image.

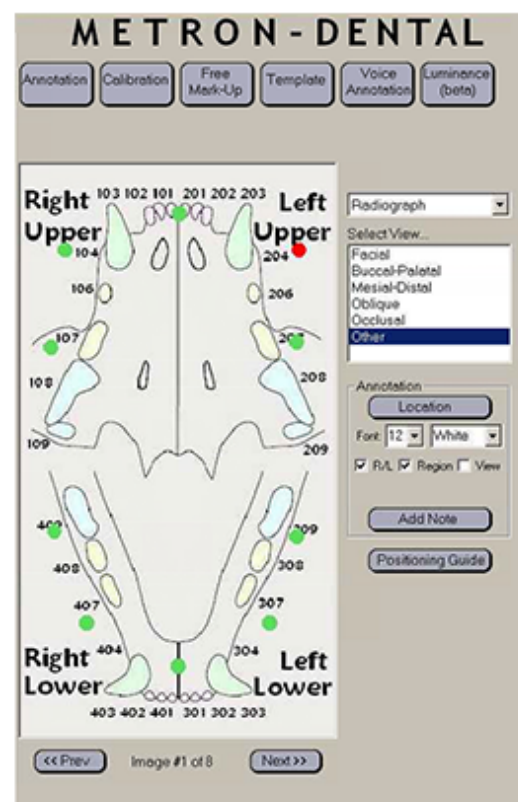


METRON
SUPPORT SERVICES

Dental Positioning Guides for Canine and Feline

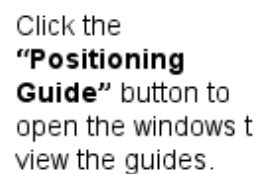

5

Dental Positioning Guides for Canine and Feline

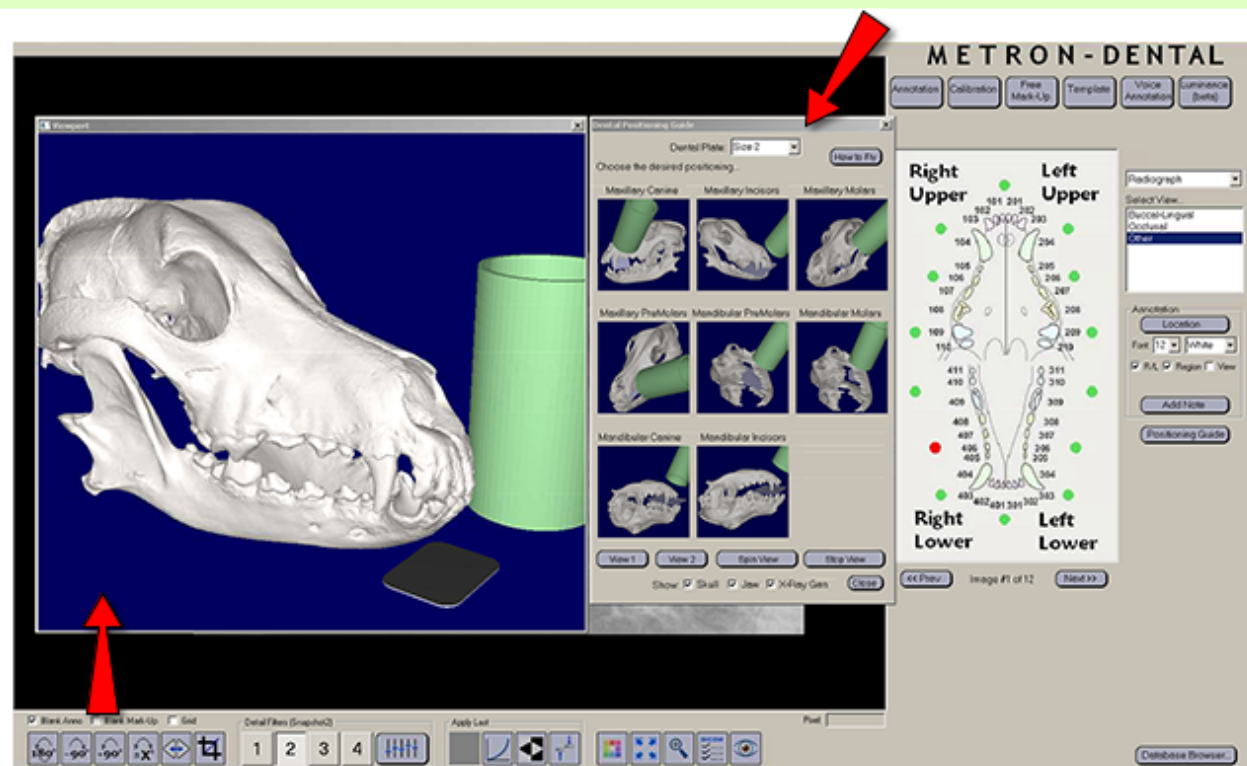


Once you have clicked on the 'Positioning Guide' button, two windows will show on the screen.

METRON
SUPPORT SERVICES

6

Dental Positioning Guides for Canine and Feline

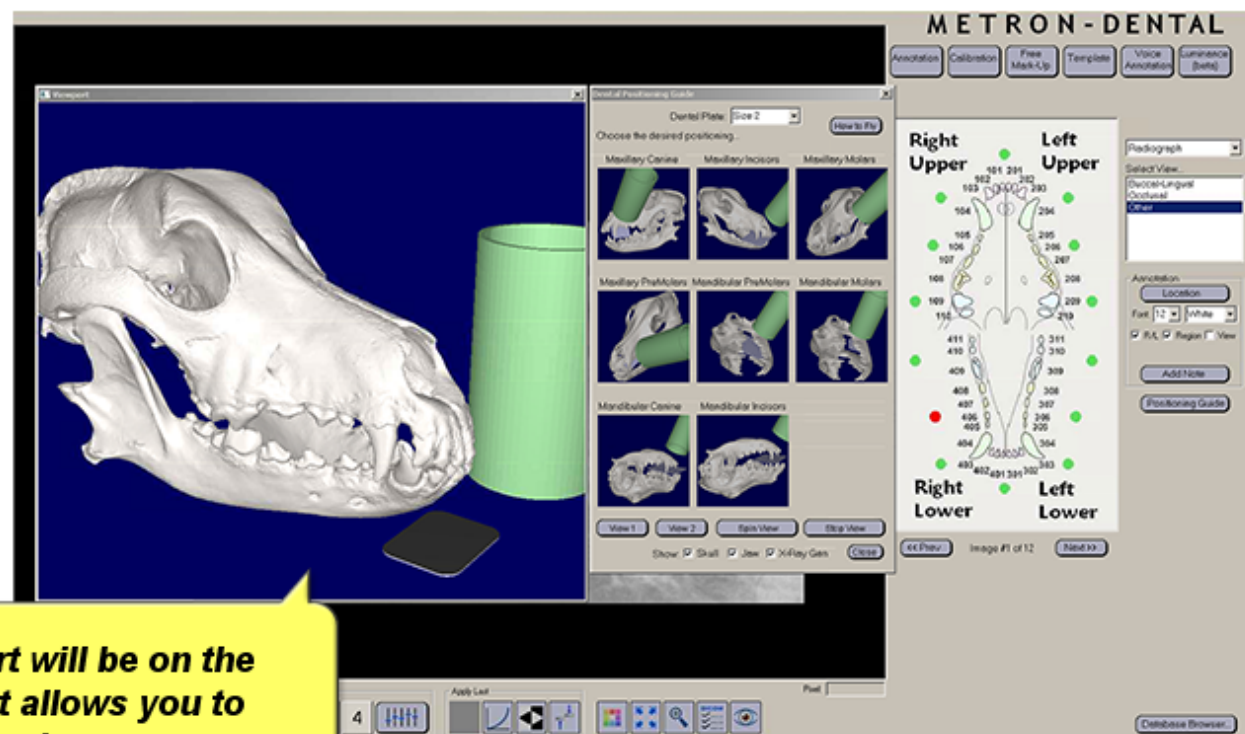


The Viewport will be on the far left and it allows you to view, spin, and move a particular position image.

METRON
SUPPORT SERVICES

7


Dental Positioning Guides for Canine and Feline

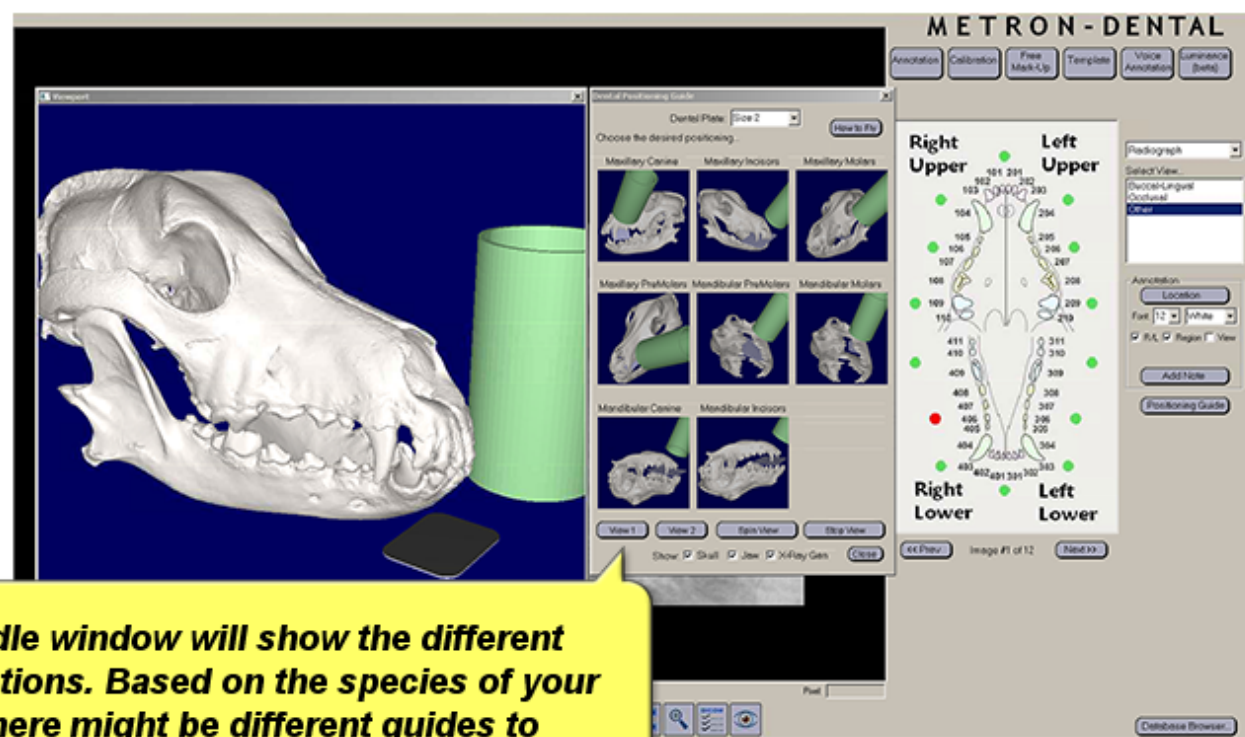

The middle window will show the different guide options. Based on the species of your image, there might be different guides to choose from for viewing.

METRON
SUPPORT SERVICES


Dental Positioning Guides for Canine and Feline

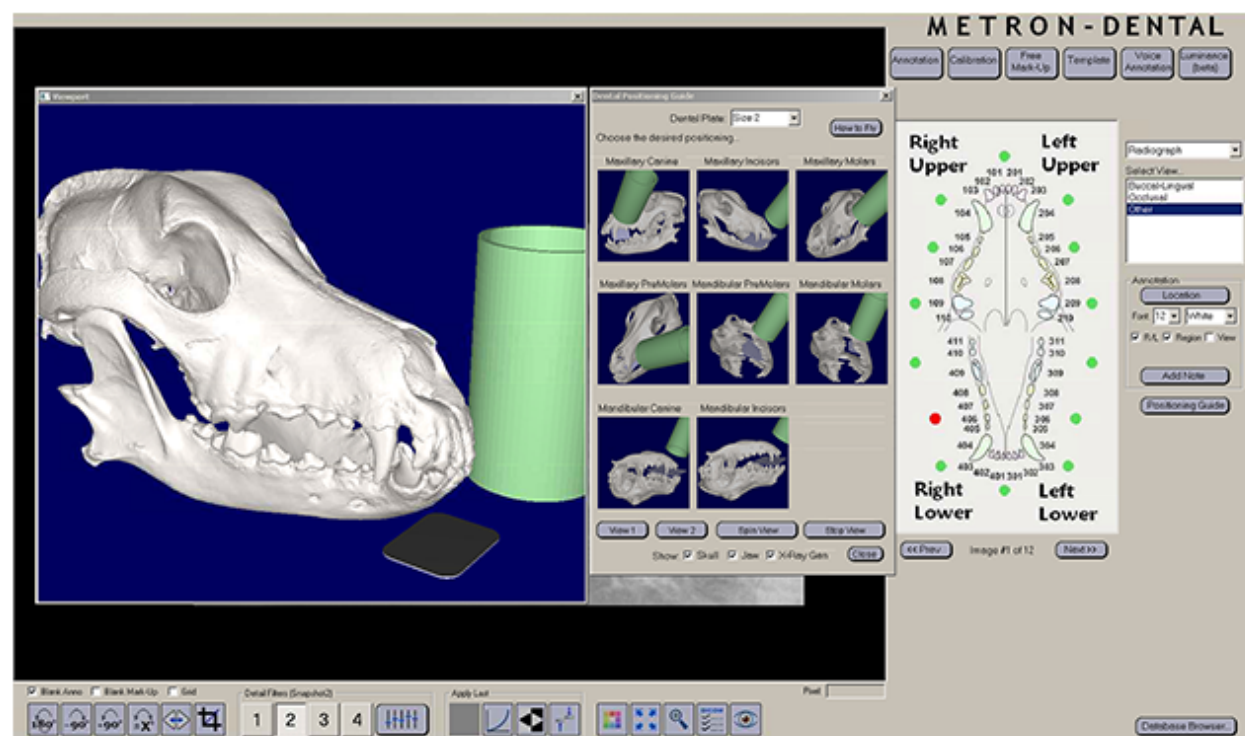


For this training example, we are viewing the Canine Positioning Guides.

METRON
SUPPORT SERVICES

Dental Positioning Guides for Canine and Feline

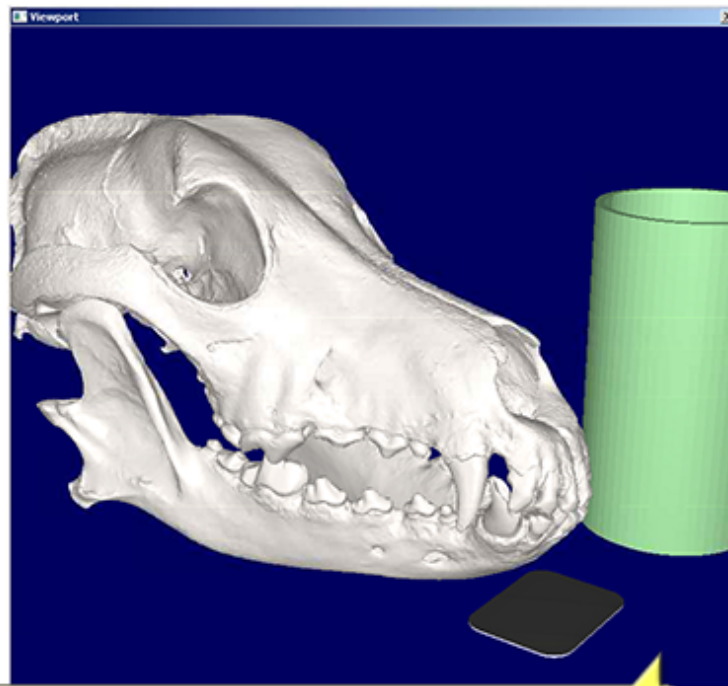

When you first enter the Positioning Guide feature, you will notice the Viewport shows the full image and the X-Ray source in green with no position initially defined.


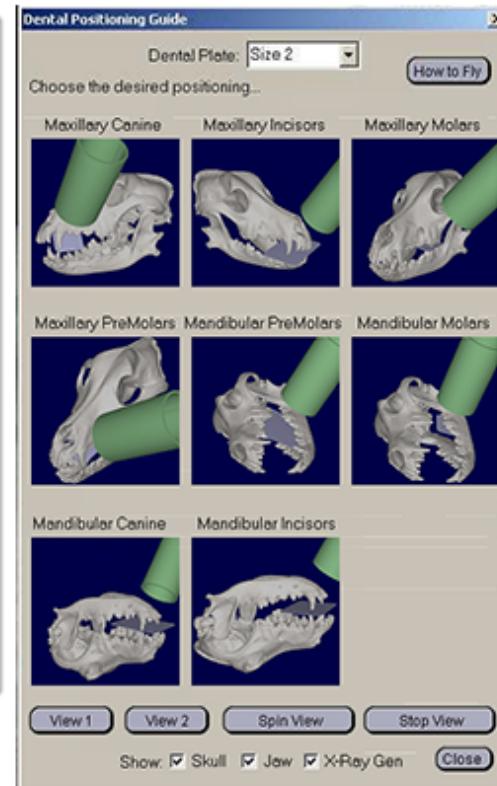


METRON
SUPPORT SERVICES

10

Dental Positioning Guides for Canine and Feline


You are able to define your plate size.

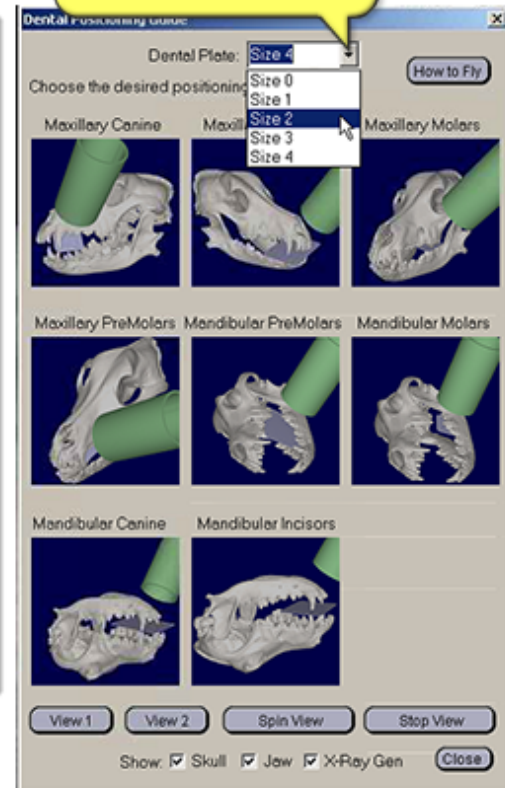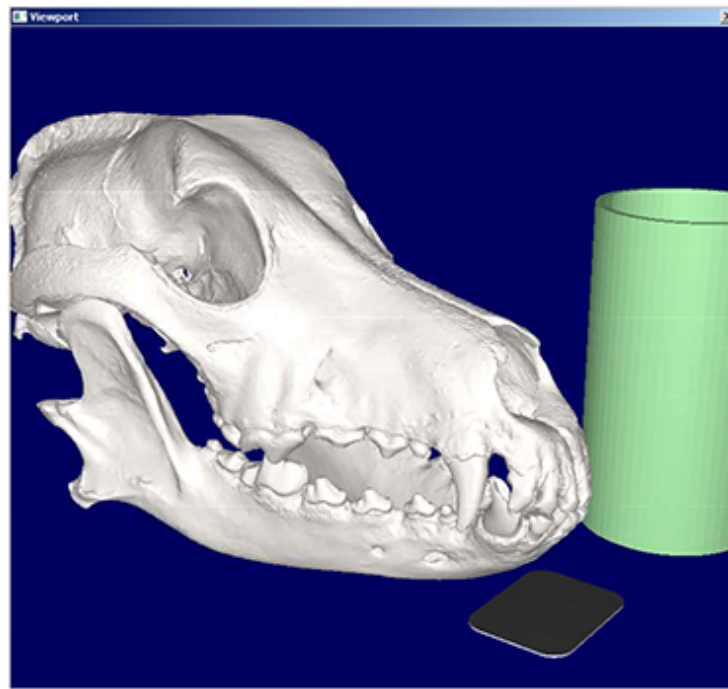

METRON
SUPPORT SERVICES

11


Dental Positioning Guides for Canine and Feline

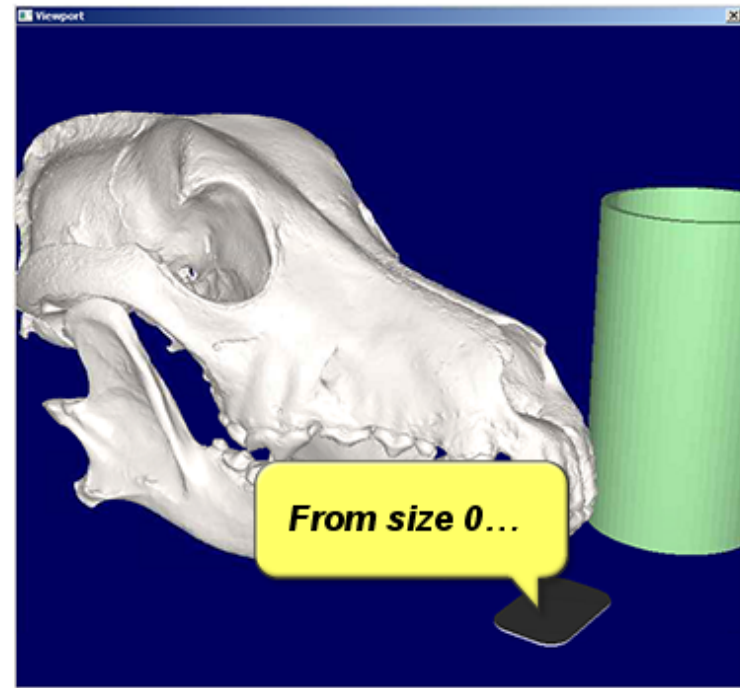

METRON
SUPPORT SERVICES

12


Dental Positioning Guides for Canine and Feline

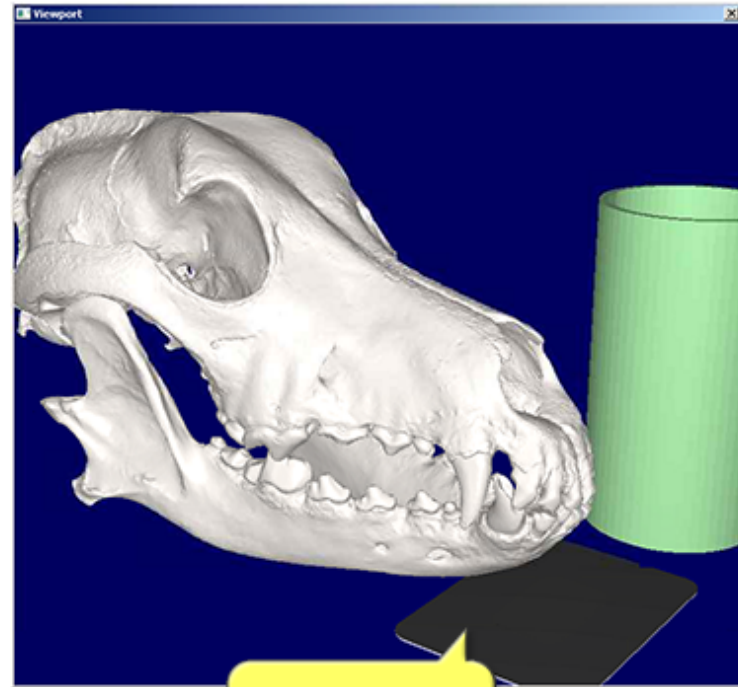


...to size 4.

METRON
SUPPORT SERVICES


Dental Positioning Guides for Canine and Feline

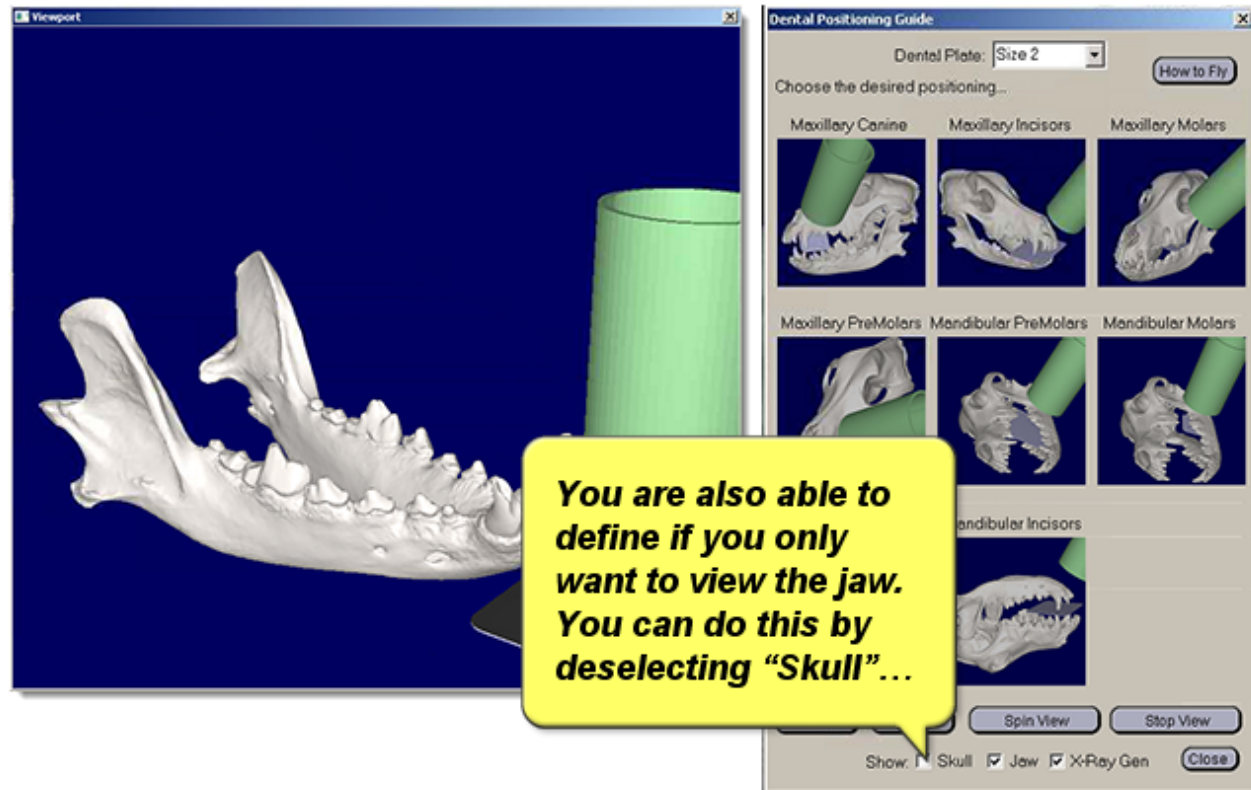

Dental Positioning Guides for Canine and Feline

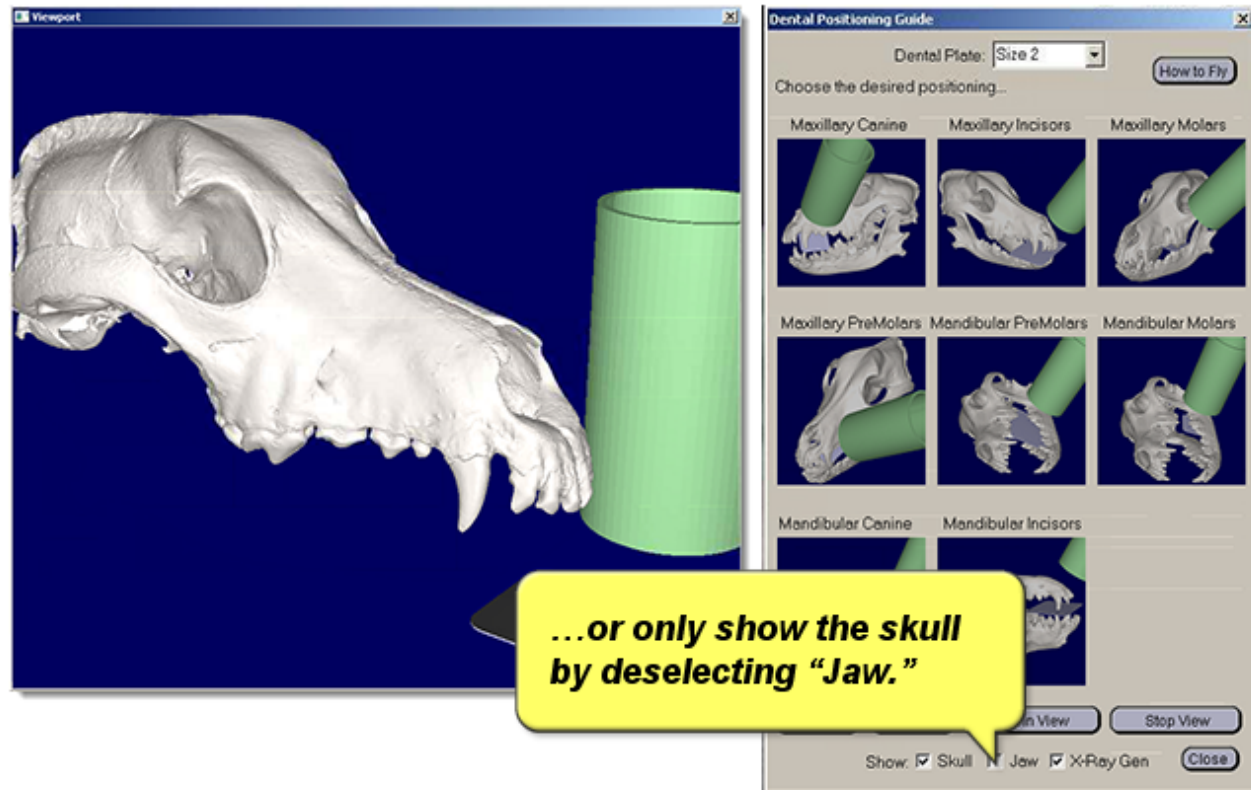


Dental Positioning Guides for Canine and Feline

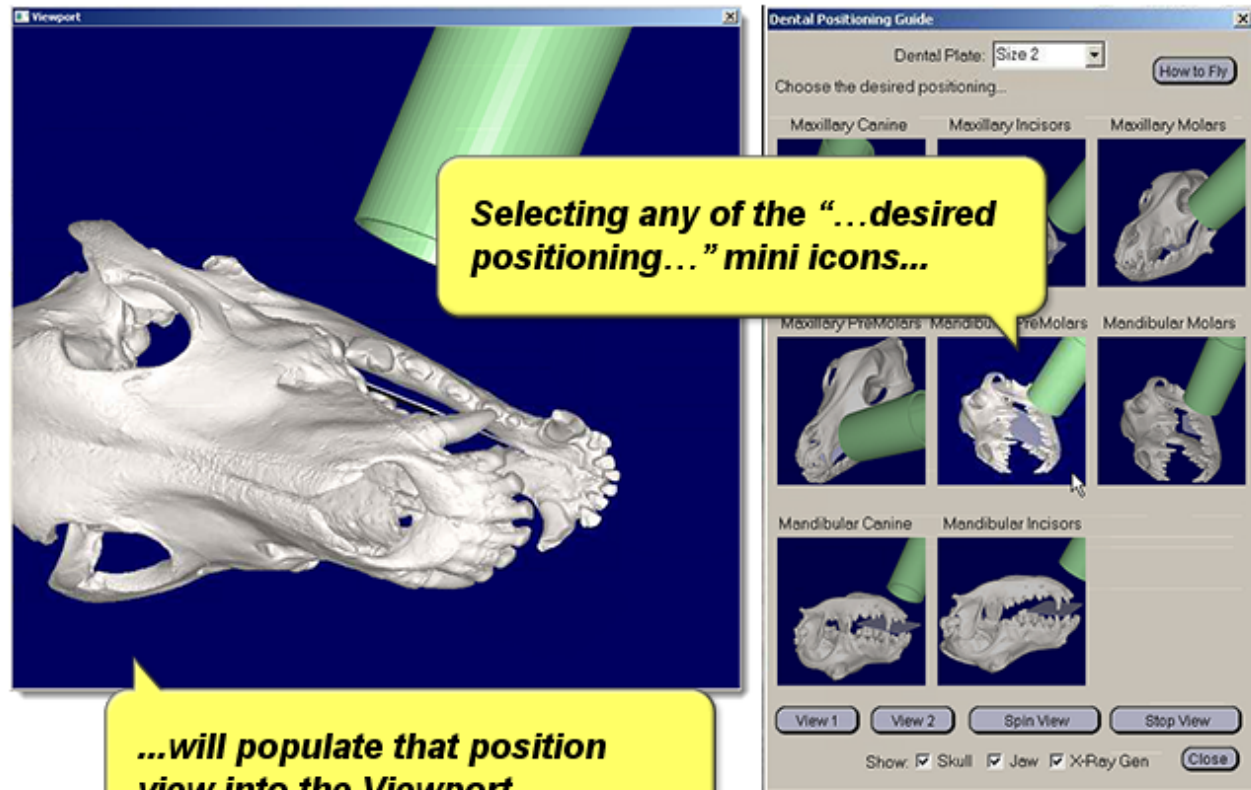


Dental Positioning Guides for Canine and Feline

You may click on any of the mini icons. The position seen in the Viewport window will change based on your selection.

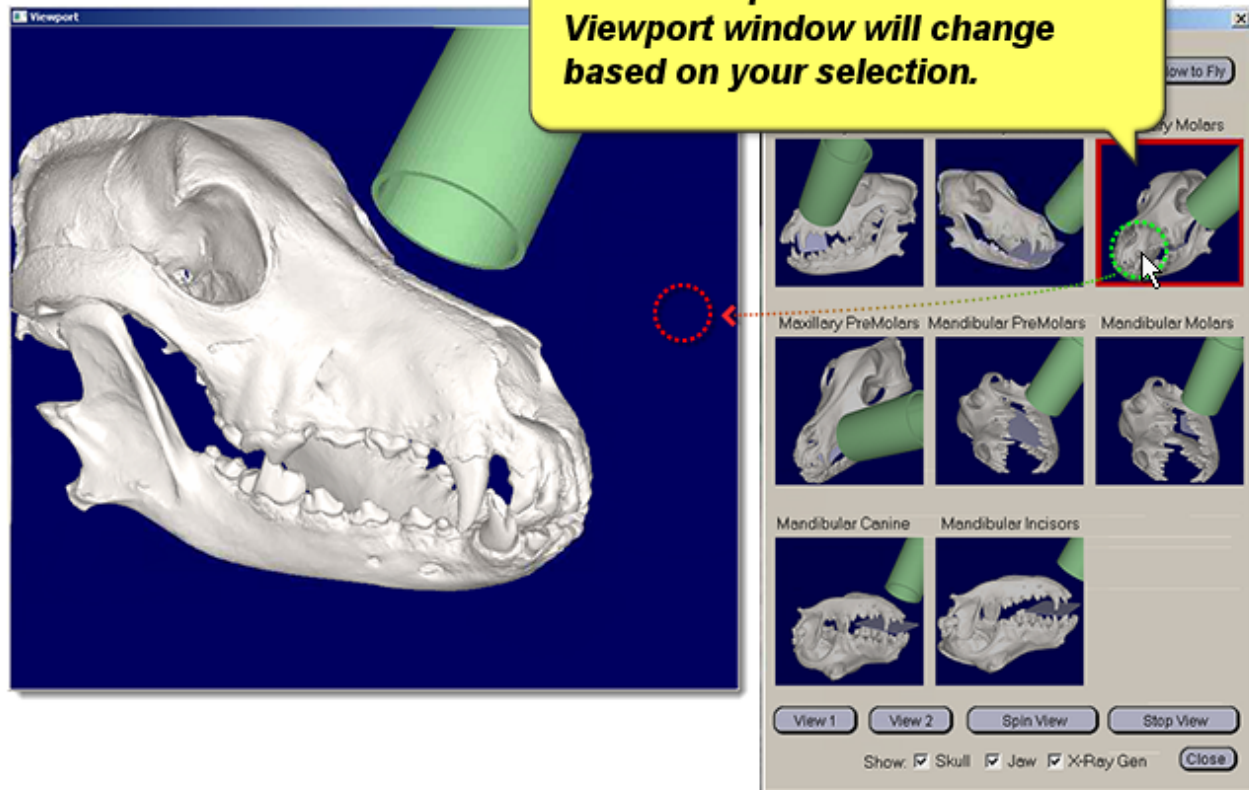


METRON
SUPPORT SERVICES

Dental Positioning Guides for Canine and Feline

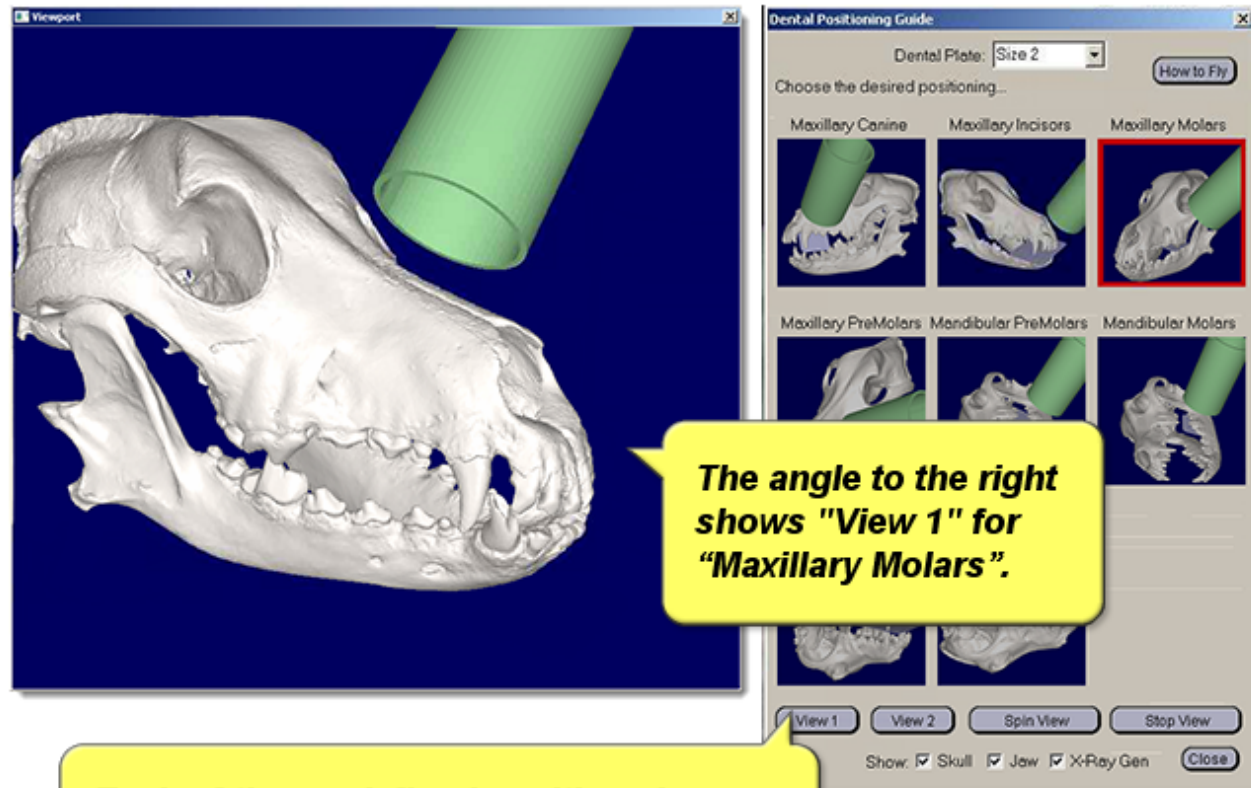

METRON
SUPPORT SERVICES


Dental Positioning Guides for Canine and Feline

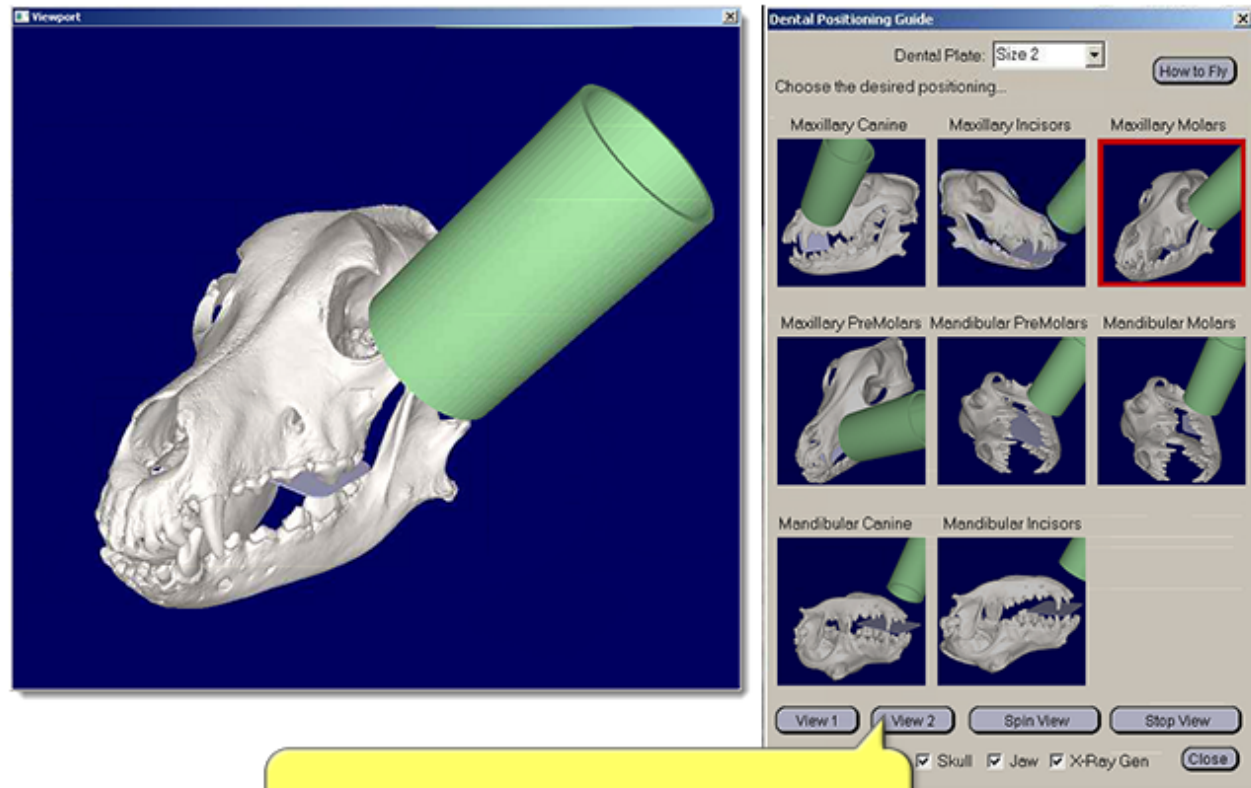


**And this example shows "View 2"
for "Maxillary Molars".**

METRON
SUPPORT SERVICES

Dental Positioning Guides for Canine and Feline

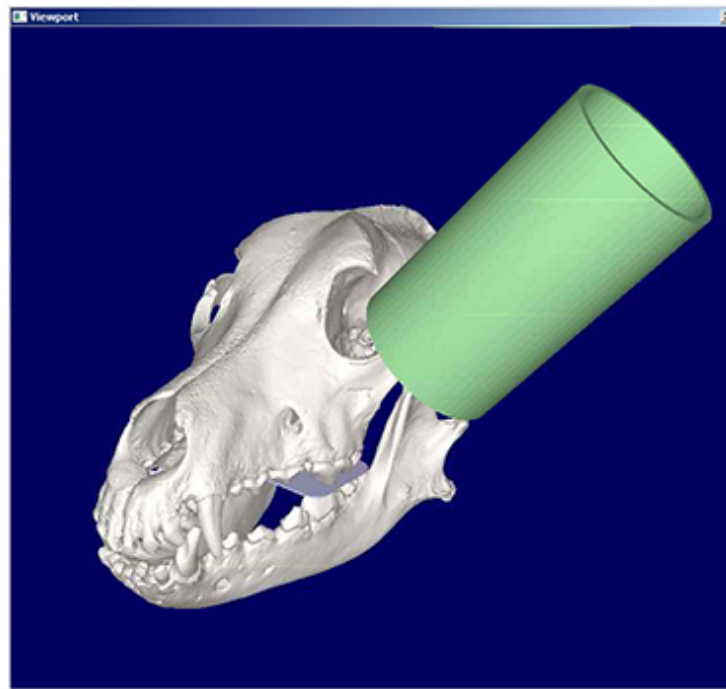


These views are stored for each of the different positions, so depending on which icon is selected these buttons display relevant views for that position.

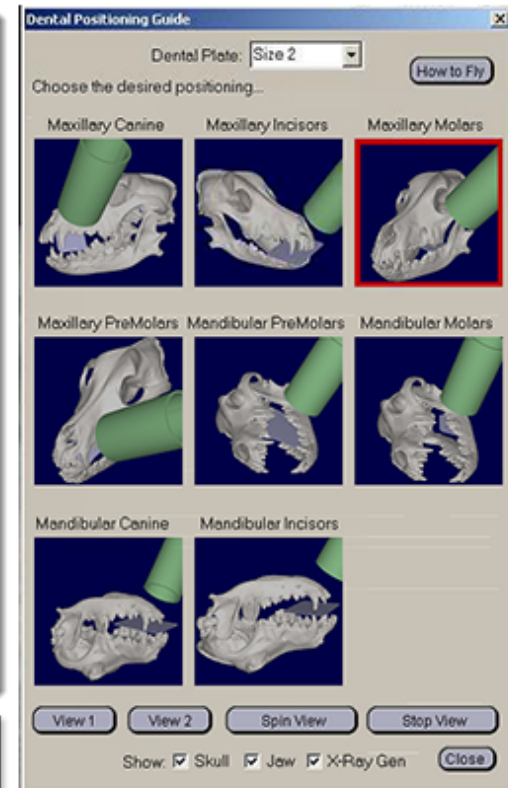

METRON
SUPPORT SERVICES

Dental Positioning for Canine and


To understand how to manually rotate any of the images in the Viewport, click the "How to Fly" button.

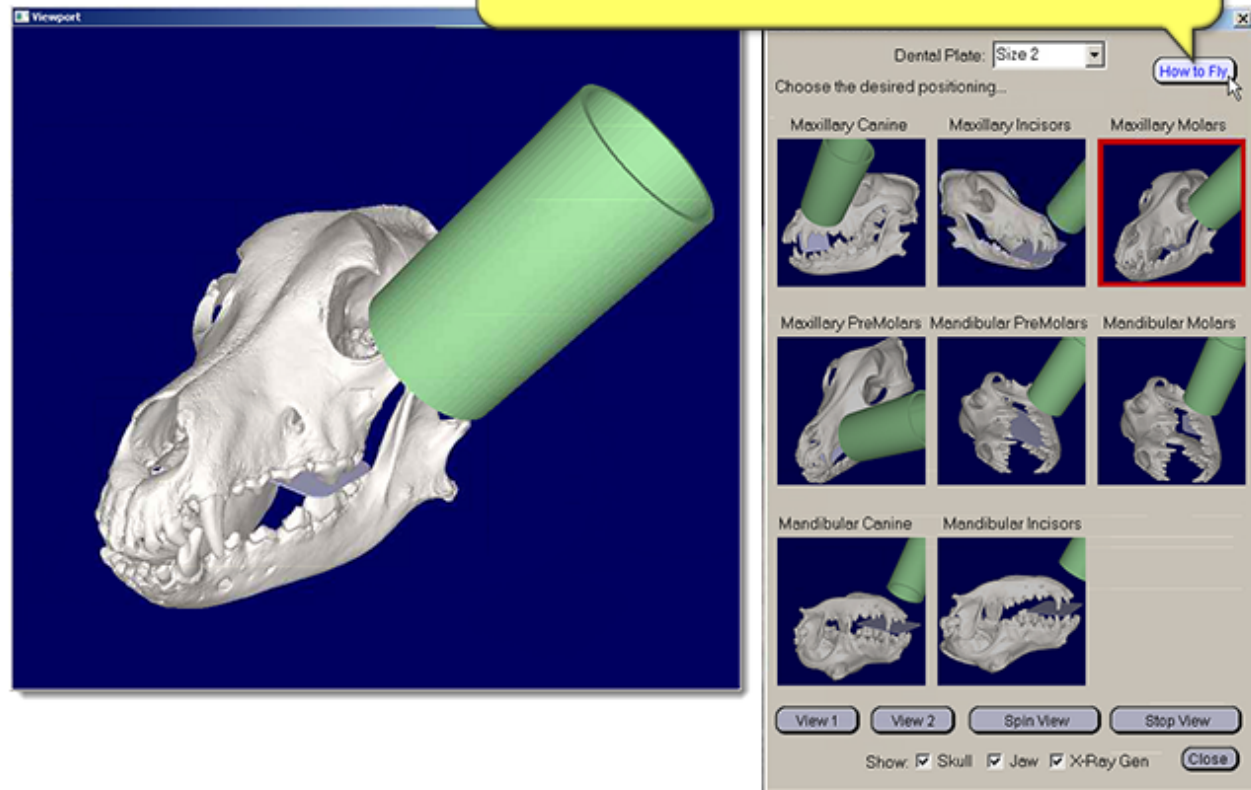

METRON
SUPPORT SERVICES

Dental Positioning Guides for Canine and Feline

This help window shows the functions that are possible using the right and left buttons on your mouse.

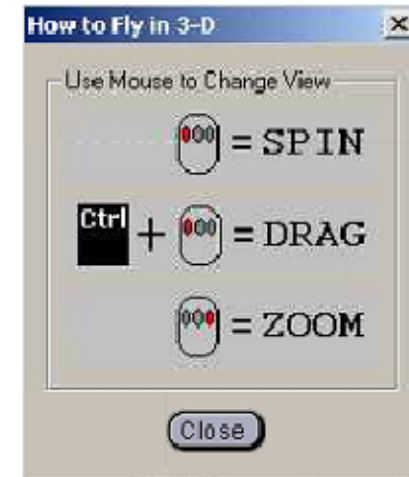

Dental Positioning Guides for Canine and Feline

The video to the right presents an example of a position being moved. First you will see a manual spin followed by the zoom option.

Please Note: Depending on your internet speed, this video may take a few moments to load.




Embedded SWF Movie


fly.swf

To preview this movie you must publish to SWF.

Dental Positioning Guides for Canine and Feline

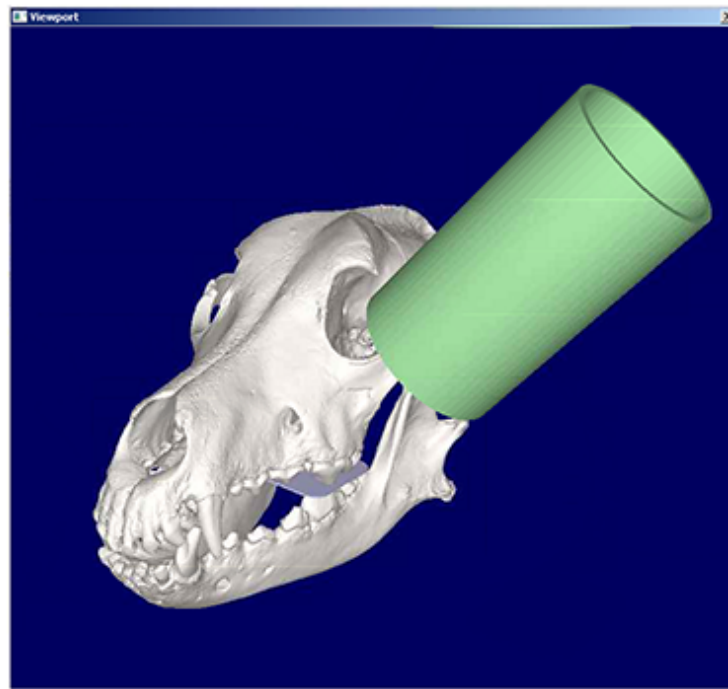

In addition to the mouse buttons, you are able to start and stop the automatic spin using these buttons.

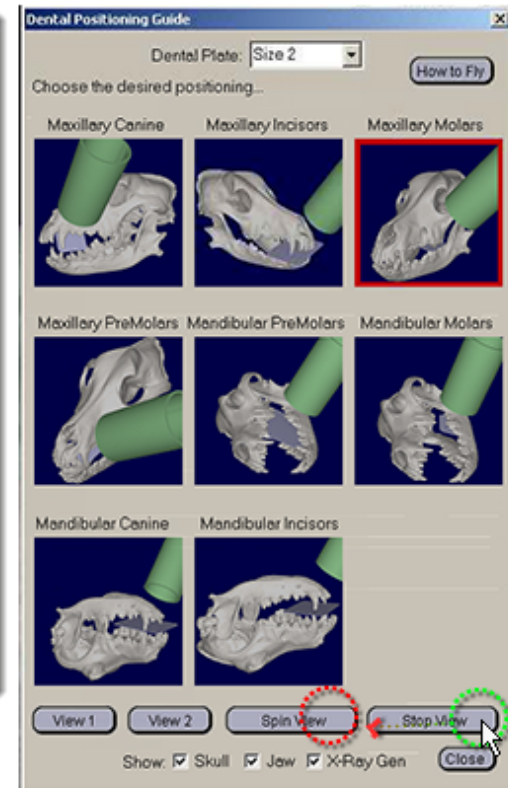

METRON
SUPPORT SERVICES

Dental Positioning Guides for Canine and Feline

Watch the video to view an example of the "Spin" button being used.

Please Note: Depending on your internet speed, this video may take a few moments to load.



Embedded SWF Movie

spin.swf

To preview this movie you must publish to SWF.



Dental Positioning Guides for Canine and Feline

This completes the "Dental Positioning Guides for Canine and Feline" online training.
Please feel free to view it again or select another online training session to review.



Click F11 on your keyboard to minimize the browser screen
and click the window closed to exit this training session.

