



Window Leveling in Metron



Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.



Window Leveling in Metron

This online training session will provide you with the additional details about "Window Leveling in Metron".




Click F11 on your keyboard for full screen display!

Dependent on your screen size, the menu buttons on the Database Browser screen may appear in one or two rows.

1

Window Leveling in Metron

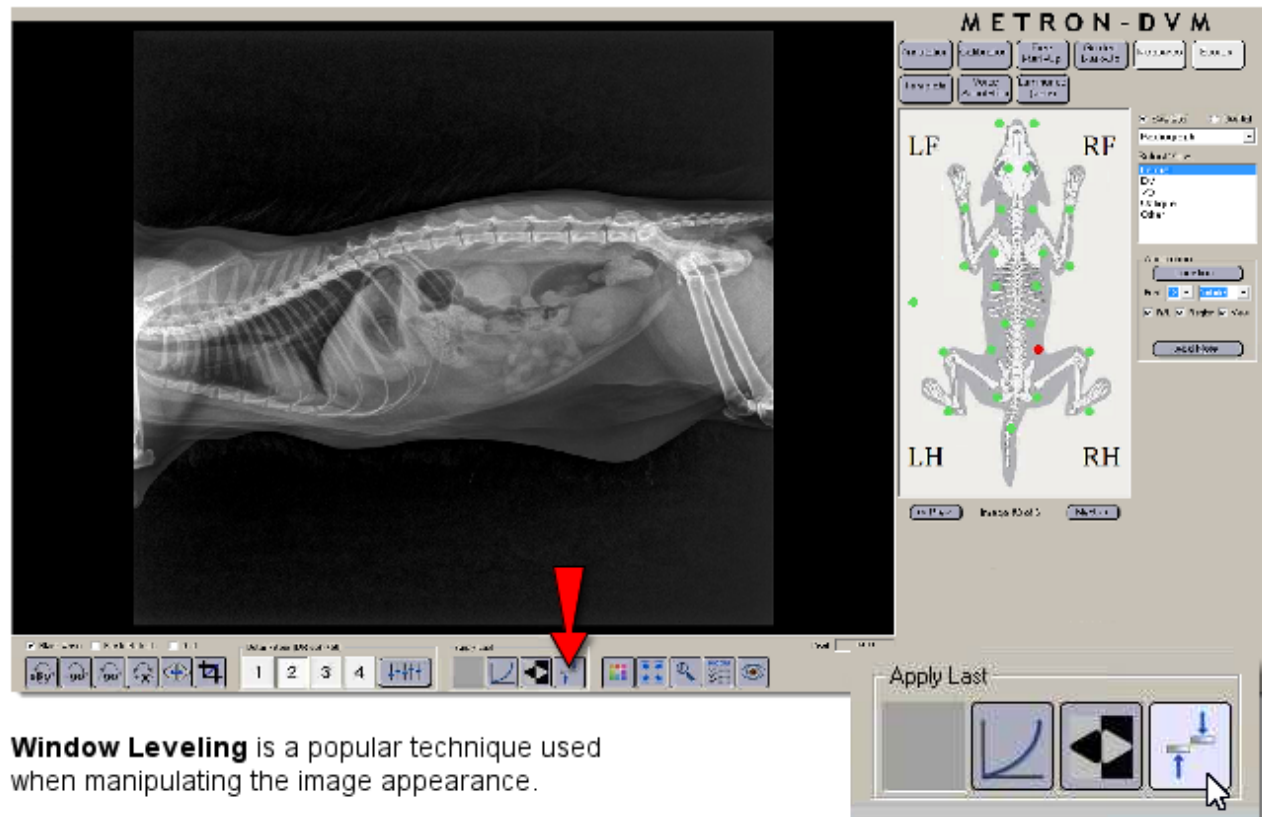


Window Leveling is a popular technique used when manipulating the image appearance.

METRON
SUPPORT SERVICES

2

Window Leveling in Metron

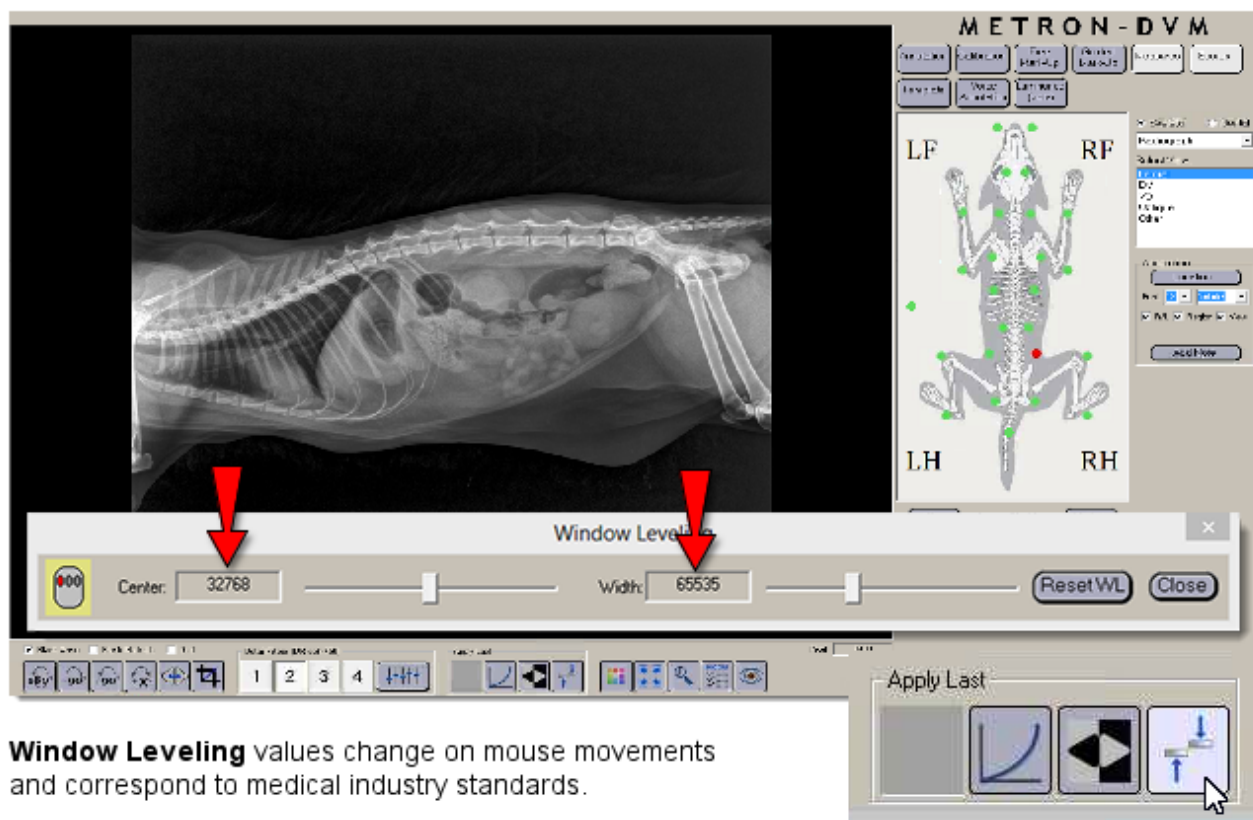


Window Leveling values change on mouse movements and correspond to medical industry standards.

METRON
SUPPORT SERVICES

3

Window Leveling in Metron

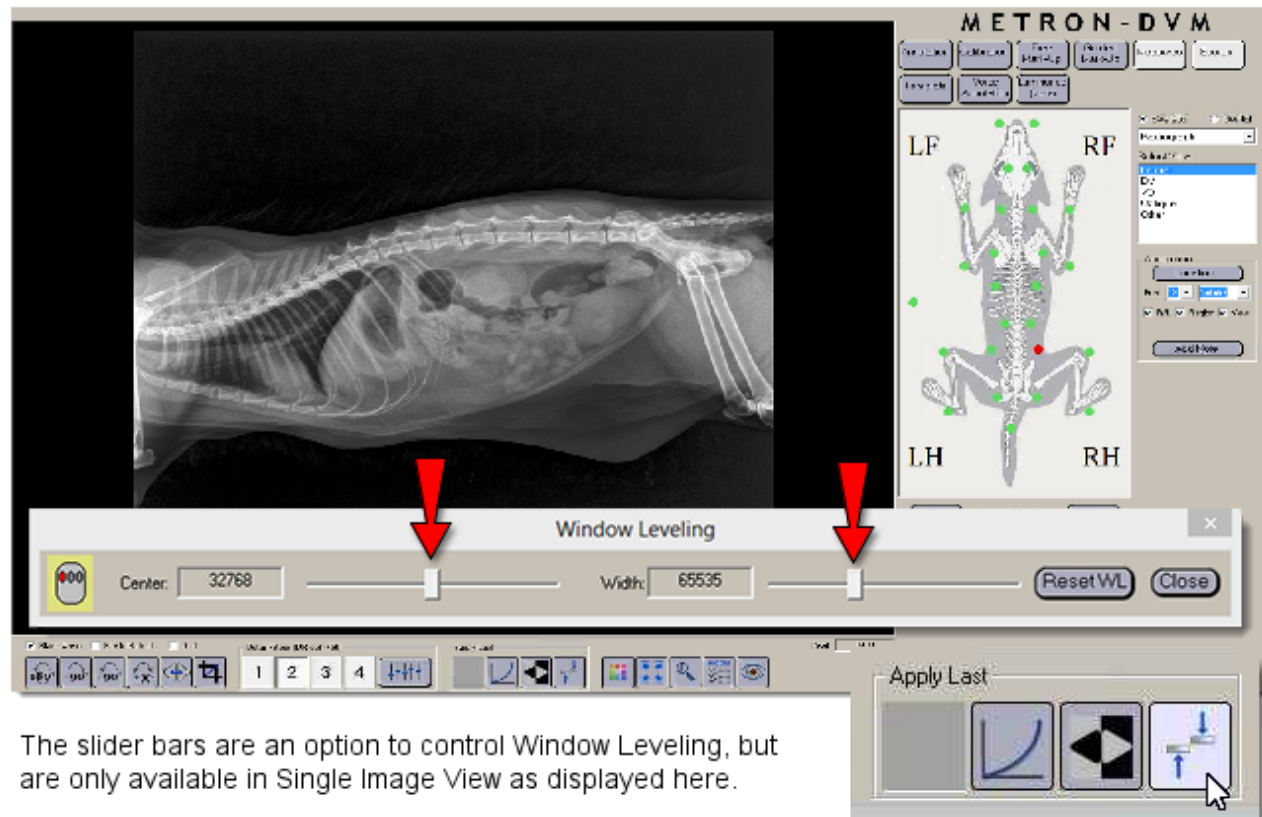

The slider bars are an option to control Window Leveling, but are only available in Single Image View as displayed here.

METRON
SUPPORT SERVICES

4

Window Leveling in Metron



Embedded SWF Movie

WinLev711-f.swf


To preview this movie you must publish to SWF.

In this example we are demonstrating the use of both the active **Window Leveling mouse movement** and **Window Leveling slider bar movement**. Just hold down the left mouse button and move to change the image.

METRON
SUPPORT SERVICES

5

Window Leveling in Metron

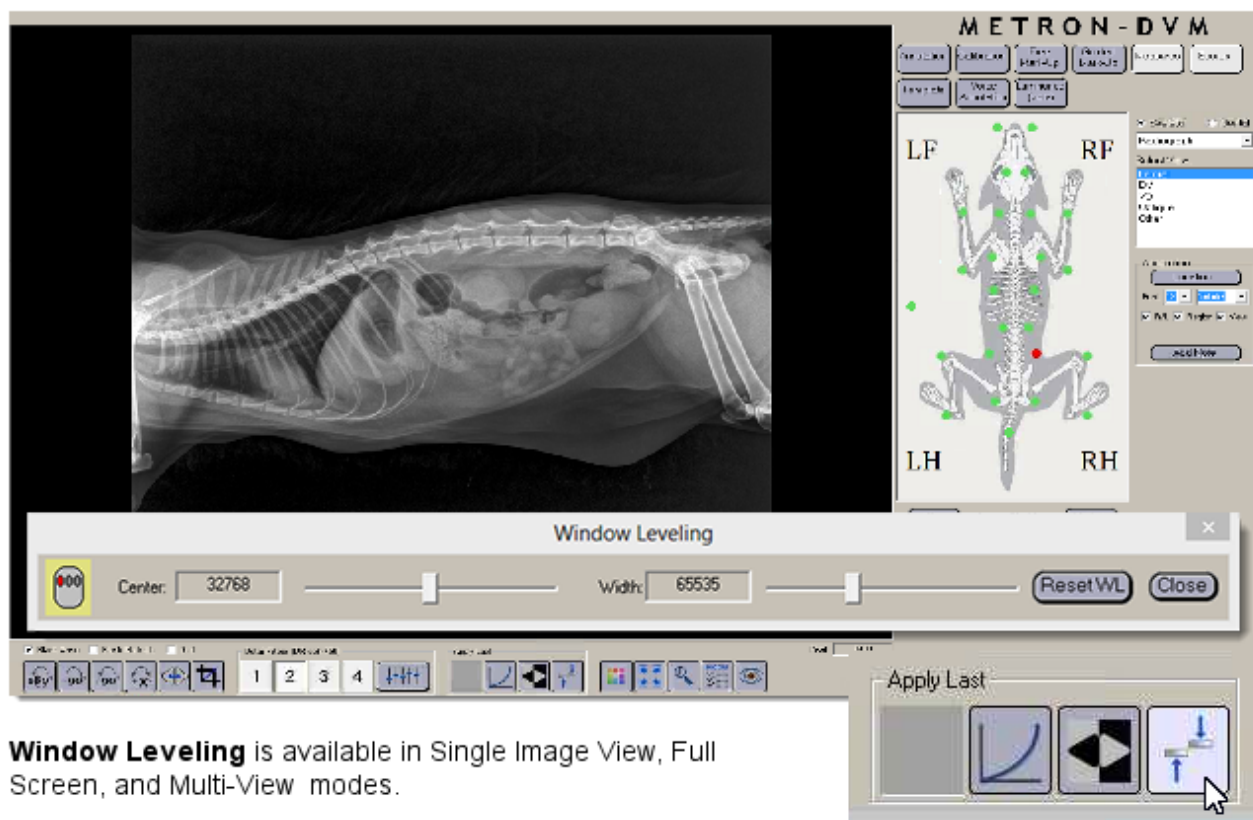

Window Leveling is available in Single Image View, Full Screen, and Multi-View modes.

METRON
SUPPORT SERVICES

6

Window Leveling in Metron


Window Leveling

is available in Full Screen and Multi-View Modes, but does not include the slider bar method as presented on the Single Image View screen.



7

Window Leveling in Metron

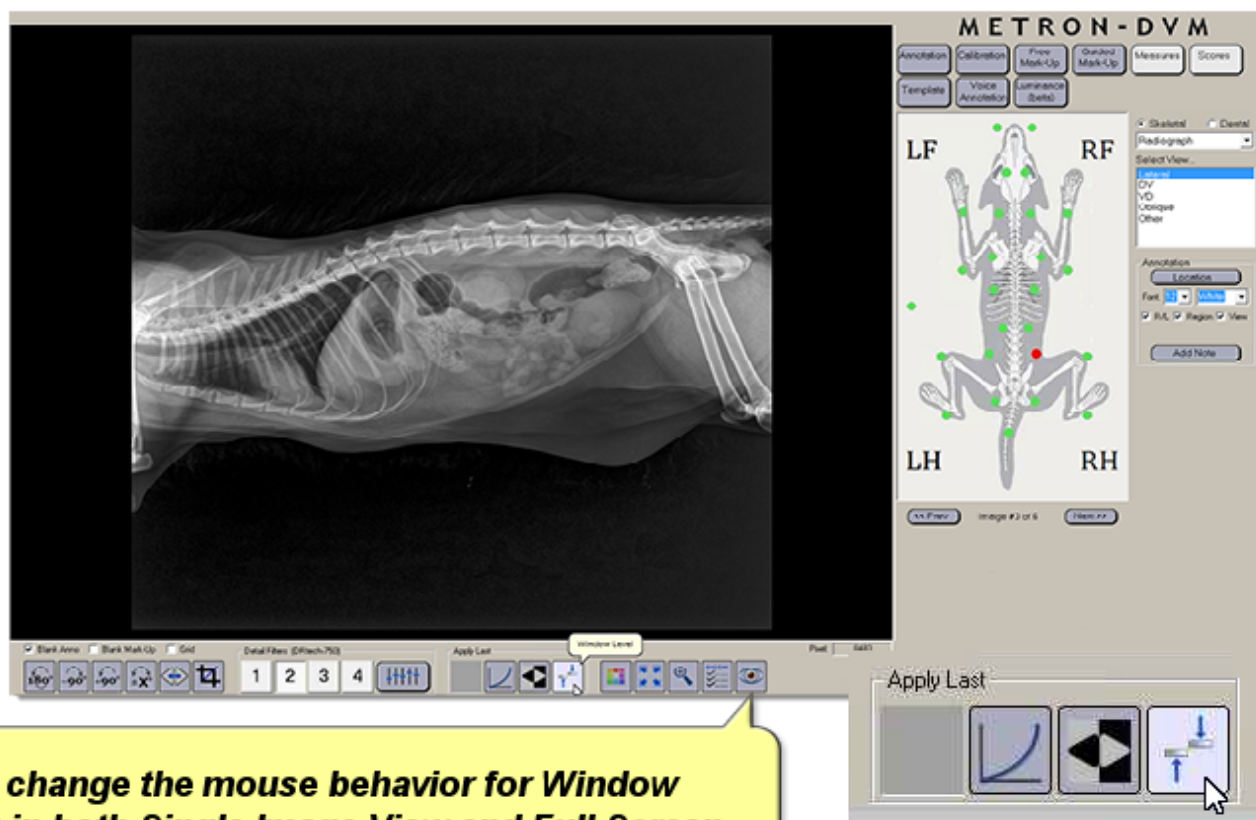


You can change the mouse behavior for Window Leveling in both Single Image View and Full Screen modes by changing the Full Screen Preferences. Select an image, then click the "View" button.

METRON
SUPPORT SERVICES

8


Window Leveling in Metron

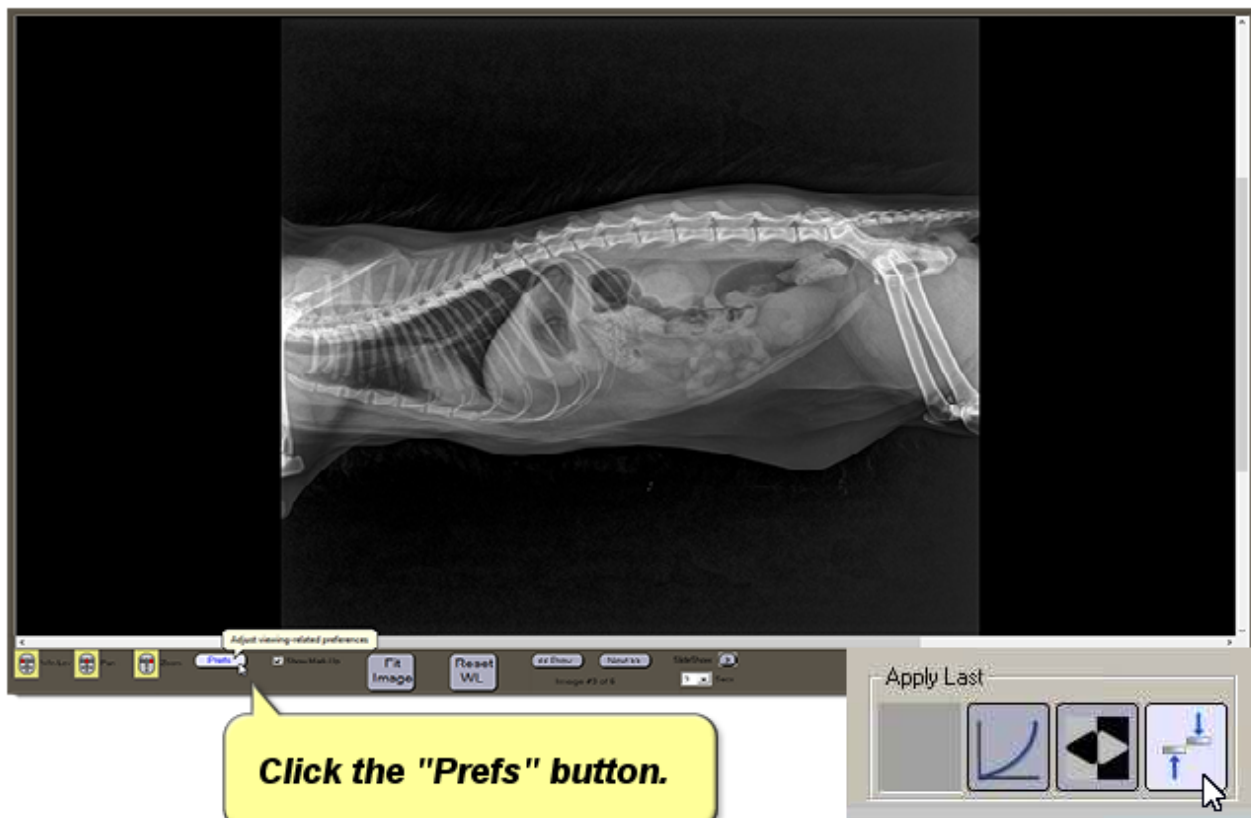


METRON
SUPPORT SERVICES

9

Window Leveling in Metron

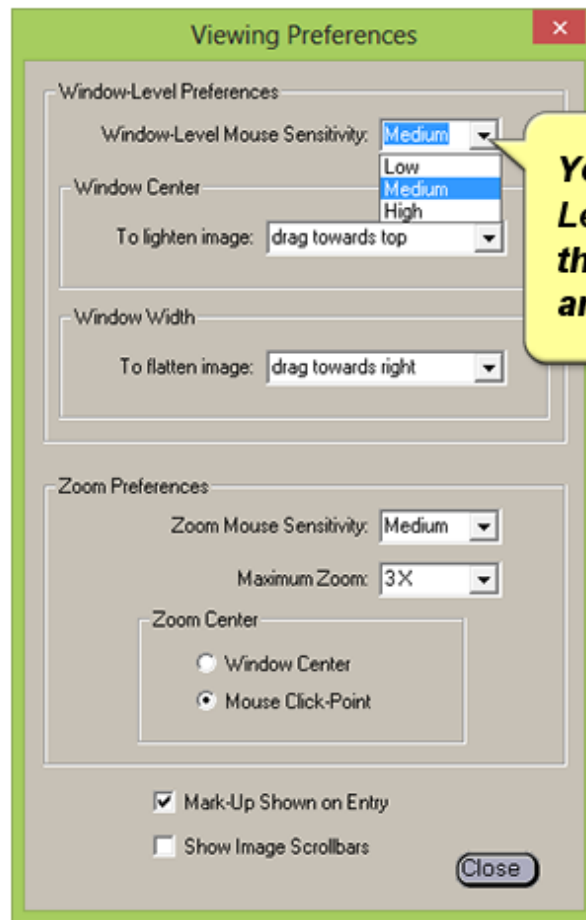


You can change the Window Leveling defaults by using the pull down arrows to view and select choices.

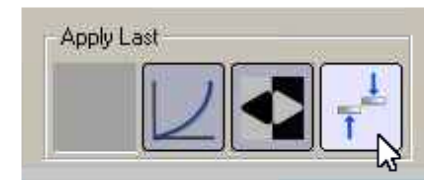

METRON
SUPPORT SERVICES

Window Leveling in Metron

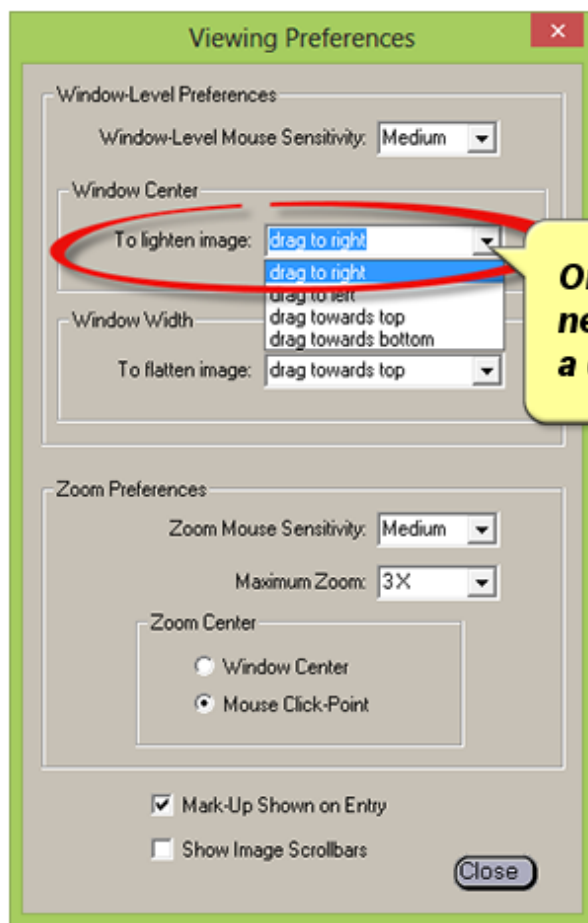


Once you have selected a new option, it will become a default setting.



METRON
SUPPORT SERVICES

Window Leveling in Metron

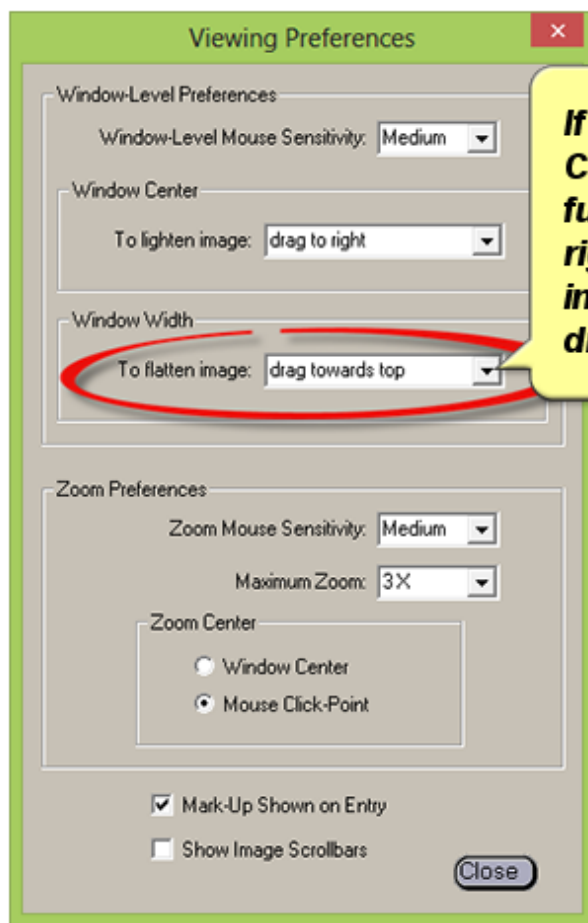


If the user prefers the Window Center (to lighten the image) to function dragging the mouse to the right...the Window Width (to flatten image) will automatically change to dragging the mouse toward the top.



METRON
SUPPORT SERVICES

12

Window Leveling in Metron

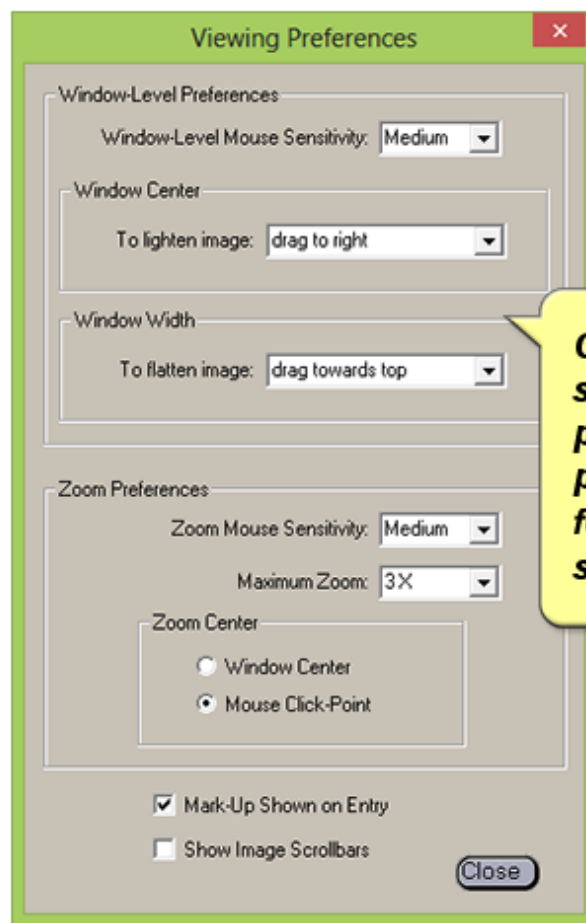


Changes you make to the default settings are based on your personal preferences only.... perhaps influenced by similar features in other imaging software programs.



METRON
SUPPORT SERVICES

13

Window Leveling in Metron

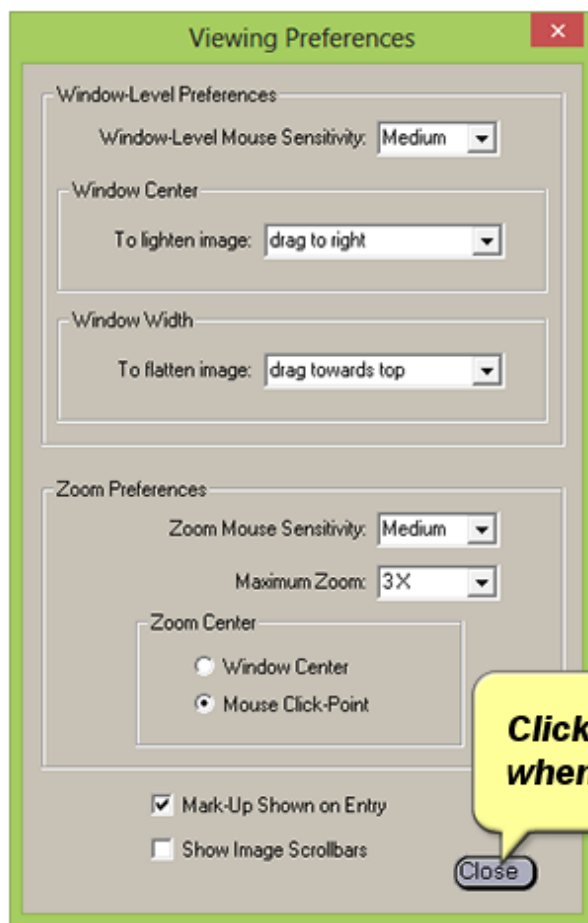

**Click "Close"
when completed.**



METRON
SUPPORT SERVICES



Window Leveling in Metron

This completes the "Window Leveling in Metron" online training.
Please feel free to view it again or select another online training session to review.



Click F11 on your keyboard to minimize the browser screen
and click the window closed to exit this training session.

